HYUNDAI

FM 605B

FM 605P

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Dear customer, thank you for purchasing our product. Please read the operating instructions carefully before putting the appliance into operation and keep these instructions including the warranty, the receipt and, if possible, the box with the internal packing.

GENERAL SAFETY INSTRUCTIONS

- The appliance is designed exclusively for private use and for the envisaged purpose. This appliance is not fit for commercial use.
- Keep it away from sources of heat, direct sunlight, humidity (never dip it into any liquid) and sharp edges. Do not use the appliance with wet hands. If the appliance is humid or wet, unplug it immediately. Do not put it in water.
- Pay careful attention to the following “Special Safety Instructions for this unit”.


SPECIAL SAFETY INSTRUCTIONS FOR THIS UNIT

- Check whether the data on the type label correspond with the voltage in your socket. The power cord plug has to be connected to the properly connected and grounded socket according.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Before replacing accessories, before cleaning or maintenance, turn off the appliance and disconnect it from the mains by pulling the power cord from the power socket!
- Do not leave the appliance plugged in power supply without supervision!
- If the power cord of the appliance is corrupted, it has to be replaced by the manufacturer, its service technician or a similarly qualified person so as to prevent dangerous situations.
- Never use the appliance if its power cord or plug is damaged, if it does not work properly, if it fell down and was damaged and it is leaking. In this case take the appliance to a special service to check its safety and proper function.
- DO NOT attempt to plug in and unplug the unit while feet are in water.
- Before plugging in or unplugging the unit, be sure the control switch is in OFF position [0].
- NEVER immerse the unit completely into water or other liquid for cleaning or other purposes
- DO NOT stand up or placing other objects in the unit, use it only when you are in sitting position.
- DO NOT leave the unit in ON position when the unit are not in use or ready for use.
- DO NOT pull the unit by the cord, or use as a handle.
- CLOSE SUPERVISION is necessary when the unit is used by, or near children, or invalid.
- DO NOT use while sleeping, or if you feel drowsy.
- USE ONLY on a flat level surface indoors, never use outdoors.
- Discontinue use if you feel discomfort, or if any pain or irritation resulted, never use for swollen or inflamed feet, or if feet have any skin eruptions. First, consult your physician.
- The appliance has a hot surface. Persons insensitive to heat must be careful when using the appliance.

- If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or a similarly qualified person in order to avoid a hazard.
- Don't put bath salts or bath oil to water in the foot masager.
- If water leaks from the appliance, the appliance should no longer be used.
- Always check the temperature of the water before placing your feet into the Foot Spa.
- Keep cord away from heated surfaces.
- Never drop or insert any object into any opening.
- Do not operate where aerosol (spray) products are being used
- Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electrical shock or injury.
- Individuals with pacemakers should consult a physician before using the product.
- In case of pregnancy, diabetes or illness consult your doctor before using the Foot Spa.

DESCRIPTION OF THE CONTROLS

- P1 Cover
- P2 Control switch
- P3 Massage balll
- P4 Air nozzle
- P5 Stimus-nodes
- P6 IR reflection
- P7 Massage wheel

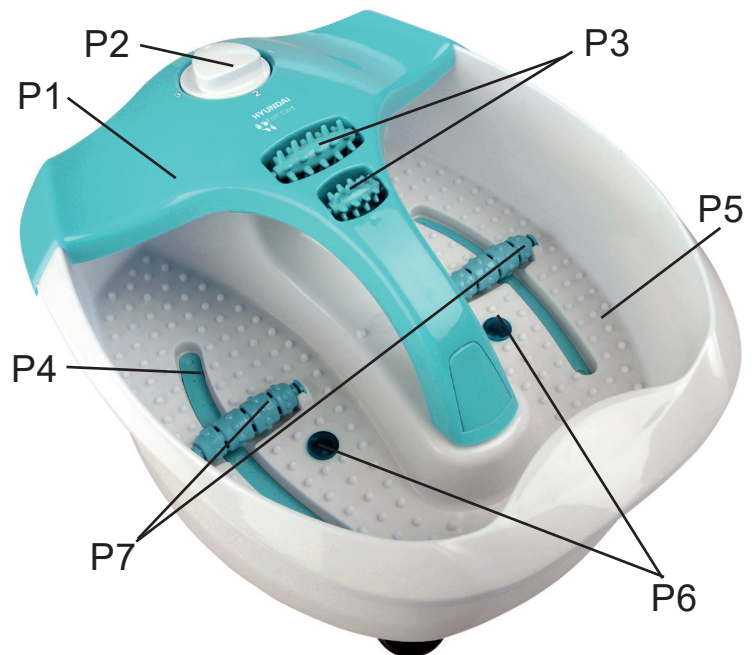INSTRUCTION FOR USE

Proper use produces a measurable increase in blood circulation in the skin of the feet which manifests itself in a pleasant and lasting feeling of warmth (rise in skin temperature). Foot massager baths at regular intervals can produce improved blood circulation in the feet and indirectly in the entire leg area as well.

There are about 200 stimu-nodes act as gentle fingers to make contact with the soles of your feet and give you most effective and invigorating massage possible. You can easily control their effectiveness by simply pressing down lightly or letting up the pressure with your feet. There are 14 air jets provide a relaxation warm bubble bath for your feet.

INSTALATION

Place the appliance on a level surface and the non-skid rubber feet will prevent it from slipping. Be sure the control switch is in "0" position. Then plug the unit into AC outlet (regular household current). Do not stand on the appliance whenever the unit is operated or not. Do not plug the appliance in while your feet are in the water.


Fill with warm water (in desired temperature) into the unit and do not fill water over the "MAX" water-level line on the keel. This unit is designed to maintain the water temperature at a comfortable level, not to raise it. To check the control switch is in "0" position. Then plug the power cord into AC outlet (regular household current).

DO NOT CONNECT OR DISCONNECT THE PLUG FROM THE POWER OUTLET WHILE YOUR FEET ARE IN THE WATER.

Seat yourself comfortably, place feet on the foot rests and select your desired setting. Using the control switch you can select the desired applications:


VIBRATION MASSAGE + IR REFLECTION (1)

Set switch to this mark if you want to massage your feet in cold water. A foot vibration massage without water is also possible with this setting. Dry massage has better result for blood circulation of your feet.

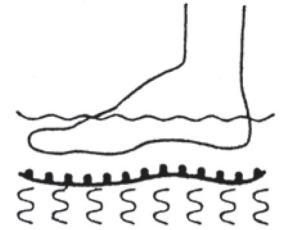BUBBLES + WHIRLPOOL WITH KEEP WARM (2)


Warmth water footbath with bubbling. Ideal for a footbath or for allowing stimulation from previous massage appliances to gradually subside.

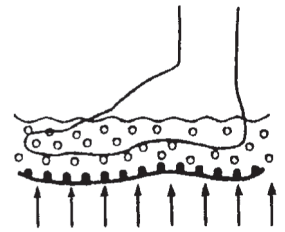VIBRATION MASSAGE + BUBBLES+ WHIRLPOOL WITH KEEP WARM + IR REFLECTION

(3)

Select this setting if you want to have the unit filled with warm water for a foot massage and bubble bath. The built-in heater delays temperature drop of the warm water while bathing. Temperature-controlled full-length heaters under each of foot rests maintain the desired temperature throughout your massager when you set the rotary switch of the unit.

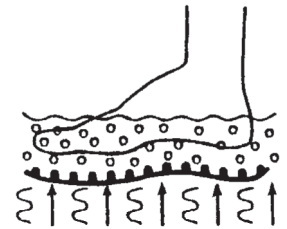OFF (O)

Turn the unit off.

MASSAGE BALL ATTACHMENT

CARE AND MAINTENANCE:

TO EMPTY

First, turn off control switch to the position **0** , then unplug the unit. Tip the unit allowing water to spill out along the drain spout. Never pour water over the control switch.

This is a health care unit, not a washing basin. Clean feet before use and use clean and soft water. Deposit will clog air jets and undermine the bubble effect. Unplug the appliance form the socket, cool the appliance before cleaning.

TO CLEAN: After using, simply rinse out with water. Never immerse the entire unit into water. If you wish, you can wipe off all surfaces with a cloth damped with water and a mild detergent, then rinse and dry. Never use sharp, abrasive or flammable cleanser. This might destroy the unit surface and/or the electric insulation.

TO STORE: Be sure the unit are unplugged, cooled and emptied of water. Then store it in a clean and dry place. Before storage just to wore the whole power cord round the wiring holder. Never hang the unit by the power cord.

TROUBLESHOOTING

The device is not working

- Check that the mains connection cable is firmly plugged in.
- Check the right function of outlet. Test outlet with other electrical appliances.
- Check the supply cord to ensure no damaged is evident.

TECHNICAL SPECIFICATIONS

- Foot massage
 - Over 200 stimulation points
 - Keep-warm function
 - Improves the skin return and blood circulation in the foot-zone
 - 3 possibilities: vibration massage + IR reflection, bubbles + whirlpool with keep warm, vibration massage + bubbles+ whirlpool with keep warm + IR reflection
 - Massage-rollers
 - Color: White-blue(605B), white-purple (605P)
 - Power drain: 80 W
 - Power supply: 230 V ~ 50 Hz
- Input in off mode is 0.00 W.

The manufacturer reserves the right to change the technical specifications and accessories for the respective models.

WARNING: Do not use this product near water, in wet areas to avoid fire or injury of electric current. Always turn off the product when you don't use it or before a revision. There aren't any parts in this appliance which are reparable by consumer. Always appeal to a qualified authorized service. The product is under a dangerous tention.

Old electrical appliances, used batteries and accumulators disposal

This symbol on the product or on its packaging indicates that this product (or battery/accumulator) shall not be treated as household waste. Instead it shall be handed over to the applicable collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. The recycling of materials will help to conserve natural resources. For

more detailed information about recycling of this product, please contact your local Civic Office, your household waste disposal service or the shop where you purchased the product.

WARNINGS AND SYMBOLS USED ON THE APPLIANCE, PACKAGING OR IN THE INSTRUCTIONS MANUAL:

HOUSEHOLD USE ONLY; DO NOT IMMERSE IN WATER OR OTHER LIQUIDS;

TO AVOID DANGER OF SUFFOCATION, KEEP THIS PLASTIC BAG AWAY FROM BABIES AND CHILDREN. DO NOT USE THIS BAG IN CRIBS, BEDS, CARRIAGES OR PLAYPENS. THIS BAG IS NOT A TOY.

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Licensed by HYUNDAI Corporation Holdings, Korea