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U2566A

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A00

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2026/03/07

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
W30

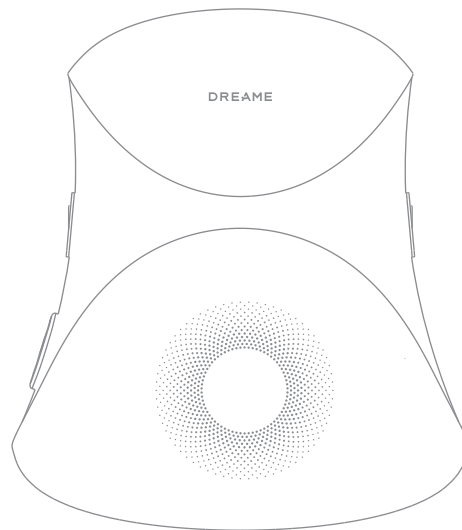Website: <https://global.dreametech.com>

Manufactured by: Dreame Trading (Tianjin) Co.,Ltd.

Made in China

A00

Product image for illustration purposes only.
Actual appearance may vary.

DREAME

THANK YOU

FOR PURCHASING THIS DREAME

Lumbar Massage Cushion

.....

The high-quality product is engineered for optimal performance.
If you have any questions on the product,
please contact us: [https://global.dreametech.com/pages/
contact-us](https://global.dreametech.com/pages/contact-us)

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Scan the QR code for the digital manual.

Safety Instructions

Please read the user manual carefully before using the appliance and keep it for future reference.

When using this product, basic safety precautions should always be followed, including the following:

Usage Restrictions & Precautions:

- This appliance is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- The appliance must only be supplied at safety extra-low voltage corresponding to the marking on the appliance.
- The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- This appliance contains batteries that are non-replaceable.
- When a power cord is required, please use one that complies with national safety standards.
- Always disconnect the power supply when leaving the appliance unattended.
- If you notice any abnormal noise, odor, excessive temperature, or malfunction, stop using the appliance immediately, turn off the switch, and disconnect the power supply.
- Do not connect electrical appliances to light sockets.
- Always ensure that the combined current load of all appliances plugged into the power cord does not exceed its rated capacity.
- Do not use this product in places where explosive or flammable materials may be present or similar hazardous environments.
- Users must not disassemble the appliance or modify internal wiring, components, or protective devices.
- Only use this appliance according to the instructions in this manual. Only use accessories recommended by the manufacturer.
- Use the product in a dry indoor environment, and avoid exposure to direct sunlight, rain, or moisture.
- Do not use the product in flammable or explosive areas, such as gas stations, high-dust workshops, or gas plants, locations subject to vibration or impact, or environments with high temperature and humidity.
- Do not use the product while driving or engaging in other high-risk activities.

Safety Instructions

- Do not apply the massager near the heart, to the head, or on areas with skin disorders.
- This product features a heat function. Prolonged use can cause low-temperature burns. Do not use while sleeping or in an unconscious state.
- For users with skin sensory impairment or sensitive skin, exercise caution when using the heat and massage functions.
- Individuals with cardiovascular or cerebrovascular diseases, skin inflammation, swelling, injuries, or other neck-related conditions should consult a physician before use.
- This product is intended for personal relaxation only and is not a substitute for medical devices.
- Pregnant women and new mothers should use with caution and only under the advice and guidance of a doctor.
- This product contains lithium batteries. Do not disassemble or modify batteries, immerse batteries in water or throw batteries into fire.

The following groups are prohibited from using this product:


- Individuals with major medical conditions or those unable to care for themselves.
- Individuals with severe cognitive disorders.
- Individuals with contractual deformities.
- Individuals in the spastic phase following a cerebrovascular accident.
- Infants, children, and individuals who are fully or partially unconscious.
- Individuals suffering from spastic paralysis.
- Do not use this product in conjunction with the following medical electronic devices to avoid serious accidents:
 - Implantable medical devices such as cardiac pacemakers.
 - Life-sustaining medical electronic equipment like artificial heart-lung machines.

Attention!


- Disassembly and repairs should be performed by qualified persons to avoid any product damages or hazards.
- Please contact Customer Service promptly if you encounter any product issues.

Instructions for Use

Main Unit

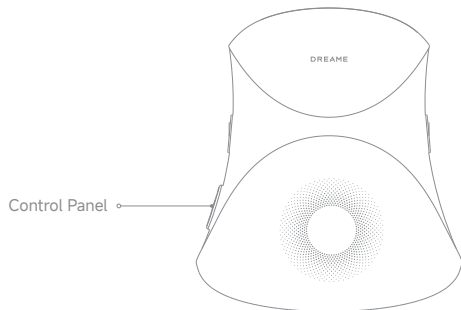Control Panel

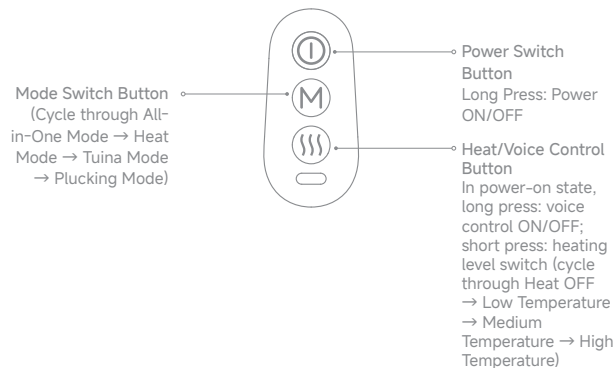Instructions for Use

Indicator Light Descriptions

Description	LED light status
Power button in on state	The indicator stays solid green if the battery is above 25% and will turn solid orange when the battery falls below 25%. The orange indicator will flash when the battery falls below 10%, and the device will automatically shut down after a voice prompt.
Power button in off state	The indicator remains off. The light will not respond to short button presses while the device is in the power-off state.
Short-circuit/Overheat protection	The orange indicator will flash 5 times and the device will automatically shut down.
Stall protection	The orange indicator will flash. If the stall is cleared, the indicator will return to the normal battery status display. If the stall persists, the device will automatically shut down.
Heat button	The red indicator cycles through "OFF → Low Brightness → Medium Brightness → High Brightness" corresponding to Heat OFF → Low Temperature → Medium Temperature → High Temperature settings.
Charging status	1. The indicator breathes orange constantly when battery level increases from 0% to 100% during charging. 2. The indicator breathes green constantly when the device is fully charged.
Massage bulbs	The brightness of the massage bulbs will automatically adjust to match the selected heating level.
Mode indicator	The blue indicator remains solid during operation. When the mode button is short-pressed, the blue light will flash twice.

Instructions for Use

Function Descriptions

Function	Mode	Function Descriptions
Voice reminder	Power on	Long press the power switch button for 1.5 seconds until the indicator lights up. The device will announce "Dreame moment, enjoy your journey."
	Power off	Long press the power switch button for 1.5 seconds. The device will announce "Massage finished, powering off..." before shutting down and turning off the indicator.
	Mode switch	The device will report the modes as you cycle through "All-in-One Mode → Heat Mode → Tuina Mode → Plucking Mode". The indicator flashes twice when switching.
	Heat adjustment	The device will report the heat levels as you cycle through "Heat OFF → Low Temperature → Medium Temperature → High Temperature". The indicator changes during switching.
	Voice control	Long press the button for 1.5 seconds in Voice Mode to turn the voice prompts on or off. You can hear the voice prompt "Voice OFF" or "Voice ON".
	Low battery alerts	When the battery falls below 25%, the device will remind "Low battery, please charge.", and the indicator turns solid orange. When the battery falls below 10%, the indicator will flash orange. The device will automatically shut down after the voice prompt: "Battery critically low, powering off..."
	Stall alert	The indicator flashes, accompanied by the voice prompt "Stuck detected. Please release."
	Overheat protection	The mode indicator flashes, accompanied by the voice prompt "Device overheating. Please try again later. Powering off now. Thank you."
Factory reset	/	While the device is powered on, long press the Heat button to restore factory settings.

Common Issues

If the product does not function normally, please refer to the table below to resolve common issues.

Common Issues	Possible Cause	Solution
Product fails to power on	Low battery level	Please charge the device using a 5V/2A power adapter.
	Low battery level; power supply not connected during charging	Check if the charging indicator is breathing.
	Incorrect startup operation	Ensure you long-press the power button for at least 1.5 seconds.
Product shuts down during use	Stall protection triggered	A voice prompt will notify you when a stall occurs. Once the device shuts down, simply restart it to resume use.
Product shuts down after repeated use while charging in low battery status	Insufficient power supply during high-load operation	If a shutdown occurs when the product is used while charging, it is recommended to wait 30 minutes before using the device again.
	Insufficient charger power	Ensure you are using a charger with a 5V output and at least 2A current.

Technical Specifications

Product Name	Lumbar Massage Cushion W30	Rated Charging Voltage	5V 
Product Model	SVP12A	Rated Input Power	10W
Battery Capacity	2200mAh/16.28Wh		

Attention!

The device should be charged using a 5V power supply with a 2A output (where the measured maximum output power is less than 15W). The power source must also comply with safety extra-low voltage standards.

Instructions for Use

Function Descriptions.

Function	Mode	Function Descriptions
Each massage cycle lasts for 10 minutes	All-in-One Mode	This is the default mode for first-time use. The device features a Memory Function, which automatically restores the settings from your previous session upon powering on.
	Heat Mode	The massage nodes will reset to their home position and the massage function will turn off. Only the heat function remains active, with the default temperature set to Medium.
	Tuina Mode	The massage nodes rotate counter-clockwise.
	Plucking Mode	The massage nodes rotate clockwise.
Heat Function	Heat OFF	The heat function is disabled. Both the mode indicator and massage bulbs are turned off.
	Low Temperature	The mode indicator and massage bulbs are set to low brightness. The heating film operates at a low temperature.
	Medium Temperature	The mode indicator and massage bulbs are set to medium brightness. The heating film operates at a medium temperature.
	High Temperature	The mode indicator and massage bulbs are set to high brightness. The heating film operates at a high temperature.
Use While Charging	/	- When the battery is above 10%, the indicator will breathe orange as you use the device while it is charging. All other operations remain the same to the non-charging mode. - It is recommended to connect the Type-C cable before powering on the device. Once the Type-C cable is connected, the indicator will breathe orange. All other operations remain the same to the non-charging mode.

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