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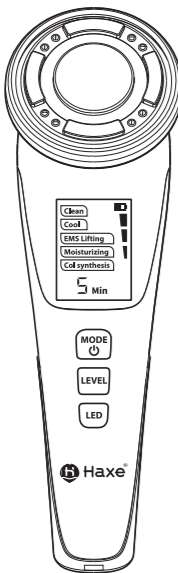
### TABLE OF CONTENTS

1. Introduction	1
2. Name of parts	2
3. Specification	4
4. Operation instructions in detail	4
5. Safety instructions	6

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### 1. INTRODUCTION

The Glove device is an innovative solution for at-home skincare. Using advanced technology, Glove enhances hydration, improves skin firmness, and helps restore the skin's natural, healthy appearance.

**Device Features:**

- **Cleaning** – opens pores and removes impurities using a warm compress.
  - **Deep Hydration** – delivers active ingredients into the deeper layers of the skin.
  - **Cryotherapy** – reduces puffiness and fatigue around the eyes.
  - **EMS Technology** – stimulates facial muscles, improving skin tightness and elasticity.
  - **Collagen Synthesis** – supports natural collagen production, smooths wrinkles, and provides a lifting effect.
  - **Photon Light Therapy** – reduces fine wrinkles and discoloration, improves circulation, oxygenates the skin, and has an antibacterial effect.
- The device is intended for home use only and facial skincare. To ensure full safety, follow the guidelines provided in the user manual.

#### 1. Introduction

Did you know that our face contains over three million pores? Unfortunately, daily cleansing and traditional skincare products only clean the surface of the skin, leaving dirt, makeup residue, and excess sebum trapped deep inside the pores. These impurities block the skin and rob it of its healthy glow. Thanks to 8000 vibrations per minute, sonic micro-vibrations effectively remove deeply embedded impurities. The skin becomes clean, smooth, and radiant, while the pores are thoroughly cleansed. As a result, the skin absorbs active ingredients from cosmetics more effectively, making your skincare routine even more efficient.

#### 2. Hot & Cold Function

Heat therapy (42°C) opens the pores, allowing for thorough cleansing and better absorption of nutrients. Cooling therapy (6°C) instantly closes the pores, smooths the skin, and locks in moisture.

The combination of heat and cold improves blood circulation, reduces dark circles and puffiness, leaving the skin firm, smooth, and radiant – just like after a professional beauty treatment.

EN

1

### 3. EMS Stimulation

EMS stimulates skin cells, boosting collagen production, which enhances skin elasticity, smooths wrinkles, and leaves the skin soft and healthy.

#### Cooling Effect

1. Supports skin firmness, maintains elasticity, and prevents allergies and other skin issues.
2. Tightens pores, locks in moisture, and smooths out wrinkles.
3. Relieves muscle tension, reduces skin fatigue, redness, and pain, while preventing the formation of acne.

#### Heating Effect

1. Warms the skin, relaxes muscles, and alleviates pain.
2. Improves blood circulation and opens pores, allowing the skin to absorb active ingredients from cosmetics more effectively, leaving it nourished and hydrated.
3. Supports deep cleansing by removing toxins, blackheads, and impurities from the pores.

#### 4. LED Light Therapy


1. **Red Light:** Wavelength of 620 nm (± 10nm) – 18 LED lights penetrate the skin (830 nm), promoting blood circulation and stimulating collagen regeneration, improving skin elasticity.
2. **Blue Light:** Wavelength of 470 nm (± 10nm) – 18 LED lights penetrate the skin (480 nm) and have anti-inflammatory and antibacterial effects. The light reaches a depth of 2.5 mm, eliminating bacteria that cause acne. It helps reduce inflammation, acne scars, and improves blood circulation and skin condition.
3. **Yellow Light:** Wavelength of 600 nm (± 10nm) – 18 LED lights enhance oxygen supply to cells, stimulate lymphatic drainage, and improve skin texture. It is particularly effective in reducing redness and smoothing rough skin.

### 2. NAME OF PARTS

1. Hot and cool panel
2. LCD display
3. Charging jack
4. USB charging cable
5. Power & Mode switch
6. LED photon button
7. Level of EMS and Temperature
8. Circle of makeup cotton
9. Base

EN

2

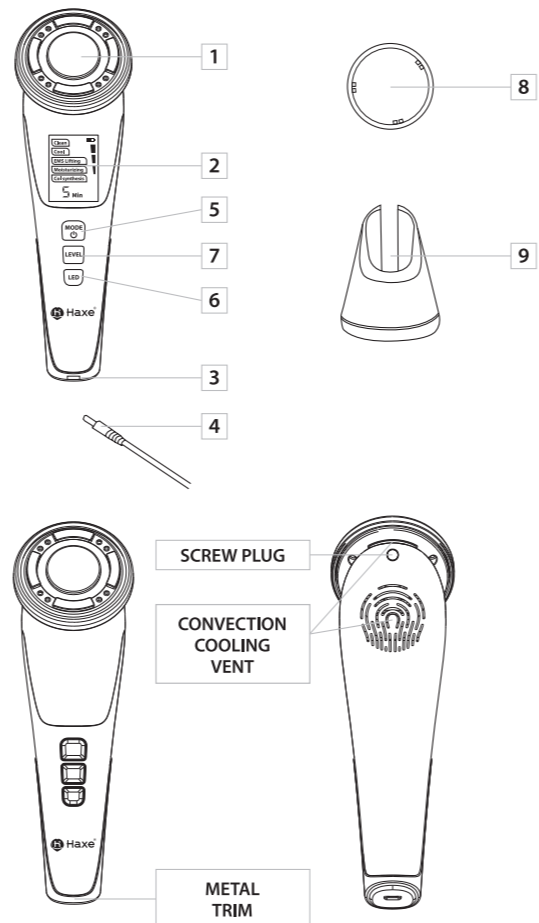EN

3

### 3. SPECIFICATION

1. Name: Needleless mesotherapy Glove
2. Input: DC 5V
3. Vibration: 8000 (10) rpm
4. Charging current: 500mA
5. Battery capacity: 1000mAh
6. Net weight: 145g
7. Machine size: 175\*65\*67mm
8. Model: HX808
9. Temperature: 6°C-42°C (5°C)
10. Storage temperature: -20°C-60°C
11. Charging time: 2 hours
12. Power: MAX 8.5W
13. Standby time: 45 days
14. Use time: 5min/10s or 2 min/10s
15. Box size: 215\*150\*80mm

### 4. OPERATION INSTRUCTIONS IN DETAIL

1. Please take full charge before the first using. The charging time is about 2.5 hours. Long press with 3 seconds on power key 5 to start machine, defaulting clean mode in 1 level, then the beauty instrument work with vibrate and heat. According to personal skin needs, temperature can be adjusted by press key 7. You can adjust LED Photon light by press key 6 in each mode (5min smart timing for mode 1).
2. Press the power key to the Cool mode, the beauty instrument work with vibration and cools. According to personal needs, you can press 7 to adjust the temperature (2min smart timing, after power off the fan still works for 30S to cool down)
3. Press the power key 5 to the EMS lifting mode. The beauty instrument vibrates and production high EMS micro current. According to individual needs, you can adjust the EMS level by press key7 (5min smart timing for mode 3).
4. Press key 5 to switch the moisturizing mode. The beauty instrument will vibrate & heat and sent EMS micro current to skin. According to individual requirements, you can adjust the EMS intensity by press key 7 (5min smart timing for mode 4).
5. Press the power key 5 to switch the Collagen Synthesis mode. The beauty instrument vibrates and cools while emitting EMS micro current to skin. (Same 2 min smart timing as model 2)

EN

4

According to individual requirements, you can adjust the EMS intensity by press key 7 (2min smart timing for mode 5).

### 6. Using Method as following pic and text

#### CLEAN MODE

Need to use with remove makeup cotton. Take off the Circle, attaching the cotton, pour with makeup remover water. Then turn on the clean mode, warm and open the pores, shatter the dirt and export the deep layer of rubbish, clear pores.

#### COOL MODE

Put the mask on the face and turn on the Cool mode. The cool function can effectively shrink the pores and lock the skin moisture.

#### EMS LIFTING MODE

After cleansing, apply skin care products on the face and turn on the EMS mode. EMS stimulates the muscles and makes the skin more energetic and elastic. It improves fine lines and remove crow's feet and wrinkles on face, and restores the skin's youthfulness, smoothness and tenderness.

#### MOISTURIZING MODE

After cleansing, apply skin care product such as lotion and essence to your face, and activate moisturizing mode to promote double absorption of skin care products! (On each mode you can select different light according to your skin needs)

#### COLLAGEN SYNTHESIS MODE

After the skin care products were absorbed, turn on the collagen synthesis mode, cool to shrink pores, lock the nutrition. The EMS micro-current stimulates the skin cells to synthesize collagen and rejuvenate the muscles anti-wrinkle. (On each mode you can select different light according to your skin needs)

EN

5

### 5. SAFETY INSTRUCTIONS

#### 1. Usage Safety

1. Read the entire user manual before using the device for the first time.
2. Use the device only as intended and as described in this manual.
3. Do not disassemble or modify the device. For repairs, contact the seller.
4. Do not use the device if you: have a pacemaker, metal implants, epilepsy, or skin conditions, are pregnant.
5. The device is not suitable for children or individuals who cannot use it safely. Adult supervision is required for users under 18.
6. Do not use the device near eyes, mouth, ears, breasts, thyroid gland, or genital areas.
7. Ensure the skin is moisturized with gel or serum during use to avoid irritation and ensure proper EMS function.
8. Stop using the device and consult a dermatologist if skin irritation or redness occurs.

#### 2. Electrical Safety

1. Do not touch the device with wet hands.
2. Use only the power adapter included in the package. Ensure the device is completely dry before plugging it in.
3. Do not use the device if the power cord or socket is damaged.
4. Do not store or use the device in damp places (e.g., bathroom).
5. Keep the device away from water, fire, and high-temperature areas.

#### 3. Storage and Maintenance

1. Store the device in a dry, cool place, away from heat sources and flammable materials.
2. Do not use aggressive chemicals to clean the device.
3. Turn off the device after use and before cleaning.
4. If the device becomes hot during use, allow it to cool for 10 minutes before reusing.

#### 4. Usage Guidelines

1. Usage time – do not exceed the recommended operating time specified in the manual.
2. The cooling function consumes more energy. Allow the device to cool down before using it again.
3. Always use appropriate moisturizing products (gel or serum) during treatment to protect the skin.
4. Check that all parts of the device are in good condition before use.

EN

6

### 5. Key Warnings (Summary)

1. Do not use the device if you have health conditions (e.g., pacemaker, pregnancy).
2. Keep the device away from water and heat sources.
3. Use a gel or serum for safe and effective treatments.
4. Turn off the device after use.
5. Stop use if skin irritation occurs.

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