Panasonic

Operating Instructions (Household use) Sonic Vibration Toothbrush Model No. EW-DA37

Series 600

English	3	Norsk	184
Deutsch	23	Suomi	204
Français	44	Polski	224
Italiano	64	Česky	244
Español	84	Slovensky	264
Nederlands	104	Magyar	284
Português	124	Română	305
Svenska	144	Türkçe	325
Dansk	164		

Operating Instructions

(Household use) Sonic Vibration Toothbrush

Model No. EW-DA37

Contents

Safety precautions.....	4
Important information	7
Parts identification	8
Before cleaning your teeth	9
Features.....	9
Charging.....	10
Basic usage.....	11
Tips on usage	12
After use	13
Changing the brush head (Sold separately).....	15
FAQ	16
Troubleshooting	17
Battery life	20
Removing the built-in rechargeable battery.....	20
Specifications	21

Thank you for purchasing this Panasonic product.
Before operating this unit, please read these instructions completely and save them for future use.

WARNING

- Be sure to switch off the appliance and remove the power plug from a household outlet before performing any maintenance. Failure to do so may result in electric shock.
- Do not immerse the USB power adaptor and USB cable in water or clean them with water.
Doing so may cause electric shock or fire due to a short circuit.
- Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and user maintenance shall not be made by children without supervision.
- Children shall not play with the appliance.
- Do not use anything other than the supplied USB power adaptor or USB cable for any purpose. Also, do not use any other product with the USB power adaptor or USB cable. (See page 8.)
- The following symbol indicates that a specific detachable power supply unit is required for connecting the electrical appliance to the supply mains.
The type reference of power supply unit is marked near the symbol.

Safety precautions

Make sure to follow these instructions.

To prevent accidents, injuries or damage to property, please follow the instructions below.

- **The following chart indicates the degree of damage caused by improper operation.**

DANGER

Indicates a potential hazard that will result in serious injury or death.

WARNING

Indicates a potential hazard that could result in serious injury or death.

CAUTION

Indicates a potential hazard that could result in minor injury or property damage.

WARNING

► Power supply

- Do not use in a way exceeding the rating of the household outlet or the wiring.
 - Exceeding the rating by connecting too many power plugs to one household outlet may cause fire due to overheating.
- Fully insert the power plug of the USB power adaptor securely.
 - Failure to do so may cause fire or electric shock.

► USB power adaptor and USB cable:

- Do not place them over or near water filled sink or bathtub.
- Never use the appliance if they are damaged or if the power plug fits loosely in a household outlet.
- Do not damage them.
 - Do not deface, modify, forcefully bend or pull them. Also, do not place anything heavy on, twist or pinch them.
 - Doing so may cause electric shock or fire due to a short circuit.
- Do not connect or disconnect them with a wet hand.
 - Doing so may cause electric shock or injury.
- Disconnect the USB plug when cleaning the main unit.
 - Failure to do so may cause electric shock or injury.
- Regularly clean them to prevent dust from accumulating.
 - Failure to do so may cause fire due to insulation failure caused by humidity.

Remove the USB power adaptor and wipe it with a dry cloth.

► In case of an abnormality or malfunction

- Immediately stop using the device and disconnect the USB power adaptor.
 - Failure to do so may cause fire, electric shock, or injury.

<Abnormality or malfunction cases>

- The main unit, USB power adaptor or USB cable is deformed or abnormally hot.
- The main unit, USB power adaptor or USB cable smells of burning.
- There is an abnormal sound during use or charging of the main unit.
 - Immediately request inspection or repair at an authorized service centre.

► This product

- This product has a built-in rechargeable battery. Do not throw into fire, apply heat, or charge, use, leave, or store in a high temperature environment.
 - Doing so may cause overheating, ignition, or explosion.
- Do not modify or repair.
 - Doing so may cause fire, electric shock, or injury.
 - Contact an authorized service centre for repair (battery change, etc.).
- Never disassemble except when disposing of the product.
 - Doing so may cause fire, electric shock, or injury.

WARNING

► Preventing accidents

- Do not store within the reach of children or infants. Do not let them use it.
 - Doing so may cause an accident or injury due to accidental ingestion of removable parts.

CAUTION

► Preventing damage to teeth and gums

- Use the toothbrush with a general type of toothpaste only. Do not use specially formulated products such as ones enhanced for tar removal.
- Do not apply the brush head strongly against teeth or gums, or use it continuously against same spot.
- Do not apply the plastic part of the main unit or brush head to teeth.
- Do not let anyone who cannot operate the brush by themselves use it.
- People with low sensitivity in their mouth should not use the brush.
 - Doing so may cause injury or damage to teeth and gums.
- People who may have gum disease or who are concerned about other symptoms in their mouth should consult a dentist before use.
 - Failure to do so may cause injury or damage to teeth and gums.

► This product

- Never use for a purpose other than cleaning your teeth.
 - Doing so may cause injury.
- Do not share the brush head with your family or other people.
 - Doing so may result in infection or inflammation.

► Note the following precautions

- Do not drop or hit against other objects.
 - Doing so may cause injury.
- Do not allow foreign objects like metallic pins or lint on the USB power adaptor's power plug or USB plug.
 - Doing so may cause electric shock or fire due to a short circuit.
- Do not charge when the appliance socket is wet.
 - Doing so may cause burns or fire.
- Disconnect the USB power adaptor from the household outlet when not charging.
 - Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.

Handling of the removed battery when disposing

DANGER

- The rechargeable battery is exclusively for use with this appliance. Do not use the battery with other products.
- Do not charge the battery after it has been removed from the product.
- Do not throw into fire or apply heat.
- Do not hit, disassemble, modify, or puncture with a nail.
- Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.
- Do not carry or store the battery together with metallic jewellery such as necklaces and hairpins.
- Do not use or leave the battery where it will be exposed to high temperatures, such as under direct sunlight or near other sources of heat.
- Never peel off the tube.
 - Doing so may cause overheating, ignition, or explosion.

WARNING

- After removing the rechargeable battery, do not keep it within reach of children and infants.
 - The battery will harm the body if accidentally swallowed.
If this happens, consult a physician immediately.
- If the battery fluid leaks out and comes in contact with the eyes, do not rub the eyes, and thoroughly rinse with clean water, such as tap water.
 - Failure to do so may cause eye injury.
Consult a physician immediately.


CAUTION

- If the battery fluid leaks out and comes in contact with the skin or clothes, rinse off with clean water, such as tap water.
 - Failure to do so may cause inflammation.

Important information

- Charge the appliance correctly according to these Operating Instructions. (See page 10 "Charging".)
- This appliance contains a battery that is only replaceable by skilled persons. Please contact an authorized service centre for the details of repair.

Parts identification

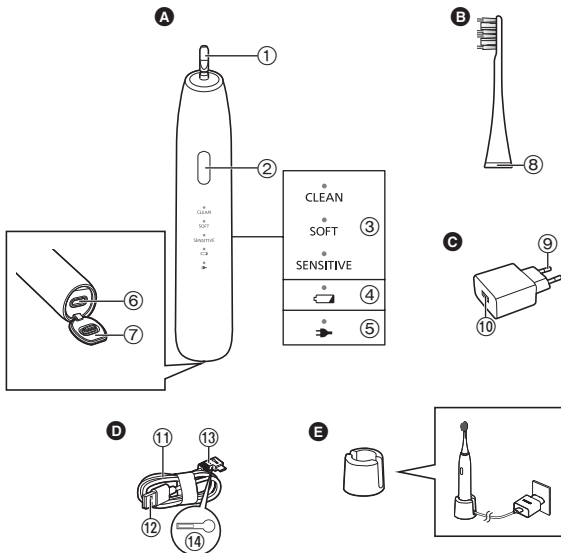A Main unit

- ① Toothbrush shaft
- ② Power switch (power off/on/mode selection)
- ③ Mode indicator (blue)
- ④ Battery capacity indicator () (blue)
- ⑤ Charge status indicator () (red)
- ⑥ Appliance socket (Type-C)
- ⑦ Socket cover

B Total care brush

- ⑧ Identification ring

C USB power adaptor (ES-RU1-02)

- ⑨ Power plug
- ⑩ USB port (Type-A)

D USB cable

- For charge only
- ⑪ Cord
- ⑫ USB plug (Type-A)
- ⑬ USB plug (Type-C)
- With overheat protection function
- ⑭ Overheat protection function mark

Accessory

E Main unit stand

Before cleaning your teeth

► This product

The Sonic Vibration Toothbrush is an electric toothbrush with sonic vibration* of approximately 31 000 brush strokes/minute.

* Sonic vibration in acoustic wave field

► Bleeding of the gums

When you first begin using the toothbrush, you may experience slight gum bleeding even if your gums are healthy. This is because your gums are receiving stimulation by the toothbrush for the first time. The bleeding should stop after one or two weeks of use. Bleeding that persists for more than two weeks may be indication of an abnormal gum condition, such as gum disease. In this case, stop using the toothbrush and consult your dentist.

Features

► About overpressure protection sensor


To avoid damaging the teeth and the gums by brushing too hard, the toothbrush will temporarily stop its brushing action when the sensor detects excessive brush contact pressure, and alerts the user with a series of rapid beeps.

- This feature will not be available in the SENSITIVE mode.
- You can choose to turn the sensor function on or off. (it will be on by default)

How to turn the sensor on or off

Press and hold the power switch for 5 seconds while the toothbrush is being charged.

- “Setting Change Successful” tones
When the sensor is turned from on → off: One long beep
When the sensor is turned from off → on: One short beep

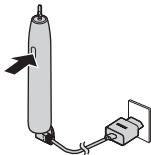► About soft start function

To get accustomed to the high speed vibration, the toothbrush will gradually reach the selected speed and vibration approximately 2 seconds after it starts operating.

Charging

1 Open the socket cover and insert the USB plug firmly into the main unit and USB power adaptor.




- Wipe any water or toothpaste residue off the appliance socket area.

2 Insert the USB power adaptor firmly into a household outlet.

- Charging is completed after approx. 2 hours.

3 Disconnect the USB power adaptor from the household outlet when the charge status indicator () turns off.

- Close the socket cover.

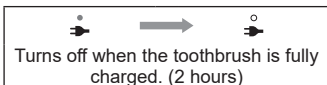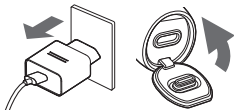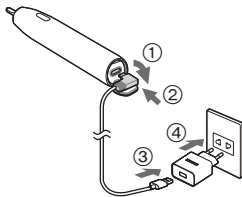► When to charge

The toothbrush will give indicator and tone signals to alert the user when the battery charge level gets low.

- When the toothbrush is switched on: the battery capacity indicator () will blink (once per second)
- When the toothbrush is switched off: A series of beeps will be given (Two beeps) (the indicator turns off after 5 seconds)

► Even when the toothbrush battery has run out of charge

- The toothbrush will be good for use for about 2 minutes after charging for 4 minutes. (Quick Charge feature)
(The available use time may vary depending on the ambient temperature and the manner in which the toothbrush is used.)

► Battery runtime on a full charge

Mode	Battery runtime
CLEAN SOFT SENSITIVE	Approx. 22 days*

(The available use time may vary depending on the ambient temperature and the manner in which the toothbrush is used.)

* Using for 2 min. 2x per day



Notes

- The charging time may vary when using for the first time, or when it has not been in use for more than 6 months. Also, the charge status indicator () may not glow for 5 to 10 minutes after charging starts, but it will glow if charging continues.
- The required charging time may vary depending on the starting charge level.
- Recommended ambient temperature for charging is 10 °C – 35 °C. When the ambient temperature is out of the recommended range, the battery may not charge properly or may take longer to be fully charged.
- If the charge status indicator () blinks quickly (twice per second), see page 17.

Basic usage

1 Attach the brush head to the main unit.

There is a small gap between the brush head and the main unit (a few mm), but this is not a malfunction.

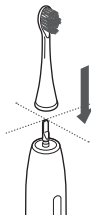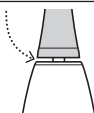2 Apply toothpaste.

- Put a small amount of toothpaste on the brush head.


3 Put the brush head in your mouth and press the power switch to turn on the main unit.

► Select the desired mode

To switch the mode, press the power switch within 5 seconds.


- The last mode used for more than 5 seconds will be stored.

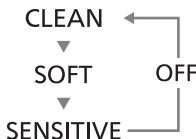4 Lightly press and move the toothbrush against your teeth.

► Example brushing sequence

- ① front of the teeth
- ② biting plane of the teeth
- ③ back of the teeth

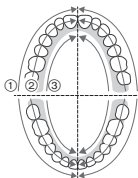► About 30-second brush timer

- The toothbrush will pause once at 30-second intervals to indicate the brushing time. It pauses twice at 2-minute intervals then starts again.

5 Turn off the main unit before removing it from your mouth.

Notes

- The appliance may not operate in an ambient temperature approximately 5 °C or lower.
- Use the appliance in the appropriate mode.

Mode	Strength	Usage or desired care
CLEAN		For everyday cleaning at standard intensity
SOFT		For first-time users or for gentle brushing
SENSITIVE		For extra gentle brushing to prevent gum irritation

- The toothbrush will temporarily stop its brush action when it senses excessive brush contact pressure against the teeth and the gums, and alerts the user with a series of rapid beeps. (See page 9.)

Tips on usage


Total care brush

This one brush head serves all your needs for gum care, plaque removal and stain off.

How to apply/move the toothbrush

Front surface and the biting surface of teeth

- Apply the brush head at a 90° angle.

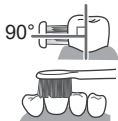Gum pocket and the boundary between the teeth and gums

- Place the brush head at a 45° angle.
- Apply the bristles lightly into the gum pocket and move the toothbrush slowly.

For the back of the front teeth or misaligned teeth

- Place the brush vertically.
- Move toothbrush to match the angle of the teeth.

After use

- Do not use thinner, benzine, alcohol, etc. to clean the appliance. Doing so may cause failure, cracking or discoloration of the parts.


Remove the brush head.

Brush head

Wash the brush head under running water.

- Do not wash with water hotter than 80 °C.
Doing so may damage the brush head.

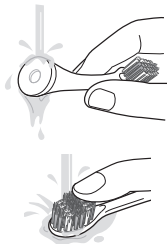Main unit


Close the cover and wash with water

Before washing the toothbrush, push the cover down to see it is tightly closed.

Wash the main unit under running water.

- Do not damage or rub the waterproof rubber seal (A) with hard materials.
- Do not leave the main unit under water.
- Do not wash with water hotter than 40 °C.
Doing so may damage the main unit.

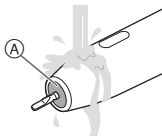Wipe off water from the brush head and main unit for storage. (This is to keep the toothbrush clean and to reduce malfunction risks.)

- Do not charge the toothbrush when it is wet.

Main unit stand

Wash and wipe off water.


Replacing the identification rings

► Removing

Slide up with your thumb.

► Attaching

Insert by covering the indentation (a) in the brush head with the inside of the identification ring (b).

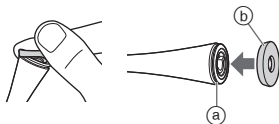If the socket cover comes off

Push the cover onto the appliance socket straight down.

Changing the brush head (Sold separately)

- Brushes are consumables.
- The brush should be replaced with a new one every 3 months.

Part number	Replacement Brush for EW-DA37	
WEW0810	Total care brush	

	Question	Answer
While using	My toothbrush temporarily stopped with rapid beeps. Why? ▶	The overpressure protection sensor has been activated. Open your mouth lightly and apply the brush to your teeth with gentle pressure. If the beeps bother you, you can choose to turn the sensor off. (See page 9.)
	Why does my toothbrush temporarily stop or reduce the brush action about every 30 seconds? ▶	The toothbrush will notify you of the elapsed brushing time every 30 seconds. (See page 12 for description of the 30-second brush timer.)
Battery	Will the battery deteriorate if not used for an extended period of time? ▶	The battery will deteriorate if not used for over 6 months, so make sure to perform a full charge at least once every 6 months.
	Can I charge the appliance after every use? ▶	The battery used is a lithium-ion, so charging after every use will have no effect on battery life.

	Problem	Possible cause	Action
While charging	The charge status indicator () does not glow when charging starts.	▶ You have just purchased the toothbrush or have not used it for a while.	▶ The charge status indicator () may not glow for 5 to 10 minutes after starting to charge it, but it will glow if you continue to charge it. If the charge status indicator () does not glow after a while, consult with an authorized service centre for repair.
		▶ The USB cable is not correctly connected with the main unit and/or the USB power adaptor.	▶ Make all connections correctly and check that the charge status indicator () glows. (See page 10.)
	The charge status indicator () is blinking rapidly. (Twice every second)	▶ The USB power adaptor is not correctly connected to the household outlet.	
		▶ The temperature is too low or high while charging.	▶ Charge the toothbrush within the recommended charging temperature of 10 °C to 35 °C. If the charge status indicator () blinks even when charging within the recommended temperature range, please contact an authorized service centre for repair.

Problem		Possible cause	Action
While charging	Toothbrush does not charge.	▶ The USB cable is not correctly connected with the main unit and/or the USB power adaptor. ▶ The USB power adaptor is not correctly connected to the household outlet.	Make all connections correctly and check that the charge status indicator () glows. (See page 10.)
	The toothbrush can only be used for a few minutes, even if it is fully charged.	▶ The wet appliance socket activated the overheat protection feature of the USB cable to suspend the charging action.	Allow the appliance socket to dry before recharging.
		▶ The rechargeable battery has reached its operating life. (Approx. 3 years)	Replace the battery. Contact the authorized service centre where purchased.

	Problem	Possible cause	Action
While using	The vibration increases gradually when pressing the power switch.	▶ The “Soft start function” is in operation.	▶ To avoid exposing the user to rapid mechanical vibration all at once, the toothbrush takes 2 seconds after power-on to increase its vibration action to the brushing level. Afterward, the toothbrush will be ready for use.
	I can feel the toothbrush shaft move when squeezed by hand.	▶▶▶▶▶▶▶▶▶▶	This is required for operation, and is not a malfunction.
	I cannot turn the overpressure protection sensor on or off.	▶ The battery charge level is low.	▶ Some setting changes may be unavailable when the battery charge level is low. Charge the toothbrush for a while and try again.
	The battery capacity indicator () turns off without blinking.	▶ The rechargeable battery has reached its operating life. (Approx. 3 years)	▶ Replace the battery. Contact the authorized service centre where purchased.
	Toothbrush does not operate even after charging.	▶ The appliance is used with approximately 5 °C or lower.	▶ Do not use the appliance with approximately 5 °C or lower.

If the problems still cannot be solved, contact the store where you purchased the unit or a service centre authorized by Panasonic for repair.

Battery life

Approximately 3 years. The battery may have reached its operating life if the usage frequency decreases drastically, even after a full charge.
(The operating life of the battery may differ depending on usage or storage conditions.)

Removing the built-in rechargeable battery

Remove the built-in rechargeable battery before disposing the toothbrush.

The battery is to be disposed of safely.


Please make sure that the battery is disposed of at an officially designated location if there is one.

This figure must only be used when disposing the toothbrush, and must not be used to repair it. If you dismantle the toothbrush yourself, it will no longer be waterproof, which may cause it to malfunction.

- Press the power switch to turn on the power and then keep the power on until the battery is completely discharged.
- Follow the procedure below, using a screwdriver to disassemble the appliance.

Perform steps ① to ③ and lift the battery, and then remove it.

- ② Insert the tip of a flat screwdriver into the groove (B) in the bottom cap, and pry the cap open with the screwdriver.
 - ③ While pushing at the toothbrush shaft (See page 8.), push the two battery holder tabs (C) inward with the tip of the screwdriver to disengage the tabs. Afterward, remove the battery.
- Take care not to short circuit the positive and negative terminals of the removed battery, and insulate the terminals by applying tape to them.

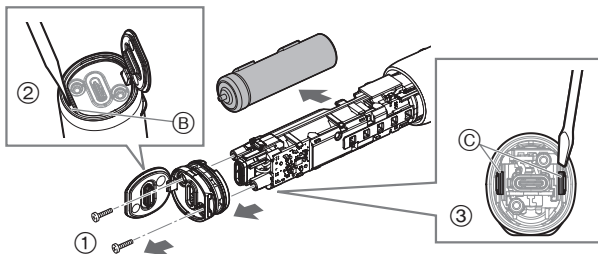For environmental protection and recycling of materials

This toothbrush contains a Li-ion battery.

Please make sure that the battery is disposed of at an officially assigned location, if there is one in your country.

Specifications

Power supply	See the name plate on the product. (Automatic voltage conversion)
Charging time	Approx. 2 hours (Room temperature: 10 °C – 35 °C)
Airborne Acoustical Noise	60.0 (dB (A) re 1 pW)

This product is intended for household use only.

For more information about the energy efficiency of the product, please visit our website, www.panasonic.com, and enter the model number in the search box.

Disposal of Old Equipment and Batteries

Only for European Union and countries with recycling systems

These symbols on the products, packaging, and/or accompanying documents mean that used electrical and electronic products and batteries must not be mixed with general household waste.

For proper treatment, recovery and recycling of old products and batteries, please take them to applicable collection points in accordance with your national legislation.

By disposing of them correctly, you will help to save valuable resources and prevent any potential negative effects on human health and the environment.

For more information about collection and recycling, please contact your local authority.

Penalties may be applicable for incorrect disposal of this waste, in accordance with national legislation.

Note for the battery symbol (bottom symbol)

This symbol might be used in combination with a chemical symbol. In this case it complies with the requirement set by the Directive for the chemical involved.

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