# FIZJOLOVE

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# 1. INTRODUCTION


EN

Thank you for choosing the FIZJOLOVE massager. The device combines graphene heating technology, vibration massage, and pneumatic compression (air cushions) to provide relief and promote recovery for your knees.

## 2. PACKAGE CONTENTS

- Main massage module (knee pad)
- Magnetic control controller
- Extension strap
- USB-C charging cable
- Herbal pouch (maxa)
- User manual

## 3. CONTROL PANEL

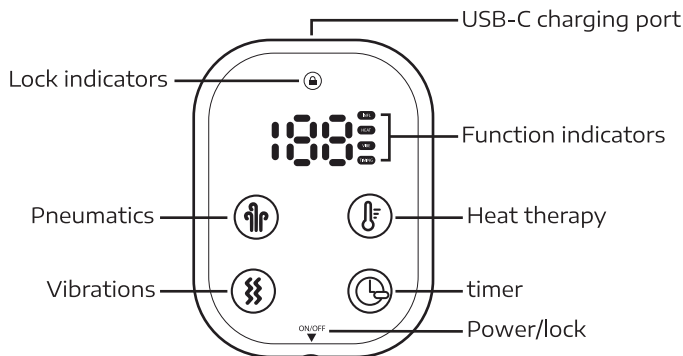- **Pneumatics:** Control of air cushions (3 intensity levels).
- **Vibrations:** Mechanical massage (3 operating modes).
- **Heat therapy:** Temperature adjustment with 5 levels (45°C to 65°C).
- **Timer:** Working time setting (default 30 minutes).
- **Power/Lock:** Turn on and lock the screen to prevent accidental touches.

## 4. INSTRUCTIONS FOR USE

### 1. Attachment:

Place the massager on your knee and adjust it using the Velcro straps. If the circumference is too small, use the extension strap.

### 2. Power On:

Press and hold the power button for 2 seconds. All icons will briefly light up, and the device will start the default 30-minute operating mode.

### 3. Pneumatic Function:

Levels: 1 (kneading), 2 (squeezing), 3 (full massage), and off mode (0).

### 4. Heat Therapy:

Range: 45°C, 50°C, 55°C, 60°C, 65°C, and off mode (0). Each press of the button cycles through the temperature levels.

Note: If using the herbal pouch (Moxa), do not exceed 60°C.

### 5. Vibrations:

The vibration function has 3 modes: 1, 2, 3. Each press changes the mode.

### 6. Timer:

When the device is turned on, the default working time is 30 minutes. Each press of the button changes the time to 30, 60, or 90 minutes. The remaining time is displayed in the device's lock mode.

### 7. Screen Lock:

The device automatically locks after 20 seconds of inactivity (lock icon turns red). To unlock the panel, press the power button briefly (icon turns white).

### 8. Charging:

Connect the USB-C cable to the controller port. The display will show battery status (0–99%). The message “FU” indicates full charge. When the battery level drops below 20%, the indicator will flash red and you will hear two beeps every 30 seconds. To check the charge level, briefly press the power button while the device is off. The device does not operate while charging. Use a standard 5V=2A power adapter for charging.

### 9. Error Alarm:

If the controller is not connected to the main module, the device will emit 6 beeps and then turn off after 5 seconds.

## 5. SAFETY PRECAUTIONS

1. The product is intended for daily care use only and is not a medical device.
2. Use the device only for its intended purpose and when dry.
3. Keep the product out of reach of children.
4. Do not use the device if it has been exposed to water or any other liquids.
5. Avoid contact of the device's surface with chemicals (e.g., perfumes, deodorants).
6. Protect the product from mechanical damage and sharp objects.
7. Do not use the device in humid, flammable, or explosive environments.
8. Use only a standard charger with appropriate specifications (5V=2A). Using fast chargers may cause the device to overheat.
9. Stop using immediately if you experience any discomfort.
10. People with sensitive skin should use the product carefully or over a layer of clothing.
11. Avoid overheating the device during use.
12. Do not leave the product unattended while in operation.
13. Do not use near flammable materials.
14. Do not use the device if any visible damage is present.

## 5. MAINTENANCE AND CARE

1. Do not place heavy objects on the device to prevent deformation or damage.
2. Clean with a soft, slightly damp cloth.
3. Do not use chemicals, solvents, or detergents.
4. For heavier dirt, clean gently without excessive rubbing.
5. Do not immerse the device in water or any other liquids.
6. Before using again, ensure the product is completely dry.

The device should be stored in a dry place, away from moisture, high temperatures, and direct sunlight. Keep it out of reach of children and protect it from contact with water or other liquids. If the device will not be used for an extended period, it is recommended to turn it off and periodically recharge the battery to maintain its performance.

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