FIZJOLOVE

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FIZJOLOVE

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
Thank you for purchasing a FIZJOLOVE product. Before using it, please read the user manual carefully to ensure proper use of the product. While reading, please pay special attention to the safety precautions and keep the manual in a safe place for future reference.

Content

1. Product description	1
2. Operation instructions	2
3. Usagemethod	3
4. Attentions	4
5. Maintenance and storage	5
6. Product parameters	5
7. Troubleshooting	6
8. Palm acupuncturepoints	7
9. Warranty card	8

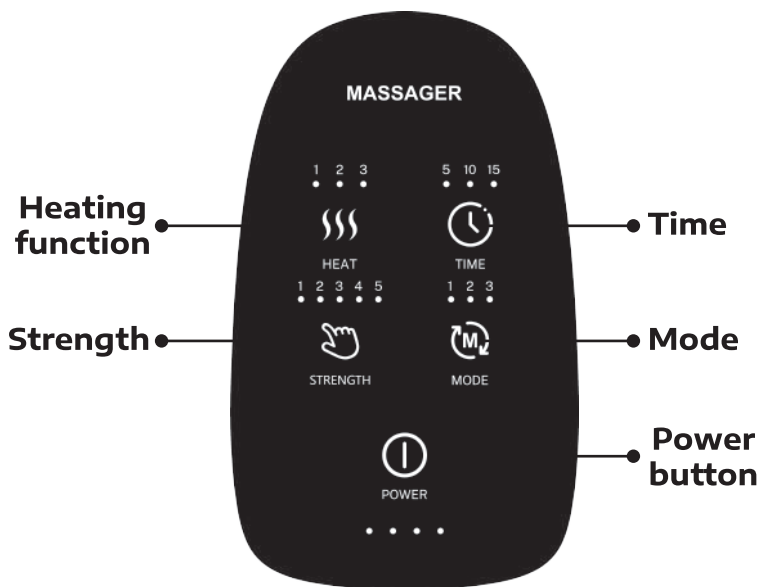1. Insert the back of your hand into the device.
2. Place your palm on the circular raised area in the center of the massage pad.
3. Press the **power** button on the panel – all button indicators will light up.
4. Press the **strength** button on the panel to select different intensity levels according to your hand size.
5. Press the **mode (M)** button to choose the massage mode: gentle, soothing, or strong – from low to high. There are 3 different intensity levels and massage modes in total.
6. To turn on the heating function, press **heat** to activate the low temperature range, then press again for medium, then high. Press once more to turn off the heating.
7. Press the **time** button to select the massage duration: 5, 10, or 15 minutes.
8. To stop the massage or turn off the device, press the **power** button.

Tips:

- ★ Try not to use it while charging.
- ★ Apply hand cream fully to your hands before using this device, which can help improve hand lines and skin.
- ★ Do not use it in heating mode for a long time to prevent low-temperature burns.

1. Use Type-C plug connection. (Note the shape of the insertion port)
2. Before use, please make sure that you have no jewelry (watch, ring, bracelet ring, etc.) in your hands, so as to avoid sharp objects piercing the airbag or squeezing the airbag and causing jewelry damage.
3. It is forbidden to use for a long time. Normally, massage time should not more than 30 minutes.
4. Please do not use this product while driving.
5. Please do not use it after getting drunk or taking sleeping pills.
6. Please use our charger (optional) or DC5V1-4A charger.
7. This product is specially designed as household electronic health care products. It is not suitable for diagnosis, treatment and other professional or commercial purposes.
8. Please keep this product away from water, fire and corrosive environment.
9. Please do not put sharp objects into this product.
10. Please do not use when hands are wet.
11. If the power cord is damaged, stop using it immediately. Contact the product service center for maintenance. Do not disassemble the product at will.
12. If you feel unwell, please stop using it in time.
13. Please do not wrap products with towels, quilts, clothes, etc. when you use it.

Please follow the doctor's advice for the following patients;

The patients who have had hand surgery, swelling and inflammation, congestion, severe arthritis, trauma and other hand diseases;
People with heart, brain and lung diseases;
Blood diseases with bleeding tendency; osteoporosis and;
soft tissue diseases; stroke; infants and pregnant women;


1. Use a towel to remove dust and other foreign matter and place it in a dry and ventilated place.
2. Please do not place it in places with high temperature, humidity and direct sunlight.
3. When the storage period exceeds 6 months, the product should be unpacked, inspected, and powered on.

Product Name	Intelligent Hand Massager		
Product No	ZH-023	Rated power	7.5W
Rated voltage	5V  1000mA	Protection time	15min
Battery Capacity	2000mAh	Power supply name	Lithium battery
Internal material	Antibacterial and deodorizing fabric	Massage way	Air bag massage
Outer size	250x154x111mm	Product weight	750g

If there is a malfunction during the use of this product, please check it as follows. If the fault can not be eliminated or the product still cannot run, please turn off the power and contact the product service center for maintenance.

Basic Faults	Possible causes	Troubleshooting Measures
The product can't start working	Check the remaining power of the battery The power may be weak.	Charging
The product cannot be recharged	Is the power adapter connected to the data line correctly?	Re-plug until LCD power display flickers
	Is the data line and USB socket intact?	Check the socket If there is any damage, please contact the manufacturer
	Are the specifications of adapters and data lines in line with operational requirements?	View specifications and proofread the technical parameters provided in the instructions
The product can be turned on, but can't work	The remaining power may be weak.	Charging

Do not disassemble or attempt to repair by yourself.

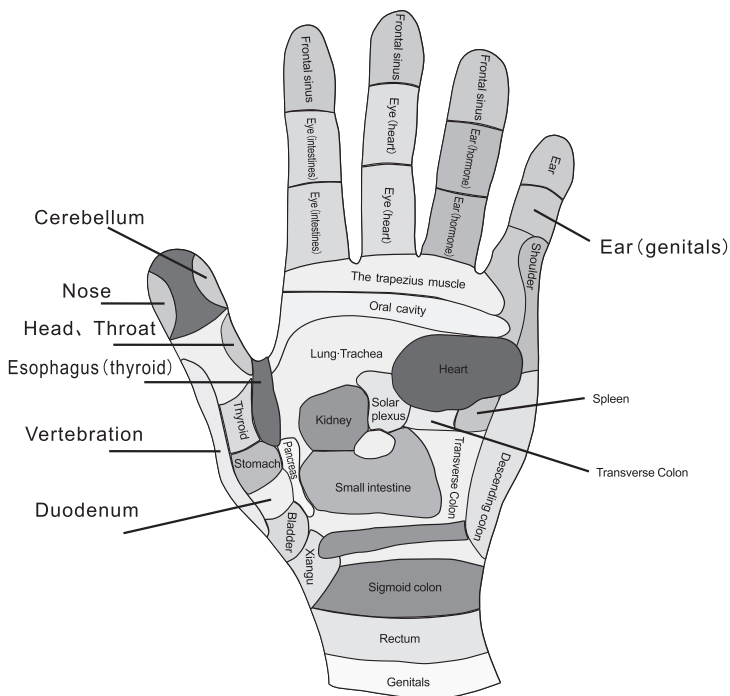Warranty card

Product name		Model	
User name		Contact number	
Address			
Post code		Purchase date	

Distributor

Address			
Post code		Contact number	

After-Sales Service Instructions

1. Users are entitled to a one-year free warranty from the date of purchase (the product's appearance is not covered by the warranty).
2. Please keep this warranty card and a valid proof of purchase for repair services. During the warranty period, both the warranty card and the proof of purchase must be presented at the same time.

3. The following conditions are not covered by the warranty:

- (1) Any damage caused by human factors, including use in an improper operating environment, incorrect operation, or improper storage not in accordance with the instructions.
- (2) The user has disassembled, repaired, or modified the device on their own, or the device has been repaired by a service provider not authorized by our company; absence of a warranty card or valid proof of purchase.
- (3) Damage caused by improper transportation or other force majeure factors (such as floods, fires, lightning strikes, or earthquakes) after the purchase of the device.
- (4) Our company reserves the right to amend and interpret the above terms.

www.fizjoloove.pl