# Oral-B®

## Professional Care

D 15525

D 15513

D 15511

powered by **BRAUN**

Type 4729

English	3, 18, 22
Français	4, 18, 22
Español	5, 19, 22
Português	6, 19, 22
Polski	7, 19, 22
Magyar	8, 20, 22
Český	9, 20, 22
Slovenský	11, 20, 22
Türkçe	12, 22
Русский	13, 21, 22
عربي	16, 15, 22


- (GB)** If you should require further assistance, please call our Consumer Service line.  
**0 800 73 11 792** (free call)
- (IRL)** **Helpline**  
**1 800 509 448** (free call)
- (F)** Appelez le Service Consommateurs Oral-B :  
**0 810 365 855** (numéro azur)
- (B)** **02-711 92 11**
- (LU)** Vous avez des questions sur ce produit? Appelez:  
**4 00 50 51**
- (E)** Servicio al consumidor para España:  
**901 11 61 84**
- (P)** Serviço ao Consumidor para Portugal:  
**808 20 00 33**
- (H)** **Információs vonal**  
**+36 (1) 801-3800**
- (CZ)** **Zelená linka**  
**00420-800-11-33-22**
- (SK)** **Zelená linka**  
**00420-800-11-33-22**
- (TR)** **Danışma Hattı**  
**0212-473-75-85**
- (RUS)** **Телефонный помощь**  
**8-800-200-1444**


**Internet:**  
[www.braun.com](http://www.braun.com)  
[www.oralb.com](http://www.oralb.com)

3-729-167/00/VI-05/M  
GB/F/E/P/PL/H/CZ/SK/TR/RUS/Arab  
Printed in Germany

D 15511

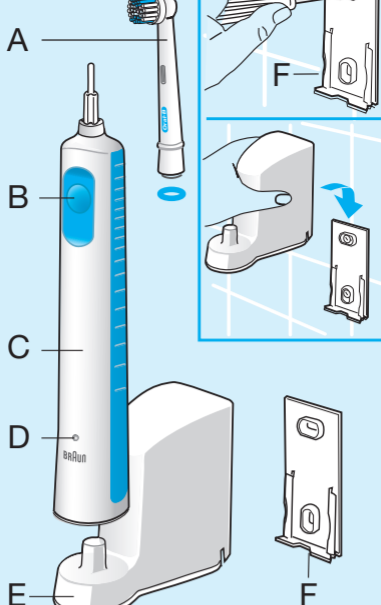D 15525

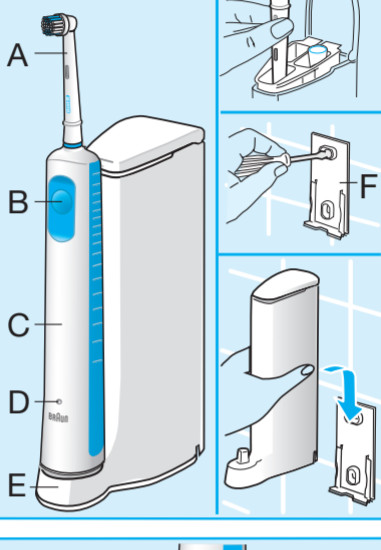
D 15513

I


a

b

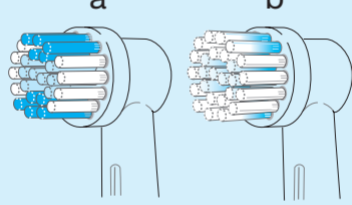II

III

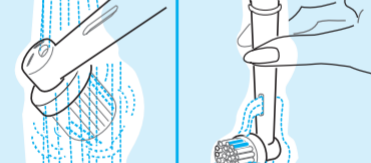IV

outsides

1

2

3


insides

4

tops

5
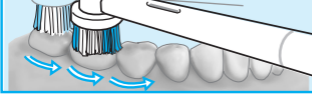
## English

Please read the use instructions carefully before first use.

### Important

Periodically check the cord for damage. Should this be the case, take the charging unit (E) to a Braun Service Centre.

A damaged or non-functioning unit should no longer be used.

## Description

- A Brushhead
- B On/off switch
- C Handle
- D Charge indicator
- E Charging unit
- F Wall holder (D 15511 and D 15525 only)

## Specifications

Power supply: see bottom of charging unit

Handle voltage: 1.2 V

## Connecting and charging

The toothbrush is electrically safe and can be used in the bathroom without hesitation.

- Plug the charging unit into an electrical outlet. Place the handle on the charging unit.
- A full charge takes approx. 12 hours and will provide an operation time of approx. 30 minutes.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the toothbrush at full power. Overcharging is impossible.

### Battery maintenance

To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and discharge the handle by regular use at least every 6 months.

## Using the toothbrush

The toothbrush can be used with any leading brand of toothpaste. However, we do not recommend toothpastes with whitening agents.

- To avoid splashing, guide the brushhead to your teeth before switching the appliance on.
- Slowly guide the brushhead from tooth to tooth, spending a few seconds on each tooth surface.
- Brush the gums as well as the teeth, first the outsides, then the insides and finally the chewing surfaces.
- Do not press too hard or scrub, simply let the brush do all the work.
- To thoroughly remove plaque, brush for at least two minutes.

### Memory timer

An incorporated timer memorises the elapsed brushing time, even when the handle is briefly (max. 30 sec) switched off.

After two minutes, the timer signals by a short stutter.

### Pressure Sensor

If you apply too much pressure, the oscillating movement will continue, but the pulsation will stop. There is a noticeable change in feeling and sound when the Pressure Sensor is activated.

During the first days of using the appliance, your gums may bleed slightly. In general, bleeding should stop after a few days.

Should it persist after two weeks, please consult your dentist or hygienist.

Dentists recommend that you replace any toothbrush every 3 months.

The INDICATOR<sup>®</sup> bristles help you monitor and improve your brushing performance.

With proper brushing, twice a day for two minutes, the colour on the INDICATOR<sup>®</sup>

bristles will disappear half-way down within 3 months, signaling it is time to change the brushhead.

### Replacement brushheads

Available at your retailer or Braun Service Centres.

## Cleaning

After use, rinse the brushhead for several seconds under running water with the handle switched on (II). Then switch off the handle and take off the brushhead. Rinse both parts separately under running water (III) and wipe them dry.

From time to time, also clean the charging unit with a damp cloth (IV).

## Environmental notice

This product is provided with a nickel-hydrate battery which does not contain any harmful heavy metals. Nevertheless, in the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at a Braun Service Centre.

However, if you prefer to dispose of the battery yourself, please refer to the chapter «Battery removal» on page 17.

Subject to change without notice.

This product conforms to the European Directives EMC 89/336/EEC and Low Voltage 73/23/EEC.

At the end of the product's useful life, please dispose of it at appropriate collection points provided in your country.

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## English

### Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose. This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Braun Customer Service Centre.

### For UK only:

This guarantee in no way affects your rights under statutory law.

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**Hong Kong**

Audio Supplies Company,  
Room 506,  
St. George Building,  
2 ICE House Street,  
Hong Kong,  
☎ (852) 2524 9377

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Tehran Bouran Company  
145 Ghaem Magham Farahani  
Tehran  
☎ +982 18 31 27 66

**Ireland (Republic of)**

Gavin's Electronics,  
83-84, Lower Camden Street,  
Dublin 2,  
☎ 1800 509 448

**Israel**

S. Schestowitz Ltd.,  
8 Shacham Str.,  
Tel-Aviv, 49517,  
☎ 1 800 335 959

**Italia**

Servizio Consumatori Braun  
Gillette Group Italy S.p.A.,  
Via G.B. Pirelli, 18,  
20124 Milano,  
☎ 02 / 6678623

**Jordan**

Interbrands  
Wadi Al-Sir, Adwiat alhekma  
St. Mekadabi bld.  
Amman  
☎ +692 582 75 67

**Kenya**

Radbone-Clark Kenya Ltd.,  
P.O. Box 40833,  
Mombasa Road,  
Nairobi,  
☎ 2 82 12 76

**Korea**

Gillette Korea Ltd.  
144-27 Samsung-dong,  
Kangnam-ku,  
Seoul, Korea,  
☎ 080-920-6000

**Kuwait**

Union Trading Company,  
Braun Service Center,  
P.O. Box 28 Safat,  
Safat Code 13001, Kuwait,  
☎ 04 83 32 74

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**Philippines**

Gillette Philippines Inc.,  
Corporate Corner Commerce Avenues  
20/F Tower 1, IL Corporate Centre  
1770 Muntinlupa city  
☎ 02-771071 02-06/-16

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**South Africa (Republic of)**

Fixnet After Sales Service,  
159 Queen Street, Kensington South,  
P.O. Box 751770,  
Johannesburg 2094,  
Kensington South,  
☎ 11 615 6765

**St. Maarten**

Rupchand Sons n.v. (ram's),  
P.O. Box 79,  
Philipsburg,  
Netherlands Antilles,  
☎ 52 29 31

**St. Thomas**

Boolchand's Ltd.,  
31 Main Street,  
P.O. Box 5667,  
00803 St. Thomas,  
US Virgin Islands,  
☎ 340 776 0302

**Suomi**

Gillette Group Finland Oy,  
P.O. Box 9,  
Niittykatu 8, PL 9, 02200 Espoo,  
☎ 09-45 28 71

**Sverige**

Gillette Group Sverige AB,  
Dept. Sweden,  
Stockholm Gillette  
Räsundavägen 12,  
Box 702,  
16927 Solna,  
☎ 020-21 33 21

**Syria**

Ahmed Hadaya Company  
Hadaya building  
Ain Keresh  
Unisyria, P.O. Box 35002,  
Damascus,  
☎ 963 011-231433

**Taiwan**

Audio & Electr. Supplies Ltd.,  
Brothers Bldg., 10th Floor,  
85 Chung Shan N Rd., Sec. 1,  
Taipei (104),  
☎ (886) 02 2523 3283

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**Yemen (Republic of)**

Saba Stores for Trading,  
26th September Street,  
P.O. Box 5278,  
Taiz,  
☎ 4-25 23 88

**Yugoslavia**

BG Elektronik,  
Bulevar kralja Aleksandra 34,  
11000 Beograd,  
☎ 11 3240 030

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Driva S.A.,
Gillette de Venezuela S.A.,