- PL** Instrukcja obsługi
- EN** User manual
- DE** Benutzerhandbuch
- FR** Manuel de l'Utilisateur
- IT** Manuale d'uso
- ES** Manual de usuario
- PT** Manual do usuário
- CZ** Uživatelský manuál
- SK** Používateľská príručka
- HU** Használati utasítás
- RO** Manual de utilizare
- NL** Handleiding

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INSTRUCTION MANUAL


Dear Customer!

Thank you for purchasing this acupressure mat. We encourage you to browse our complete range. Please read this instruction manual carefully and retain it for future reference.

The HYPERMOTION Team

Compliant with REACH (EC) 1907/2006 – no SVHC above 0.1%,
Tested by SGS: PAHs (AfPS GS 2019:01 PAK), REACH SVHC

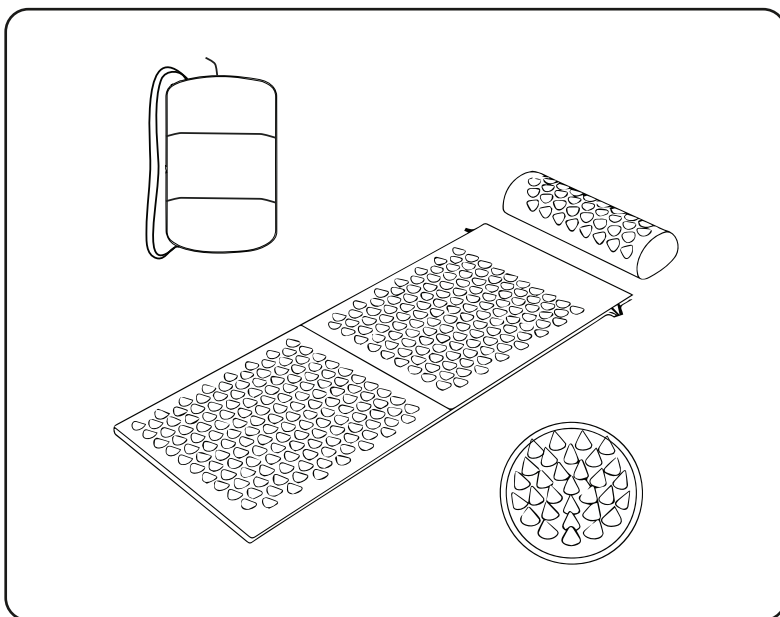THE IMAGE OF THE ACUPRESSURE MAT IS PROVIDED FOR GUIDANCE
PURPOSES ONLY. ACTUAL ITEM MAY DIFFER FROM THE IMAGE.

**WARNINGS!**

1. This product is intended only for adults and adolescents over 16 years of age. Do not use on children. This product is not a toy.
2. Do not use:
 - on damaged, irritated, or diseased skin
 - during pregnancy without medical advice
 - in the case of skin conditions, circulatory problems, epilepsy, or pacemakers
 - do not use without supervision on individuals with impaired sensory perception (e.g., neuropathy), impaired consciousness, or intellectual disabilities
 - do not use the mat on heated surfaces or in combination with a hot water bottle, heat lamp, or UV rays. The combination of heat and temperature may cause skin overheating or burns
 - in the case of active acne, psoriasis, or eczema
 - with advanced osteoporosis
3. Do not exceed the recommended use time (up to 30 minutes).
4. This product is not a medical device – do not use as a substitute for medical therapy. Intended for relaxation and wellness purposes only.
5. For home use only.
6. Do not use while sleeping (may cause skin irritation or damage).
7. Do not use while standing unsupported (risk of fainting or falling).
8. In rare cases, skin contact with ABS plastic may cause irritation. In the event of an allergic reaction, discontinue use immediately.
9. The manufacturer and importer are not liable for damage resulting from improper use of the product, contrary to the instructions.

Product Specifications:

- **Acupressure Mat:**
Cover: 100% cotton
Filling: 18D polyurethane (PU) foam
Dimensions: 130 × 50 × 2 cm
Spikes: ABS plastic - 414 round spiked elements (diameter 3.2 cm)
- **Cushion:**
Cover: 100% cotton
Filling: 18D polyurethane (PU) foam
Dimensions: 39 × 15 × 10 cm
Spikes: ABS plastic - 51 round spiked elements (diameter 3.2 cm)
- **Carrying Bag:**
Material: Black nylon

Benefits:

With nearly 8,000 pressure points stimulating the body's meridians, the set (mat + pillow) supports:

- Release of endorphins („happy hormones“)
- Improved blood circulation
- Reduction of pain and muscle tension
- Relaxation and stress reduction
- Increased energy
- Improved sleep quality

How to use?

- Lie on the mat directly on bare skin for best results.
- If you are a beginner, you can use a thin sheet or wear light clothing.
- Step on and off the mat very carefully, as the points are sharp and can cause intense pressure.

Material:

- Cover: 100% cotton
- Spikes: ABS plastic
- Filling: PU foam

Note:

- You may feel some discomfort at first – this is normal!
Breathe deeply and relax – this feeling will subside after a few minutes.
- If pain persists, discontinue use and consult a physician.

Application Areas:

- Neck and Back: Relieves Pain and Tension
- Legs and Feet: Activates Reflexology Points
- Abdomen and Chest: Gentle Massage and Relaxation
- Arms and Buttocks: Reduces Tension

Recommended Use Time:

- Beginners: 5–10 minutes per day
- Experienced Users: Up to 30 minutes per day

What Not to Do:

Do not machine wash – may damage the spikes or deform the foam.

Do not use bleach, fabric softeners, or abrasive cleaners.

Do not tumble dry – high temperatures may deform the ABS plastic.

Do not expose to prolonged sunlight or high temperatures.

Cleaning:

Remove the cover from the mat and/or cushion (if removable, check the zipper or Velcro). Gently hand wash in lukewarm water (up to 30°C) with a mild detergent. If you can't remove the cover, wipe with a damp cloth and a mild detergent. Do not rub the sharp spikes to avoid damaging them.

Drying:

Dry flat, in the shade, in a well-ventilated area. Avoid contact between the spikes and other fabrics when drying (they may damage them).

General care:

Store in a dry place, away from sources of heat and moisture. Avoid folding the mat in sharp places (e.g., in half) to avoid damaging the foam. Air the product regularly, especially if used directly on skin.

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