MERIDEN SONIC+ PROFESSIONAL

Sonic toothbrush with five brushing modes
and USB travel case

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TABLE OF CONTENTS:

Your new toothbrush	3
Contents	3
Additional accessories	3

Important tips for beginners

Safety tips	3
Important tips concerning batteries and their disposal.....	4
Product details.....	5
Turning the toothbrush on	6

Application

Cleaning and care.....	7
Removal of defects	7
Technical data	14
Warranty card.....	15

Information and answers to the most frequently asked questions concerning our products can be found on the website: www.meridenglobal.com

YOUR NEW SONIC TOOTHBRUSH

Dear Customer!

Thank you for purchasing the Meriden Sonic+ Professional sonic toothbrush with five brushing modes and an USB charger. In order to ensure the optimal use of the toothbrush please read the instruction manual and proceed in accordance with the following guidelines.

Contents:

(the configuration may vary depending on the model)

- sonic toothbrush
- charger
- USB cable
- 2 toothbrush heads
- USB case
- EU power adapter
- operating instructions, service instructions

IMPORTANT TIPS FOR BEGINNERS

Safety tips

- The purpose of this instruction manual is to familiarize the user with the various functions of the product. Please keep the instruction manual readily available.
- Improper use of the product may cause damage. CAUTION! Risk of personal injury!
- Do not open the device and do not repair it on your own.
- Please handle the product with appropriate care. The product can be damaged by shocks, impacts and falls from a small height.
- Do not store the device in damp and extremely hot premises.
- Do not immerse the device in water or other liquids.
- In the event of bleeding from gums lasting for an extended period of time, please consult your dentist.
- Before connecting the device, make sure that the voltage indicated on the device is compatible with the voltage in the electrical network.
- Do not leave the charger connected to a power outlet.
- Only use the toothbrush in accordance with the intended purpose described in this instruction manual. Do not use any accessories that are not recommended by the manufacturer.
- Do not attempt to open, repair or modify the product without authorization.
- If the toothbrush is dropped, the head of the toothbrush should be replaced before next use, even if there are not visible signs of damage.
- Do not use the charging base and the cable outdoors or near sources of heat.
- If the device does not work, do not attempt to disassemble it or to replace any part of the toothbrush as this could cause the risk of electric shock. Contact your technical support hotline.
- Keep the charger in a dry and safe place when it is not used
- In the event of a serious dental surgery, contact your dentist before using the toothbrush.
- Contact your dentist in the event of heavy gum bleeding after using the toothbrush, or if minor

bleeding of the gums persists for more than a week after you start using the toothbrush.

- The Meriden Sonic Professional toothbrush meets the safety standards for electromagnetic devices. If you have an implanted pacemaker or another medical device, please contact your physician or the manufacturer of the implanted device.
- Do not use the device if the tip of the toothbrush is damaged or if the bristles are already worn out. Damaged bristles may fall out during brushing. Replace the head of the toothbrush every 3 months or more frequently, if there are visible signs of wear.
- The device should not be used by persons with reduced physical, sensory or mental capabilities, and by persons who do not have any experience or knowledge in the use of this type of equipment.
- The Meriden Sonic Professional toothbrush is designed for cleaning your teeth, gums and the tongue. Do not use it for other purposes. The Meriden Sonic+ Professional toothbrush is a personal hygiene device, which is not intended to be used by a larger number of people.
- The product can be used by children under the supervision of adults, provided they understand the principles of operation of the device and the risks arising from improper use. The cleaning and maintenance of the device should not be performed by unattended children.
- Children must be supervised when using the toothbrush and must be warned against tampering with it. Children are not allowed to play with the toothbrush.
- We reserve the right to make technical changes and correct errors.

IMPORTANT TIPS CONCERNING BATTERIES AND THEIR DISPOSAL

Do not dispose of batteries with normal household waste. After use, discard the battery into an appropriate container. The battery can also be returned to a waste treatment point. The battery is permanently attached to the product and cannot be replaced. If the battery is stored for a longer period of time, we recommend that it is only charged after reaching 30% capacity. Avoid solar radiation at the place of charging - the ideal temperature is 10-20 degrees Celsius.

PRODUCT DETAILS

TURNING THE TOOTHBRUSH ON

- Place the charger on a flat and stable surface. Connect the charging cable with the USB charger. Place the toothbrush on the charging base. The charging process will start automatically and will be signaled by the LED indicator lighting up - after about 16 hours the LED indicator will turn off.
- Put the head of the toothbrush on the toothbrush handle.
- Moisten the head of the toothbrush with water and apply an appropriate amount of toothpaste.
- Turn the toothbrush on by pressing the function key.

SONIC TOOTHBRUSH MODES

The multifunction button is used to select the subsequent functions.

White: gently but effectively removes superficial surface discoloration and plaque.

Removes stains from tea, red wine, cigarettes and tobacco within a week of use.

Daily: for daily brushing


Polish: brightens and polishes your teeth; thanks to the high frequency mode, it increases the shine, and then brightens your teeth.

Massage: gently massages your gums.

Study: the entirely new Study mode ensures unique cleaning - from delicate to high frequency - by offering an adaptive process to new users.

The toothbrush should be moved slowly, using small movements and should not be pushed against the teeth too hard. A brief pause after 30 seconds indicates the transition to the next step of the brushing.

2 minutes times
30 seconds reminder

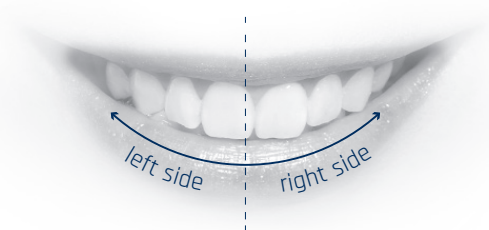CLEANING & CARE

- Turn the toothbrush off. After each use it is necessary to rinse the handle and the head of the toothbrush under running water, and to gently dry them with a cloth.
- The charger should be gently cleaned with a moistened cloth.
- The toothbrush should be kept in a dry place
- The head of the toothbrush should be replaced every three months or sooner, depending on the frequency of use.

THE WARRANTY DOES NOT COVER

- Damage resulting from random events (electrical fault, fire, inundation, flooding, etc.),
- Mechanical damage caused during the use of the device by the holder of the warranty, and any defects caused by them,
- Damage resulting from the user's fault,
- Damage resulting from the use of consumables not recommended by the manufacturer in the instruction manual,
- The breach of the warranty seals, or any other interference of persons other than authorized service personnel,
- Transportation and handling.

The warranty does not cover the failure of the product due to changes in the properties (reduction in quality) of components undergoing natural wear. Warranty repairs do not include the device maintenance and cleaning activities described in the instruction manual.

Disposal of unnecessary equipment by home users in the European Union.

The presence of this symbol on the product or its packaging indicates that it cannot be disposed of in the same way as household waste. As a result you are responsible for the disposal of waste equipment and are required to deliver it to an authorized facility recycling unwanted electrical and electronic equipment. The sorting, removal and recycling of waste equipment contributes to the conservation of natural resources and ensures that recycling is carried out in accordance with the principles of respect for human health and the environment. For more information about used equipment collection points, please contact your local authorities or your local household waste recycling point.

CE This product complies with the requirements of the EU Directive on electromagnetic compatibility and electrical safety.

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