
# AquaFloss Station

## EN - INSTRUCTIONS FOR USE

**Product description:**

- Storage for replacement heads
- Water tank
- Control wheel ON/OFF
- Rubber suction pads
- Lid
- 360° rotational socket
- Release button
- Water On/Off switch
- Base

### Instructions:

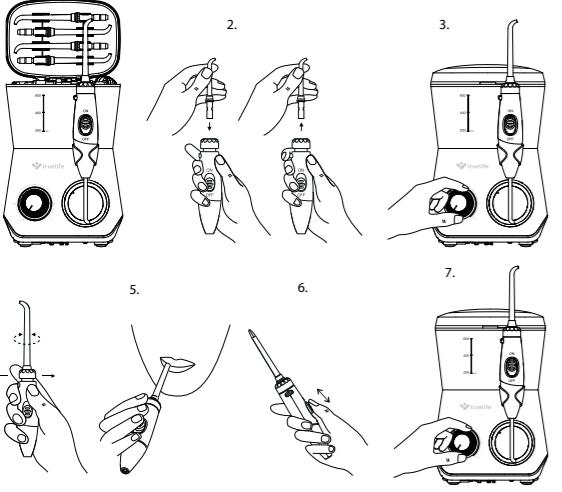**1. Preparing the water tank**  
Remove the water tank from the base and fill it with lukewarm water. Insert back into the base and connect it to the socket. Turn the On/Off control knob to activate the device.

**2. Inserting and removing the heads**  
Insert the head into the center of the dental nozzle handle. It should click into place. To remove once more, press and pull the head out.

**3. Nozzle pressure settings**  
Turn the On/Off pressure control knob to the desired setting or that recommended by your dentist.

**4. Usage**  
Lean over the wash basin and place the nozzle into your mouth. Direct it to your teeth. At the top of the handle, slide the Water On/Off switch up and commence cleaning.

**5. Recommended technique**  
Direct the nozzle at a 90° angle to the gums, gently squeeze your lips to prevent water spraying out and thereby allowing it to flow into the sink. For better results, start at the molar area (back teeth) and continue towards the front. Move the head over the gums, pause at the interdental spaces. Continue until all areas around and between the teeth have been cleaned.

**6. Water On/Off Switch**  
You can temporarily interrupt water flow from the nozzle by sliding down the switch.

**7. Termination**  
Switch off the device by turning the On/Off control wheel on the base, remove from the socket and pour out the remaining water.

**Package contents:**

- Appliance
- Manual
- Standard Jet

**Accessories:**

- Standard Jet – Rinses hard-to-reach places where a classic toothbrush does not reach.
- Orthodontic Jet\* – The special orthodontic tip is designed to clean and rinse hard-to-reach areas around dental braces.
- Periodontal Jet\* – The periodontal jet directs a stream of water or mouth wash below the gumline through a narrowed silicone tip.
- Dental Plaque Jet\* – Removes dental plaque effectively by providing easy access to affected areas.

\* not included

**Maintenance:**

- After the placement of specific liquids (e.g. mouthwash) in the water tank, the container should be rinsed out after each use.
- Fill the tank partially with clean, warm water, direct the nozzle into the sink, turn it on and allow it to run until the tank is completely empty. Regular flushing prevents clogging and prolongs product life.
- When necessary, clean the device with a fine cloth and non-aggressive cleaners.
- Before cleaning, please disconnect the appliance from its power source.
- Prior to potential exposure to frost, dry the appliance thoroughly and dry off any residual water.
- For hygienic reasons, replacement heads must be replaced at least once every 2-3 months.

**Potential problems and troubleshooting:**

Problem	Cause	Solution
Water leaking at the connection of the handle and the head	The head hasn't been properly connected	Remove the head and attach once again
Periodontal head is damaged	The head is worn-out	Replace the head
Water is leaking from the hose/handle	The hose/handle is damaged	Warranty
The tank is leaking	The tank's valve is missing, or has been fitted incorrectly	Position the valve correctly. Contact your dealer should the valve be missing
Low water pressure	The tank is incorrectly fitted	Ascertain that the outlet isn't blocked. Check the electric socket is functioning by testing it on other home appliances
Water is leaking from the appliance	The motor is switched off	Empty the tank and clear any water from the hose

The water smells	The hose, tank or brush head is blocked	Clean the tank, rinse the head with clean, warm water and change the hose
------------------	---	---

### Manufacturer's responsibility:

Truelife products are under warranty to cover any manufacturing defects. Eleme6 s.r.o. is obliged to repair or replace any part or product with a comparable component or product at its sole discretion. Conversion, alteration, or other unauthorized manipulation of the product exceeds the scope of this warranty. This applies even if defects, negligence or other errors in the product result in further loss or damage. The warranty does not apply to damage resulting from accidental or normal wear and tear. The company eleme6 s.r.o. assumes no responsibility for any accidents, injuries, death, loss or other claims related to or arising out of the use of this product. Eleme6 s.r.o. shall, under no event, be made liable for any incidental or consequential damage resulting from, or due to, the use of this product or any of its components. Repayment or replacement of spare parts and/or products may be subject to charges for transportation, handling, replacement and/or replenishment.

This product should be stored out of the reach of children. Children should be supervised by an adult.

### SAFETY WHEN USING THE PRODUCT:

- Read the safety instructions carefully.
- For optimal use of the product, please refer to the current instruction manual which is available for download at: [www.truelife.eu](http://www.truelife.eu).
- Printing errors reserved.

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