oromed
DBAJ O ZDROWIE

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DO ZĘBÓW

MODEL: ORO-SÓNIC BASIC

Spis treści

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Model produktu: ORO-SONIC BASIC

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Korpus szczoteczki: 182*19 (mm)

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Before operating this product, please read these instructions completely and save this manual for future use.

SAFETY INSTRUCTIONS

Please read the safety instructions before using this product. The safety instructions provided here are help you use this product safely and correctly, and also to prevent you from injuring yourself or others.


- Be sure to switch off the product and remove the plug from the electrical outlet before performing any maintenance. Do not touch the plug with wet hands.
- Never operate this product if it has a damaged cord, or if the cord is warm. Failure to do so may result in electrocution, short-circuit and/or fire.
- Do not damage, break or modify the cord. Do not bend, pull or twist the cord with unnecessary force. Do not place heavy objects on the cord or allow the cord to get pinched between objects.

- Do not modify, disassemble or repair the product at any time except when disposing of it.
- Do not use the charger in the bath. Do not wash the charger.
- Use this product only for its intended use as described in this manual. Keep the product out of reach of infants. Use the accessories recommended by the manufacturer.

FOR YOUR PERSONAL DENTAL CARE

Thank you for purchasing the ORO-SONIC BASIC which is designed for those who wish to keep their teeth clean and healthy. This product should be used personally only. People with tooth or gum disease or people receiving oral treatment should consult their dentist before use.

PARTS IDENTIFICATION

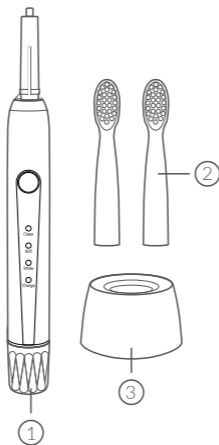1. Handle
2. Brush Head
3. Charging base

CHARGING

Place the toothbrush upright in the charger station. The charger indicator will light up red.

- The toothbrush charges whenever it is placed in the charger.
- The indicator light will change to green when fully charged.
- Charging time is 5-6 hours about.

Caution:

Do not leave metal objects, such as coins or paper clips, on the top of the charge station, they may become very warm.

BRUSHING MODES

The ORO-SONIC BASIC offers different brushing modes for varying oral care needs:

- Clean – Brush head swing in high frequency, remove tartar and bacteria on teeth.
- Soft – For sensitive teeth and gums, brushing your teeth in gentle frequency.

- White – Clean and Polish mode switching continuously and whitening teeth.

HOW TO USE

1. Place the brush head firmly onto the toothbrush shaft. There is a small space between the brush and the main unit. This is normal and is necessary for the head to vibrate properly.
2. Place the toothpaste on the head and wet it.
3. Place the toothbrush in your mouth. Turn the brush on.
4. Guide the brush head slowly from tooth to tooth, the ORO-SONIC BASIC will stop for a while each 30s for changing brushing area. After 2 minutes, ORO-SONIC BASIC stop automatically.

HOW TO SWITCH BETWEEN BRUSHING MODES

1. Power on/off and brushing mode button. Pressing this button within 5s, the operating sequence is power on, in clean, soft, white mode, and power off.
2. The ORO-SONIC BASIC automatically remember the last used mode, and starts in this mode when using again. To shift from one mode to another mode, press the brushing mode button within 5s.
3. When you press the brushing mode button, the brushing mode indicator will lighting, indicating the selected working mode.

CLEANING RECOMMENDATIONS

- After brushing, rinse your brush head thoroughly under running water.
- Remove the brush head from the handle and clean both parts separately under running water. Then wipe them dry before you put the handle onto the charging unit.

- Before cleaning your charging unit, make sure you unplug the adapter. Only use a dry or damp-dry cloth to wipe the surface of the charging unit.
- Place the charging unit onto a dry and safe place.
- Do not plug in the charging unit until it is completely dry.

DISPOSAL

Important: do not clean the charging unit under running water or immerse it in water. Place the charging station on dry and safe place.

This product contains a rechargeable battery. In the interest of protection the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposable can take place at appropriate collection points provided in your country.

HI-TECH MEDICAL Sp. z o.o., Sp. k.

ul. Ptasia 10

60-319 Poznań

Polska, Europa

www.oromed.pl

Wersja: 1_08.2019 r.