BEAUTIFLY

B-Shine

Blush

series

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B—Shine
Blush
series

3 LAMPS OF DIFFERENT DEPARTMENT INCLUDED

LONG LAMP LIFE OF 500 000 IMPULSES

5 LEVELS

WINDOW SIZE OF 3.9 CM²

CLEAR LCD DISPLAY

HANDY AND EASY TO USE

STORAGE COVER

THE POSSIBILITY OF REPLACING THE LAMP AFTER SANITATION

Multifunctional cosmetic laser B—Shine – Our IPL device uses E-IPL technology that works on hair removal, skin rejuvenation, acne treatment. Each of these functions is possible by applying the appropriate lamp. Multifunctional cosmetic laser B—Shine will not only be a stylish accessory, but also a small device with a great potential of several professional cosmetics devices. Choosing B—Shine means that you hold in your hands your private beauty salon.

Introduction

Thank you for choosing the multifunctional B-Shine cosmetic laser. Our IPL device uses E-IPL technology (filters harmful and invalid light to get effective light according to three different function) that works on hair removal, skin rejuvenation, acne treatment. Each of these functions is possible by applying the appropriate lamp. Multifunctional cosmetic laser B-Shine will not only be a stylish accessory, but also a small device with a great potential of several professional cosmetics devices. Choosing B-Shine means that you hold in your hands your private beauty salon. Congratulations and welcome to the Beautify users.

For the best results, we suggest that before the first use, read the instruction manual including all safety rules and warnings, for correct and safe using the device.

Function

B-Shine has three functions and separate lamps for each of them.

- Function 1: IPL hair removal
- Function 2: IPL skin rejuvenation
- Function 3: IPL remove acne

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Contraindications

1. The device is not intended for children under 15 years of age. Young people aged 15–18 can use the device with permission and / or under the supervision of parents or guardians. Adults over 18 years of age can use the device freely.
2. Do not use if the skin in the area to be treated is tanned or burned. Using the device on tanned or burned skin can cause temporary or even permanent discoloration of the skin.
3. Do not use if you are pregnant, try to get pregnant, you suspect you may be pregnant or you are breastfeeding. The B-shine device has not been tested for this occasion, therefore, do not use it in these cases. Hormonal changes occurring during pregnancy or breastfeeding may increase the risk of skin damage or skin sensitivity.
4. Do not use the device if you have a disease such as epilepsy or or light-induced seizures or migraine. Energy flashes can cause severe headaches and migraines.
5. Do not use on areas covered with tattoos, dark brown or black places (such as dark freckles, birthmarks, or warts). Tattoos or dark spots can absorb too a lot of light energy that can cause swelling of the skin, temporary or permanent blisters, discoloration or whitening of the skin.
6. Do not use the device if you have a tendency to form keloids.

Do not use on the scalp, ears, neck, nipples or on mucous membranes (lips

7. vulva smaller, around the vagina and anus). At those points of the bikini zone where the skin has stronger pigmentation, use a lower level of light intensity. Use in these places may cause swelling of the skin, blisters or temporary discoloration or skin bleaching.

8. Do not use if any active inflammation of the skin occurs in the treatment area, including acne, psoriasis, pimple, burns, herpes, abrasions or epidermis rash. These factors can make the skin more sensitive to light and make the procedure work unpleasant or cause swelling, blisters or infections on the skin. If the skin is enough healed, before using the device, consult a doctor.
9. Do not use the device if you are taking medicines that sensitize the skin to light including non-steroidal anti-inflammatory agents (e.g. aspirin, ibuprofen, cetaminophen), tetracyclines, phenothiazines, thiazides, diuretics, sulfonyleureas, sulfonamides, DTICs, fluorouracil, vinblastine, griseofulvin, alpha-hydroxy acids (AHAS), Beta-Hydroxi Acids (BHAS), Retin-A, Accutane and / or local retinoids.
10. Do not use if you have a disease related to photosensitivity, e.g. porphyria, polymorphic light eruption, solar urticaria, lupus etc. These conditions make it the skin is more sensitive to light and the treatment can be unpleasant or causing swelling of the skin, blisters or infections.

Operating warnings

1. The product will not turn on without inserting the lamp, insert the lamp before connecting to power. Do not block the intake air openings and the and ventilation holes during operation so that it does not adversely affect heat generation.

2. Wear protective goggles when using the device.

3. Clean the area of skin to be treated before use. It is not allowed to apply any cosmetics to avoid the destroy of instrument and affect result.

Adjust the energy level from low to high level. On condition of safety, the

4. higher the energy level is, the better result it will be. After treatment, it is normal to have skin redness, suggest to use ice pad after treatment.

5. After the phototherapy treatment, use moisturizing preparation and sun block product.

6. After treatment then power off, please use wet cotton swab or cloth to clean treatment window. This can avoid hair and sebum of carbon sticking glass and affecting energy of shooting light.

7. Please don't let system and accessories touch liquid or put it in humid environment. Please store instrument in ventilated and dry environment and keep far away from children.

8. Please don't put metal object into heat emission hole. Otherwise there will be risk of electric shock.

Before starting the therapy

1. Clean the cartridge light output window with wet cloth or cotton swab.
This is done in order to remove hairs and dust from the treatment window.
2. Please hold lamp above and below, then insert or pull out parallelly. (figure 1)
3. Pull out lamp. Install the lamp.
4. Connect the power cord to the B-Shine device.
5. Place the B-Shine device near an electrical outlet. Away from all kinds of liquids.
6. Make sure the lamp is properly installed before turning it on.
After connecting the power cord to the wall electrical outlet, status indicators of device's two sides will flash in green slowly.
7. Press the power button above 2 seconds to get device power on/off. After turn on, then the fan starts and LED lights. Lamp type, remaining pulse, energy level (default is level 1), thermometer icon, fan icon will be shown on the LED (as figure 2).
8. The status indicator flashes which means the device is at the ready condition. After ready the indicator goes out. When the treatment window touches skin, the indicator turns green which means the device is ready for emitting lights.

(Figure 1)


Pull out lamp

Install lamp

(Figure 2)

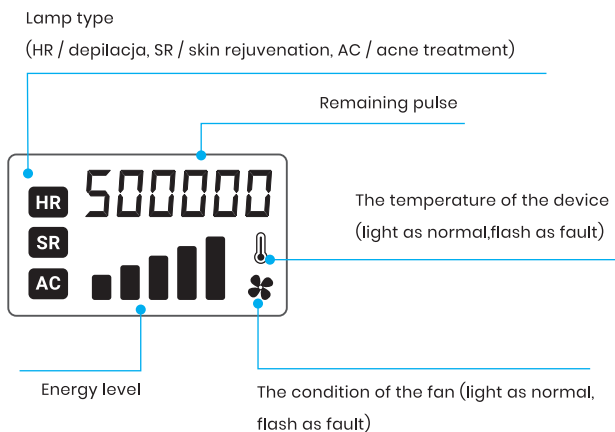9. Press the power button in a little while to adjust the energy level from low to high, level 1 is the lowest and level 5 is the highest.
10. Let the treatment window touch on treated area tightly and vertically. (figure 5)
11. The device will power off automatically if there isn't any operation in 5 minutes.

(Figure 3)

The device has two automatic and manual modes. Manual mode is activated as the default mode. To start the automatic mode, press the flash button for a few seconds. If you want to return to manual mode, press the button again for a few seconds.

Perform skin sensitivity test

If you are using the device for the first time or in a new area of the body, select a small area near the area you intend to undergo to test the sensitivity of the skin. Use the lowest level of the device. Follow the instructions for the procedure below, but emit no more than 2-3 flashes. After the skin test, wait 60 minutes and check that the skin does not react. If after 60 minutes the skin appears normal and there is no redness or other skin reaction, you are ready to perform the full procedure.

Choose the right intensity

Choose the suitable energy according to your skin tolerance and increase energy from one to five level gradually. If there isn't unsuitable or reaction, you could increase energy gradually.

Don't hold the applicator window still over one spot on the skin and flash multiple pulses on the skin. Doing so may cause too much heat and may result in skin swelling, blisters, or temporary or permanent changes in the color of your skin. After each flash, move treatment window slowly to a neighboring skin spot. Before returning to a previously treated skin spot, wait for at least 20 seconds.

Operation procedure

Hair removal

(Iamp depilation / detailed depilation)

Application

- Lip hair, armpit hair, body hair and legs, hair on places that affecting looks such as hair line on forehead and bikini area.
- The effect is different according to skin color. It is not suitable for white and light color hair.

Skin and hair color comparison table

| skin color \ hair color | white | Ivory tint | beige | light brown | brown | dark brown and black |
|-------------------------|-------|------------|-------|-------------|-------|----------------------|
| black | ✓ | ✓ | ✓ | ✓ | ✓ | ✗ |
| dark brown | ✓ | ✓ | ✓ | ✓ | ✓ | ✗ |
| brown | ✓ | ✓ | ✓ | ✓ | ✓ | ✗ |
| golden | ✓ | ✓ | ✓ | ✓ | ✓ | ✗ |
| red | ✗ | ✗ | ✗ | ✗ | ✗ | ✗ |
| white | ✗ | ✗ | ✗ | ✗ | ✗ | ✗ |
| gray | ✗ | ✗ | ✗ | ✗ | ✗ | ✗ |

✓ applicable ✗ not applicable

How often to perform depilation treatment:

Usually, it is about 7-9 treatments. The first 3 treatments take place once a week. In the next 4-9 treatments, do one treatment every 2-3 weeks.

During the period of maintaining the depilation effect and preventing the hair roots from reactivating, regular treatments should be repeated so that the skin remains smooth. If the treatments are not continued, the hair roots will reactivate and the hair will grow back. In order to preserve the results, it is enough to repeat the treatments every 4-8 weeks, which will ensure that the hair roots will not activate and the hair regrowth will remain inhibited. The intervals between treatments may vary depending on the individual rate of hair regrowth.

How to perform depilation treatment:

- Step 1:** Before using the B-shine device, prepare the skin by shaving the selected area to be treated. The skin should also be washed.
- Step 2:** Put on safety glasses to avoid eye stimulation with IPL light.
- Step 3:** Select the desired intensity (see "Choose the right intensity") and start from a low level.
- Step 4:** Let the treatment window touches treated areas vertically, and treat the areas by turns. Do not repeat again and again in the same place, nor omit any area.
- Step 5:** Start in the initial area and repeat step 4, 2-3 times.
- Step 6:** After the treatment it is recommended to use a moisturizer after depilation.
- Step 7:** Turn off the device and use a damp cloth or swab to clean the treatment window.
- Step 8:** Place the device and accessories in a case or box and store them in a safe place.

Photorejuvenation of the skin

(Iamp wrinkle reduction)

Do you want to keep the skin young and healthy for years? We use our skin throughout the year by subjecting it to various atmospheric conditions. The exposure to the sun is particularly effective. Intense tanning usually leaves marks on the face or neckline, which is often difficult to mask with makeup. Similar discolorations also appear on the backs of the hands over time. How to deal with them?

People who want to visibly rejuvenate the skin and refresh its appearance, even out the color and reduce the stumbling discoloration, we offer photorejuvenation treatments with a modern B-Shine cosmetic laser. Laser therapy B-Shine allows for significant skin smoothing, complete removal of age spots or blushes.

The purpose of rejuvenation of selected parties are short and safe light pulses directed towards the surface of the skin. The light wave stimulates skin cells to produce collagen and elastin, which accelerate the regeneration of the skin. The effect of the therapeutic effect of IPL radiation is the improvement of blood supply and oxygenation of the skin, increase in the number of collagen and elastic fibers, better absorption of moisture, regulation of the dyes and keratinization process.

Effects:

- Elimination of signs of skin aging
- Reconstruction of collagen fibers
- Firm and elastic skin
- Reducing the pores
- Smoothing wrinkles
- Improving skin hydration
- Lightening discoloration

How often to perform photorejuvenation treatment:

To obtain the desired effects, a series of treatments should be carried out. It is recommended to carry out a series of 8 treatments, one treatment each week.

How to perform photorejuvenation treatment:

- Step 1:** Before proceeding with the treatment, thoroughly clean the skin with makeup and cosmetics. Also, wipe the lamp.
- Step 2:** Put on safety glasses to avoid eye stimulation with IPL light.
- Step 3:** Select the desired intensity (see "Choose the right intensity") and start from a low level.
- Step 4:** Allow the device window to touch selected areas vertically and treat them areas in turn. Do not repeat again and again in the same place, nor omit any area.
- Step 5:** Start in the initial area and repeat step 4, 2-3 times.
- Step 6:** After the procedure, it is recommended to use a strong moisturizer.
- Step 7:** Turn off the device and use a damp cloth or swab to clean the treatment window.
- Step 8:** Place the device and accessories in a case or box and store them in a safe place.

Acne treatment

(Iamp acne treatment)

The treatment of acne with IPL is intended for people with moderate acne lesions. IPL light destroys rapidly multiplying bacteria responsible for inflammatory changes in acne. The effect is similar to antibiotics that are used topically or generally. However, they can cause a number of side effects by destroying the liver and stomach.

Laser therapy B-shine, acts directly on the causes of acne. It works on synthesized by anaerobic bacteria, molecules of the so-called porphyrins. These make the bacterium particularly sensitive to light. When porphyrin is exposed to visible light, it becomes chemically active and achieves a state of arousal. This causes the formation of singlet oxygen, which connects to the cell membrane and destroys the anaerobic bacteria responsible for the formation of acne lesions and sebaceous gland cells.

IPL light therapy can not be combined with antibiotic therapy. You should also wait a few weeks after the end of such treatment

How often to perform acne treatment:

The treatment period lasts about 8 - 16 weeks (2-4 months), and the treatments are performed every 1-2 weeks. It is best to perform two exposures of all areas attacked at a distance of about 10 minutes. Starting from the second week, the first signs of improvement are visible (reduction of the degree of inflammation of the affected areas and a slight reduction of the symptoms themselves).

After a series of treatments:

Almost acne spots are almost eliminated.

However, those that remain - they are no longer active - disappear within the next weeks.

After curing, it is recommended to apply Skin Rejuvenation therapy, which is able to stimulate collagen production in order to reduce the amount of acne scars.

During the entire treatment it is recommended to use appropriate skincare cosmetics for oily and acne skin.

How to perform an anti-acne treatment:

- Step 1:** Before proceeding with the treatment, thoroughly clean the skin with makeup and cosmetics. Also, wipe the lamp.
- Step 2:** If there is inflammation in the area that you want to treat, you should cure the skin before using the device.
- Step 3:** Wear eye protection to avoid eye stimulation with IPL light.
- Step 4:** Select the appropriate intensity (see "Choose the right intensity") and start from a low level.
- Step 5:** Allow the device window to touch selected areas vertically and treat them areas in turn. The interval between pulses should be at least 20 seconds.
- Step 6:** After the treatment it is recommended to wash the face with a moisturizer.
- Step 7:** Turn off the device and use a damp cloth or swab to clean the treatment window.
- Step 8:** Place the device and accessories in a case or box and store them in a safe place.

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The crossed-out wheeled bin symbol indicates that in the European Union after the end of the product's useful life it must be disposed of in a separate, dedicated collection point. Do not dispose of these products together with unsorted general household waste

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Beauty Factor Sp. z o.o.
Plac Bankowy 2, 00-095 Warszawa

English

Hereby, Beauty Factor Sp. z o. o. declares that the device equipment Beautifly B-Shine Blush is in compliance with Directive 2014/30/UE.

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TERMS OF WARRANTY

1. The product warranty is valid for a period of 24 months.
2. The warranty covers the period from the date of receipt of the device from Beauty Factor Sp. z o. o. and covers defects resulting from causes inherent in the received device.
3. The device will be repaired as soon as possible, not exceeding 21 working days from the date of acceptance of the device for repair by an Authorized Service Center.
4. The customer has the right to apply for replacement of equipment free from defects, if the Authorized Service finds that the removal of the defect is impossible or requires excessive costs.
5. If only part of the product is defective and can be detached from the part of the product working properly, the Customer's rights resulting from this warranty limit only to the defective part of the product.

In order to benefit from the warranty conditions, please send a complaint notification together with the warranty card to the e-mail address serwis@rgmedia.pl

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.....
NR DOWODU ZAKUPU
NUMBER OF EVIDENCE OF PURCHASE

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