






medivon.

User manual

Collar Simple

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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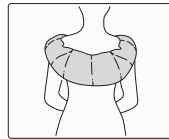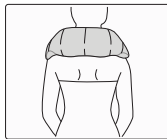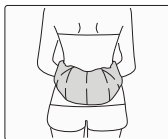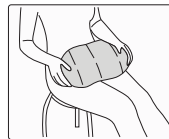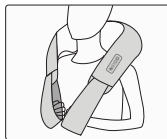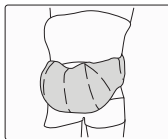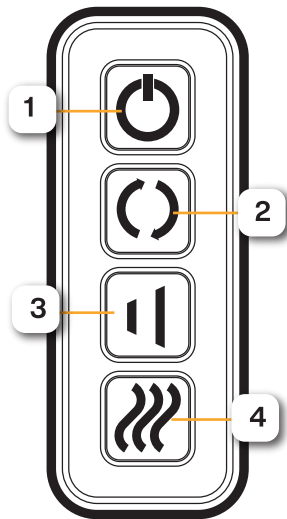Controller support:

1. Power button

Press to turn on the massager, press again to turn off the massager.

2. Directional button

Press while the massager is operating to change the direction of the massage heads.

3. Intensity button

Press while the massager is operating to set the intensity (high or low).

4. Heating button

Press while the massager is running to turn the heating function on or off.

Note: When you turn on the massager, the default massage mode will start. The massage heads will rotate clockwise at a high intensity level along with a warm-up function.

The device will turn off automatically after 15 minutes.

Safety instructions

The massager meets safety rules and technical requirements. Do not wet the device or remove the cover. The massager is not a toy, therefore children and disabled people must be under adult supervision when using it. Do not leave the massager unattended after it is plugged in. Repairs may only be performed by authorized personnel.

Improper use and unauthorized repairs may void the warranty and compromise safety. Avoid contact of the device with water, high temperatures and sunlight. Use only undamaged cables, plugs and sockets. Do not cover the ventilation opening of the device. If the plug or cord is damaged, consult the manufacturer or a specialist. Immediately disconnect the device from the network in the event of a failure. Do not use the massager if you have skin problems, wounds, swelling or inflammation.

The recommended massage time is no more than 15 minutes at a time, after this time the device will turn off automatically. Avoid massage during pregnancy in case of bruises, skin damage, inflammation and swelling. Do not use the massager while driving or sleeping.

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005.

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Short Declaration of Conformity

The manufacturer declares under its sole responsibility that the Collar Simple product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility, Directive 2014/35/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU as regards the list of restricted substances.

The full text of the Declaration of Conformity is available at the following website: www.medivon.pl/deklaracje

Importer: Armare S.A.
ul Grzybowska 87,
00-844 Warszawa, Poland
Made in China
Marka: Medivon
E-mail: serwis@medivon.pl
www.medivon.pl
SN:01265908275917007

Producent: XIAMEN COMFIER TECHNOLOGY CO.,LTD
2-5/F, No.65 Sunban South Road, Jimei District, Xiamen, China.
E-mail: sales@comfier.com

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Shortened warranty terms

This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860. We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: <http://medivon.pl/warranty>