medivon.

User manual

Cosy Shiatsu

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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Failure to observe the above instructions may constitute misuse of the product and can cause serious injury or burns.

- Place the massager on a suitable seat or backrest
- Use the elastic straps and Velcro fasteners at the back massager to attach the massager to the seat / support.
- Connect the adapter output / car adapter output to the socket of the massaged device.
- Turn on the massager with a wired controller.

NOTE: After pressing the ON button / OFF After the first use, the massage will be automatically moved to the lowest body position before starting the massage program. The massage mechanism always returns to this position when the device is turned off. If the power supply is interrupted, the mechanism will also be restored to the lowest position when power returns.

1. Safety Instructions

Please read the following instructions carefully BEFORE putting your massage device into use to ensure trouble-free operation and optimum efficiency. Please retain these operating instructions for further use!

- This product is not intended for medical use, nor shall it be used as a substitute for medical attention. The product is solely intended to provide back massage tusers.
- This product is not to be used by children or invalids without adequate adult supervision.
- Consult your doctor prior to using this product if:
 - you are pregnant,
 - you are fitted with a cardiac pacemaker, implant or other aid,
 - you are suffering from diabetes,

- you have any concerns regarding your health,
- Do not use if you have skin disorders or open wounds.
- Do not use on swollen or inflamed areas.
- Do not use in bed, or while you sleep.
- Discontinue use immediately should you feel any discomfort or pain.
- Never leave the product unattended, especially if children are present.
- Do not use in cars or automobiles when the ambient cabin temperature exceeds 40°C.
- Do not use near open flames.
- Never cover the product when it is in operation.
- The massager will get switched off automatically after 15 minutes. Please rest the unit for 15 minutes before using it again.
- The device has an overheat protection system. If device fails to start, please allow 15-30 minutes for it to cool down before restarting.
- Do not use this product for more than 15 minutes at a time. Extensive use could lead to the product's excessive heating / wear and tear and shorten life.
- Remove adaptor from socket when not in use.
- To minimise risks of electric shock, do not use on wet skin.
- Keep product away from moisture. Do not use if the product or any of its components or accessories are damp or wet.
- Product is intended for indoor use only.
- Do not dismantle or attempt to dismantle or replace product components.
- Use only with the supplied accessories. Mains voltage must be compatible with the voltage indicated on the supplied power adaptor.

3. Controller instructions

1. Power button:

Press to turn the massager on or off.

2. Demo button

Quick demonstration of different cushion modes

3. Heat button

Press to turn on heating. The function is active with the massage option turned on.

4. Vibration button

Press to turn on the seat vibration.


You can switch between low, medium, high vibrations or turn it off completely.

5. Mode button

Use the buttons to select the area to be massaged, upper back, lower or all.

6. Timer button

Press to set the massage time to 5, 10 or 15 minutes.

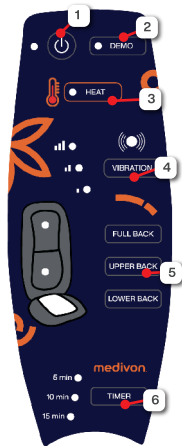4. Care and cleaning

- Do not wash in washing machine or immerse the product or the adaptor in water.
- Always switch off before cleaning and remove the adaptor from the socket with clean, dry hands.
- The massager can be wiped with a slightly damp cloth and dried thoroughly. Do not wipe the product with alcohol or corrosive chemicals such as thinner or bleach.
- Do not overstretch the elastic backstraps.
- Keep away from heat, fire or exposure to direct sunlight.
- Do not apply excessive pressure or impact to the massager/ massage heads. The use of force higher than 50 kg may damage the massager and invalidate the warranty. Refrain from using the device excessively.
- Store the massager in a cool, dry place out of reach of children. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface.

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005.

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Short Declaration of Conformity

The manufacturer declares under its sole responsibility that the Cosy Shiatsu product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility, Directive 2014/35/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU as regards the list of restricted substances.

The full text of the Declaration of Conformity is available at the following website: www.medivon.pl/deklaracje

Importer: Armare S.A.
ul Grzybowska 87,
00-844 Warszawa, Poland
Made in China
Marka: Medivon
E-mail: serwis@medivon.pl
www.medivon.pl
SN:01265908275917120

Producent: XIAMEN COMFIER TECHNOLOGY CO.,LTD
2-5/F, No.65 Sunban South Road, Jimei District, Xiamen, China.
E-mail: sales@comfier.com

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Shortened warranty terms

This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860. We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: <http://medivon.pl/warranty>