GIO-400

CE 1639

innoGIO

Glözoe

SONIC ADULT ELECTRIC TOOTHBRUSH

INSTRUCTION MANUAL

Manual Version: V 1.0

Date of Issue: 04.2020

Introduction

Congratulations on purchasing the GLOzoe sonic toothbrush and welcome to the users of InnoGIO products.

Important informations

- Before using the device's instruction manual.
It is also worth using the manual for the future.
- Use the device only for its intended purpose as described in this manual.

Danger

- The machine being used is damaged in any way (brush head, brush holder or charger), do not use it.
- If the device is damaged, connect to InnoGio Sp. Zoo.

Attention


- Mainly brushes, handles and stands must not be washed in a dishwasher.
- If you have undergone jaw or gum surgery in the last 2 months, contact your dentist before using the brush
- If you use the brush, there is heavy bleeding or if bleeding does not stop after a week of use.
- The InnoGio GIOzoe brush covers all requirements of safety standards related to electromagnetic devices. If you have a pacemaker or other device implanted, you must contact the device or device before joining.
- If you have health concerns before connecting to the InnoGIO GIOzoe brush connect to the application.
- This device is designed to handle teeth, gums and tongue.
Do not use it for any other purpose. If you start to feel discomfort or pain when cleaning with brushes, stop using the devices and consult with the devices
- InnoGio GIOzoe is a personal hygiene device that is not intended for use by a larger number of people in dental surgeries and facilities.
- Do not use the head with broken or bent bristles.
Replace the brush head every 3 months or more often if there are signs of wear. Do not use heads other than those recommended by the manufacturer.
- If the toothpaste contains peroxides, baking soda or bicarbonate (often found in whitening pastes), clean the head of the soapy water brush thoroughly after each use. Prevents plastic from cracking.

General information

Using the innovative Sonic Pulse technology, the InnoGIO GIOzoe brush makes 13,000 high intensity pulsations per minute. Equipped with an ultra-hygienic, medical grade silicone head and with micro movements it breaks and removes plaque in an extremely gentle way, without damaging tooth enamel or gums.

InnoGIO GIOzoe brush

Note: GIOzoe may be blocked before first use. To unlock, hold down the + and - buttons simultaneously. Then the lamp will light to indicate that the brush is unlocked.

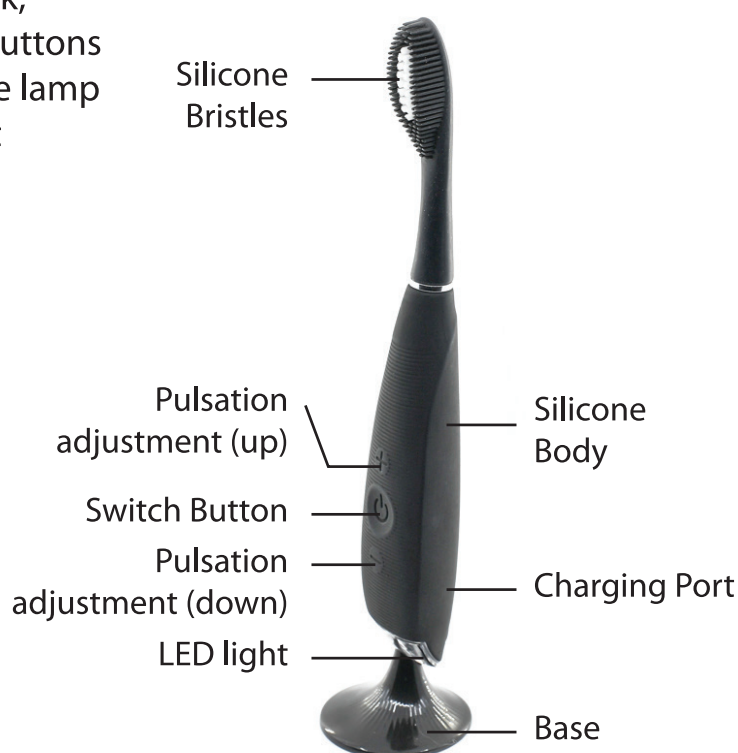Hybrid brush heads

Silicone-Nylon head
with cleaning function

Nylon head for dedicated usage

Silicone head for sensitive usage

Massage head for relaxing usage

Cleaning head
for beauty usage

Properties

ON / OFF button:

- Press once to activate and once to deactivate
- After 3 minutes, GIOzoe will turn off automatically to prevent excessive battery consumption

Pulsation intensity controller:

- 16 modes for different oral care needs: you can adjust the intensity by pressing the + or - buttons
- The device remembers your preferences when you use it again


Hybrid brush heads

- Silicone and Nylon bristles very effectively cleans and at the same time is gentle on the gums.
- Nylon bristle brush heads are intended for adults who have excessive plaque buildup and need more effective cleaning.
- Tongue cleaning heads efficiently and conveniently remove food residues and bacteria from the tongue, providing fresh breath, while ensuring a holistic

In the set you get:

- Silicone-Nylon head with tongue cleaning function
- Silicone head with tongue cleaning function
- Nylon DuPont head
- Peeling head
- Massage head
- GIOzoe main unit
- USB charging cable • The base

How to use GIOzoe

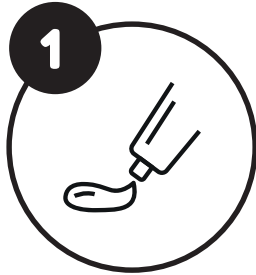Apply toothpaste.

Squeeze the right amount of your favorite toothpaste onto the GIOzoe toothbrush.

Brush.

Press the middle On/Off button to turn on the GIOzoe. To adjust the intensity of the brush pulsation, press + or -. Make sweeping, circular sweeping movements, lightly pressing the toothbrush for 2 minutes. GIOzoe will give an impulse (every 30 seconds) as a signal to move to the next part of the mouth.

Rinse.

Turn off the GIOzoe brush with the middle On / Off button. Rinse out mouth and rinse GIOzoe under running warm water. Repeat brushing twice a day to enjoy a beautiful smile.

Tongue brush – For complete oral hygiene, rotate the brush head onto the tongue brush and gently clean the tongue with vertical movements.

Anti-aging mode for beautiful eyes

Features:

1. Use the massage head.
2. Press the power button to turn on the device. The LED will light up. Press "+" or "-" to adjust the intensity of the device to your preferences (the intensity can be increased or decreased with the "+" or "-" button). We recommend using the weakest pulsing mode for eye area massage.

Recommended ways of performing a massage:

1. **Apply a cream** or serum to the skin under the eyes from the inner corner of the eye to the temple
2. **Care with circular movements:** perform circular movements with the massage head moving it from the inner corner of the eye to the temple. Repeat about 5 times for both sides.
3. **Pressure care:** using the massage head, gently press the skin under the eyes from the inner corner of the eye to the outer corner. The compression time should be about 2 seconds. Repeat 5 times for both sides.
4. **Care by moving:** move the massage head from the center of the lower eyelid to the inner corner of the eye and then to the outer corner of the eye. Repeat 5 times for both sides.

Apply eye cream or essence

Circular movement care

Point pressing care

Movement care

Recommended ways to clean your face with an exchangeable peeling head:

1. Remove makeup, wet your face and apply facial cleanser.
2. Clean your chin and cheeks in circular movements.
3. Starting from the inside, clean your forehead with gentle strokes to the outside.
4. Clean one side of the nose by moving the brush up and down, then clean the other side in the same way.
5. Gently clean the area under the eyes with the pulls out.
To clean the area under the eyes, it is recommended to use the mode of low pulsation.

Rinse and dry your face. Then you can apply your favorite skin care products.

Cleaning the GIOzoe

Clean the brush after each use by washing it with warm water, running your fingers through the bristles. Let the device dry in the open. Placing the brush on the Stand or in the sterilization box.

Note: Never use cleaning products containing alcohol, gasoline or acetone.

Replacement of GIOzoe brush heads

When replacing the head, make sure the brush is turned off. Pull the brush head up and insert the new brush head so that the gap between the brush head and the handle is not visible without twisting it.

Important safety information

The use of this device by or around young children requires close adult supervision.

Stop using if the product is damaged in any way (head, handle, charger).

Never use cleaning products containing alcohol, gasoline or acetone.

To avoid the risk of electric shock, connect the device only to power supply networks with a protective ground.

Charging

NOTE: Before charging, make sure that there is no water or toothpaste residue on the charging port and charger. Connect the charger to the USB charging port (DC 5V). One hour of charging ensures full charging, which is enough for 120 uses. When the battery level is low, it will be signaled by a flashing LED. When the LED lights up steadily, it means the device is fully charged. When the battery is discharged, the engine will stop. DO NOT charge for longer than 24 hours.

ATTENTION: The brush is powered by 3.7V 420mAh lithium batteries, which have a service life of 5 years.

FAQ

1. *Does the brush not turn on despite pressing the On button?*
The battery is dead. Charge the device
2. *The brush is blocked. Press "+" and "-" simultaneously to unlock the device.*
3. *Can't turn off Brushes and / or device buttons are not responding?*
The microprocessor is malfunctioning temporarily. Connect the charging cable to restore normal operation.
4. *The LED does not flash when the toothbrush is connected to the charger?*
The battery is fully charged. The device is ready for use.
5. *The battery is exhausted.* It takes a few minutes for the LED to start flashing while the device is charging. The charging cable is not connected properly. Check the electrical outlet / USB socket and charging port.
6. *The device fell into the water and the charging port is wet?*
Dry and continue using. Our brush is waterproof and the charging port is fully sealed.

Information on the disposal of equipment

Disposal of electronic equipment (applies to the EU and other European countries with selective waste collection systems). The crossed out wheeled bin symbol indicates that this device should not be treated as normal household waste, but should be delivered to an appropriate collection point for the recycling of electrical and electronic equipment.

Removing the battery

The battery should be removed before disposal and should not be disposed of as household refuse.

To remove the battery, remove the brush head, remove the silicone cover, open the plastic housing and remove the battery. Dispose of the battery in accordance with local environmental regulations. Remember that you must first cut off the red positive pole line and then cut off the black negative pole line on the battery.

For safety reasons, gloves should be worn to perform this operation.

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- size: 194 * 29.5 * 25 mm
- weight: 60 g
- battery: 420 mAh, 3.7V
- battery life: 120 use on a fully charged battery
- frequency: 220/120Hz
- maximum noise level: 50 dB
- charging input: usb 5V/2W
- output/voltage: 3.7V/1W

Disclaimer:

Users of the brush take responsibility for its use. Both suppliers and sellers assume no liability through injury or damage, physical or other available equipment, or indirectly from the use of this equipment.

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