SN

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
Instructions for Oclean Portable Oral Irrigator

Thanks for using Oclean Oral Irrigator.

Important Safety Notice

- This product can only be used with Oclean nozzles. We will not bear any responsibility due to the use of any third-party accessory.
- Do not use it for purposes other than oral cleaning, otherwise it may cause accidental injury when the water is sprayed to the eyes, nose, ears or throat.
- Persons who have just undergone oral surgery for treatment. It is recommended that the consumer who has mouth disease, such as seriously hemorrhaged gums, consult a dentist before using this toothbrush.
- For the first time you use it, you may experience discomfort or bleeding, this is a normal phenomenon, and as your gum health improves, the discomfort will disappear.
- If the bleeding is heavy or if it still bleeds after 3-4 weeks, please consult the dentist.
- This product shall be used only with cold or warm water (below 40 °C), or specific solutions as defined in the manufacturers instructions.
- It is not recommended for children under the age of 12. Children shall not play with the appliance.
- This appliance can be used by children aged from 12 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.
- It must only be supplied at safety extra low voltage corresponding to the marking on the product.
- As the water compartment can be unplugged, for safety reason, the product needs to be charged horizontally.
- Instructions for products containing non user-replaceable batteries shall state the substance of the following:
This product contains batteries that are only replaceable by skilled persons.
- Please make sure the charging port is dry before charging.
- Only the Oclean special Type-C charging cable shall be used for the product. Oclean is not responsible for any consequence otherwise.
- Remove the package foam in the water compartment before the first use.

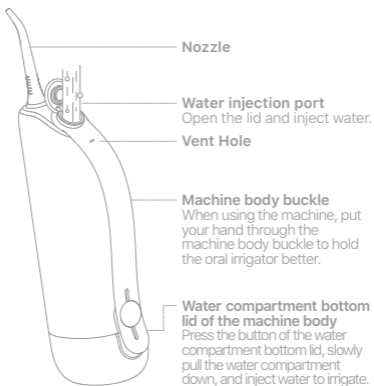Battery indicator light

Red light flashing: low battery alert;
White light breathing: charging state;
White light normally on: the battery is fully charged; the oral irrigator is in working condition

Powerful Normal Mode**Standard Normal Mode****Sensitive Normal Mode****Impulse Mode****Point Shooting Mode****Type-C charging interface**

Open the lid, insert the USB charging cable, and charge.

**Nozzle assembly button**

Press the button and insert the nozzle into the main body of the oral irrigator. Press the button again, and take out the nozzle from the main body of the oral irrigator.

Power button

Perform a short press on the power button to start the oral irrigator, the white light is normally on, and perform a short press on the power button again to shut down.

Mode button**Instruction for use****1 Install the nozzle**

Insert the nozzle into the main body of the oral irrigator.

**2 Inject water into the water compartment**

Lid and inject clean water, the water temperature is appropriate to feel comfortable (recommended not to exceed 40°C), then close the lid, and the installation is completed.

**3 Switching button**

Perform a short press on the "mode button" to switch modes.

4 Preparation before use

Before use, please bend down slightly and extend the nozzle into your mouth, aim at the teeth, gaps among teeth or gums to ensure that the water flow is perpendicular to irrigation positions, and at the same time, open your mouth slightly so that the irrigation water can flow out of your mouth.

Warm tips:

The first time you use the oral irrigator, you may experience itching or slight gum bleeding. This phenomenon usually stops within 1 to 2 weeks, if it exceeds 2 weeks, it may indicate that your gums are abnormal, please stop using the oral irrigator and consult the dentist.

5 Irrigation operation

Press the on-off button to start using this product. This product has a 2-minute timed shutdown function, if you need to continue, please press the on-off button again and make sure that there is water in the water compartment. If you need to pause during use, please press the on-off button, and if you need to continue, please press the on-off button again.

Water volume reminding function

The strength of the water flow impulse changes twice to remind the water volume change, and the water flow impulse returns to the normal value immediately after the reminding.

When the water compartment is filled with water by default: in the Powerful Normal Mode, Standard Normal Mode, Sensitive Normal Mode and the impulse mode, the oral irrigator will remind you each time when arriving at 3/4, 1/2, 1/4 of the theoretical working time of the current gear, reminding 3 times in total.

For example, in the straight irrigation mode, the oral irrigator can continue to spray water for 60s. After startup, the reminding is performed once respectively at 45s, 30s, and 15s.

Note: the "remaining water reminding" function is only enabled in the straight irrigation mode and the impulse mode, and the function is not enabled in the "point shooting mode".

Cleaning Method

Remove the nozzle from the main body wash the nozzle after each use to remove the residue.

After use, press the button of the water compartment bottom lid and slowly pull it down, and if there is residual water in the water compartment, please empty the residual water. Open the bottom water compartment and pour out all the residual water in the water compartment.

After cleaning, wipe the surface of the main body with a dry soft cloth.

Charging

During use, when the battery is lower than 15%, the battery indicator light starts to flash the red light, please charge in time. Uncover the charging waterproof lid at the front of the machine body and insert the USB charging cable to start charging. In the charging state, the battery indicator light is in the white light breathing state; when the charging is completed, the battery indicator light is in the white light normally-on state.

Nozzle instructions

Standard nozzle
For daily cleaning

Periodontal pocket nozzle
For periodontal pocket cleaning

Orthodontic nozzle
For dental brace cleaning

Tongue scraper
For tongue cleaning

Product Specifications

Product Name	Oclean W10 Portable Oral Irrigator
Model	W6010
Water Compartment Capacity	200ml (±5%)
Input Power	5V $\overline{=}$ 1A
Charging	USB Type-C charging cable
Input	DC 5V/1A
Recommended Charger Specifications	Charger with charging current of DC 5V/1A
Service Period	30Days (*)
Battery Capacity	1400mAh (**) Lithium ion battery
Charging Time	4Hours (**)
Product Size	72×55×267mm (with nozzle height)
Product Net Weight	278g
Executive Standard	IEC/EN 60335-1, IEC/EN 60335-2-52

* Note: 1. Power is input power;
2. Irrigate teeth once a day in the impulse mode, 1 compartment of water each time.

** The battery capacity theoretically has a capacity tolerance of ±50mAh.

If you need more detailed product usage information, please follow the Oclean Official Account

www.facebook.com/ocleanglobal

Important use tips

1. Please do not operate this product in a way other than this manual.
2. Please do not use this product without installing the nozzle.
3. Persons (including children) who lack relevant experience and knowledge, such as persons who do not have full civil capacity must use, clean and maintain this product under the supervision and guidance of a guardian.
4. Only the nozzles and accessories supplied, sold, and recommended by Oclean are allowed to be used.
5. Please do not disassemble this product or change the parts of this product in order to avoid damage or danger.
6. Please do not use this product for purposes other than oral cleaning.
7. Please do not press the nozzle forcefully on the gums or teeth in order to avoid injury.
8. Please do not use or operate this product in places where aerosols (sprays) are used outdoors or where oxygen is supplied.
9. Please remove all tooth decorations when using this product.
10. Please do not use iodine, salt water or concentrated essential oils that are insoluble in water in the product. The use of these liquids will reduce the performance of the product and shorten the service life; please do not inject pure mouthwash liquid into this product, if you need to use mouthwash, please dilute it according to the proportion before injection, and after using the diluted mouthwash, you need to irrigate with clean water again to clean the internal pipes.
11. Please do not use the main body completely immersed in water, and please wipe the main body with dry cloth after use.
12. If you will not use this product for a long time, please empty the water compartment to a dry state, clean up the nozzle, main body and charging cable, and store it in a cool and dry place.
13. There may be slight plastic odor during the first few uses of the new machine, this is a normal phenomenon, and after a period of use, the odor will disappear.

Please request a maintenance report from the maintenance center and keep it together with the warranty certificate.

Shenzhen Yunding Information Technology Co., Ltd., in strict accordance with the Responsibility Regulations on Maintenance, Replacement and Sales Return of Household Appliances, offers a guarantee for consumers who purchase the Oclean Oral Irrigator. In the case of any inconsistency between the service articles and the provisions of guarantee, the provisions of the guarantee and the contents of the guarantee certificate shall govern.

If the product you have purchased is found to have a fault or defect, please contact our service line at service@oclean.com

Brand Owner: Shenzhen Yunding Information Technology Co., Ltd.
#28G, Building 3, Dachong Business Center (west of Tonggu Road and north of Shennan Avenue), Yuehai Subdistrict, Nanshan District, Shenzhen, China

Manufacturer: Shenzhen Risun Technology Co., Ltd.
Building A, No.6 of Xinmu Road, Xinmu Community, Pinghu, Shenzhen, China

Website: www.oclean.com

Name of parts	Toxic and harmful substance or element						Environment Protection
	Pb	Hg	Cd	Cr6+	PBB	PBDE	
PCBA	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	10 years
Plastic parts	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	
Metal hardware	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	5 years

Name of parts	Toxic and harmful substance or element						Environment Protection
	Pb	Hg	Cd	Cr6+	PBB	PBDE	
Motor	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	10 years
Accessories	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	
Battery	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	<input type="radio"/>	5 years

Remarks

- This indicates that the content of toxic and hazardous substances in all homogeneous materials of this component is less than the limited quantity as specified in Directive 2011/65/EU(Rohs).
- This indicates that the content of toxic and hazardous substances in at least one of the homogeneous materials of this component is more than the limited quantity as specified in Directive 2011/65/EU(Rohs).
- * Except for very few materials, due to technical development reasons, most of the materials used are environmentally friendly. Please use this product with confidence during the environmental protection period.

Guarantee Certificate of Oclean Portable Oral Irrigator

User Information

Name: _____ Phone: _____ E-mail: _____ Retailer: _____ Phone: _____ Invoice no.: _____ Product type: _____

Address: _____ ZIP code: _____ Address: _____ ZIP code: _____ Sales date: _____ Manufacturing site: _____

Retailer Information

Product Information

Warranty Card of Oclean Portable Oral Irrigator

Delivery date	Fault description/reason/disposal	Equipment employed	Acceptance date	Serviceman

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Proprietario del marchio: Shenzhen Yunding Information Technology Co., Ltd.
Produttore: Shenzhen Risun Technology Co., Ltd.