
medivon®

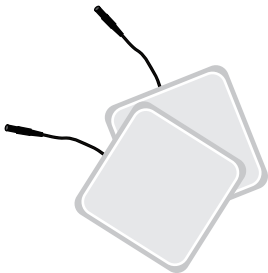Electrodes pad

User manual

PL | EN | DE | ES | IT | FR | CZ |
SL | PT | LT | RU

1. Opis produktu

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1. Product description

1.1 Structure and composition:

The product mainly consists of a conductive gel silica and carbon film and holds fasteners for connecting the electrode to main product.


1.2 Application:

The product is intended to be used in conjunction with a therapeutic device to conduct electrical impulses onto the skin, manufactured by Medivon.

2. Main technical parameters

- Dimensions of the conductive part of the electrode: 50 × 50mm;
- Conducting electrode impedance: 50-150Ω, permissible tolerance ± 10%
- Product service life: 2 years;
- Transport and storage conditions
- Temperature range: 0 ° C~45 ° C;
- Relative humidity: Range ≤80%;
- Atmospheric pressure range: 50 kPa~106 kPa.

3. Illustrating the product

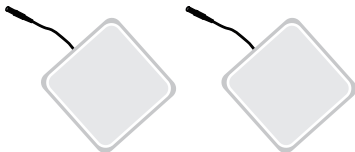4. Application method

- Connect the dedicated cable to the conductive viscose leather electrode (the short name is electrode). Connect the other end of the cable to the therapy device. Remove the protective foil from the electrode, then place the electrode near the painful area (Ashi point), turn on the therapy device and use it as directed.
- When in use, electrodes should be closely in direct contact with the skin and evenly spaced.
- The electrodes can be applied to the shoulder, legs, arms and other parts of the body, but only in areas without wounds or scars. Do not use them on these and other diseased areas of the skin.
- After use, stick the electrode to the dedicated plate and store it in a dry and ventilated place until the next use.

- The electrode surface should be kept clean, avoid any dirt, oily or sticky substances etc. Otherwise, the stickiness of the electrode patch will decrease. If the electrode surface is dirty, it can be washed with clean water and reused after drying.
- The electrode should only be used by the same person. The life of the electrode allows about 80 uses, 20 minutes each. If the electrode cannot regain viscosity after repeated cleaning, consider purchasing new electrodes from distributors or manufacturers.

5. Notes

- Electrodes should not be placed in front of or behind the heart during treatment.
- During treatment, electrodes should avoid contact with wounds and scars.
- Using electrodes close to the chest increases the risk of heart fibrillation.
- Do not use near the heart, head, mouth or perineum, skin lesions, etc.
- Consult a physician before or during use.
- When the electrode is used up, dispose of it in accordance with local environmental regulations.

6. Cleaning, maintenance and storage

- The electrode surface should be in place kept clean. Do not expose the electrode to dirt, contact with an oily or sticky substance etc. Otherwise the viscosity will decrease electrodes. Light dirt on the surface can be rinsed with clean water, use the electrode normally after drying.
- Do not disassemble, repair or modify product on your own to avoid accident or malfunction.
- Store the product in a place:
 - inaccessible to children,
 - not exposed to sunlight, high temperature and humidity,
 - dry and ventilated.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment.

Shenzhen Kentro Medical
Electronics Co., Ltd, 2nd Floor, No.
11, Shanzhuang Road, Xikeng
Village, Yuanshan Street,
Longgang District, Shenzhen City,
Guangdong Province, China

[Tel]: +86(755)33825998

[Fax]: +86(755)33825996

WellKang Ltd
(www.CE-marking.eu)
Enterprise Hub, NW Business
Complex, 1 Beraghmore
Road, Derry, BT48 8SE,
Northern Ireland

[Tel]: +44(33)3303 1126

[Tel2]: +44(20)32876300

[Email] AuthRep@CE-marking.eu

Importer

Armare S.A.
ul Grzybowska 87,
00-844 Warszawa
Made in China