
medivon®

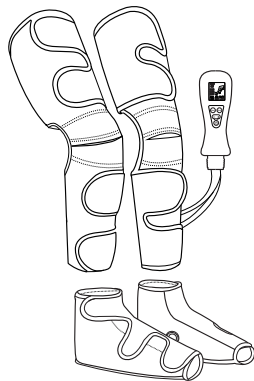
Air massager

Pressure therapy

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

1. Panel Sterowania

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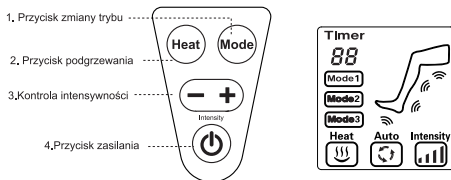
1. Control Panel

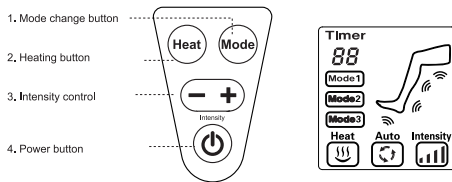2. Operation of the device

Press this button to turn on the machine when it is in standby. The device will start up with the LCD screen and backlight. The default state is automatic air compression mode, intensity level 2 and heating level 2 - the higher level. A countdown will begin (default timer: 20 minutes). The LCD display shows the current operating mode, the intensity and the remaining operating time. If the button is pressed when the device is turned on or after 20 minutes of device operation, all massage operations will be stopped and the device will be turned off and go into standby mode. When the device is turned on, the backlight will be automatically turned off 15 seconds after the last correct button press. When the backlight is off, press any button to turn it on and all function buttons will be active.

Press this button when the device is working to switch the pressure massage modes. Each time the button is pressed, the massage mode switches respectively once. The massage intensity is automatically adjusted to level 2. The LCD screen shows the massage mode and intensity. The switching sequence is: Auto> Mode 1> Mode 2> Mode 3> Auto. Mode 1 massages only the foot area, Mode 2 massages only the calf area, Mode 3 massages only the thigh area, and Auto Mode massages the entire leg.

Press this button when the appliance is running to change the status of the heating function. Each time the button is pressed, the heating status changes once. The  icon is displayed when the heating function is on. The  icon flashes at heat level 1 and lights up at heat level 2.

Press this button to increase or decrease the massage intensity. Each time the button is pressed, the air compression intensity changes to a higher or lower level, up to level 4. The corresponding intensity icon will be shown on the LCD display.

This icon indicates the area of the body that is currently being massaged

3. Dangers

1. Do not use spare parts that are not approved by the manufacturer.
2. Unprofessional service should not disassemble and replace parts of the leg massager without permission.
3. If the power cord is damaged, it must be replaced by appropriate professional personnel to avoid a hazard.
4. Do not use or store the product in a place with high humidity or high temperature.
5. After completing the massage, disconnect the power cord.
6. Do not use the product for more than 15 minutes at a time.
7. If any side effects occur, please stop using the massager immediately.
8. If the temperature is too high or the product starts to smell bad or smells burnt, stop using it immediately and disconnect the power.
9. Do not carry the product by pulling on the power cord to prevent the power cord from loosening.
10. This product is not intended for people who have no physical, sensory or mental abilities or have no experience and knowledge in this field (including children).
11. People insensitive to heat should use the heating function carefully - preferably under the supervision of a caregiver.

4. Warnings


If you have any of the following conditions, consult your healthcare practitioner before use:

- Pregnant women, children or the elderly
- People who suffer from skin infections or have reddened and burned skin
- who use a pacemaker or similar implanted medical device;
- Patients with heart disease or malignant tumors, etc.
- People who are treated or who experience discomfort
- People suffering from osteoporosis
- People with high blood pressure or abnormal blood pressure

5. Cleaning and maintenance

Turn off the product and unplug the power cord before cleaning. Wipe the surface of the product with a dry, soft cloth. Do not clean the product and accessories with harsh detergents such as benzene, thinners, and store in a cool, dry place when not in use. Also avoid exposure to direct sunlight or humid environment, otherwise it may cause color fading or product damage.

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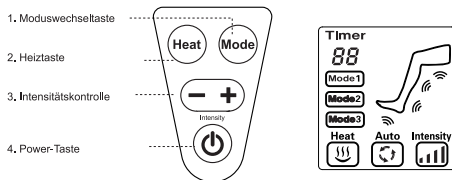The manufacturer declares under his sole responsibility that the Noble Pro device complies with Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, as well as Directive 2011 of the European Parliament and of the Council /65/EU of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to the Directive of the European Parliament and 2011/65/EU as regards the list of restricted substances. The full text of the declaration of conformity is available at the internet address: www.medivon.pl/deklaracje

Importer: Armare S.A.
ul Grzybowska 87, 00-844, Warszawa, Poland.
S/N: 03255904119280252
Marka: Medivon
Made in China
Model: Noble Pro
email: serwis@medivon.pl

Producent: Xiamen High Top Electronic Technology Co., Ltd.
Address: F2-5, No.118, Siming Park, Tong'an Industrial Zone,
Xiamen City, Fujian Province, China.
Email: sales6@xmhightop.com

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005.

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Shortened warranty terms

This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860. We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts:

24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: www.medivon.pl/warranty