
medivon.

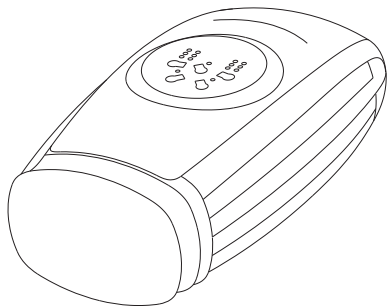User manual

Medivon Hand

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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1. Operation of the device

The device is very easy to use - one press of a button is enough to start a relaxing hand massage which will greatly improve your blood circulation and relax tense hand muscles after a hard day's work.

Put your hand into the massager and then hold the POWER button for about 2 seconds to start. The device's default timer is 15 minutes - after this time, the device will turn off automatically. You can also turn them off by holding down the POWER button again.

Changing the massage intensity - to change the massage level, just press the appropriate button on the left. The intensity will change cyclically from the weakest to the strongest L1 - L2 - L3 - L4 etc.

Changing the massage mode - the massager has defined 6 automatic massage modes. By pressing the MODE button, we can select cyclically changing programs, the one that currently suits us best.

The vibration option is intended for people who appreciate even better effects of relaxation and improved blood circulation in the hands. using the vibration button, we can change their strength or turn it off completely.

Warm-up function - the built-in heads act as a warming compress, this function is activated by the heating button on the right side of the control panel. After starting the warm-up function - the temperature in the device will start to rise gradually.

2. Dangers

1. Do not use spare parts that are not approved by the manufacturer.
2. Unprofessional service should not disassemble and replace parts of the leg massager without permission.
3. If the power cord is damaged, it must be replaced by appropriate professional personnel to avoid a hazard.
4. Do not use or store the product in a place with high humidity or high temperature.
5. After completing the massage, disconnect the power cord.
6. Do not use the product for more than 15 minutes at a time.
7. If any side effects occur, please stop using the massager immediately.
8. If the temperature is too high or the product starts to smell bad or smells burnt, stop using it immediately and disconnect the power.
9. Do not carry the product by pulling on the power cord to prevent the power cord from loosening.
10. This product is not intended for people who have no physical, sensory or mental abilities or have no experience and knowledge in this field (including children).
11. People insensitive to heat should use the heating function carefully - preferably under the supervision of a caregiver.

3. Warnings

If you have any of the following conditions, consult your healthcare practitioner before use:

- Pregnant women, children or the elderly
- People who suffer from skin infections or have reddened and burned skin
- who use a pacemaker or similar implanted medical device;
- Patients with heart disease or malignant tumors, etc.
- People who are treated or who experience discomfort
- People suffering from osteoporosis
- People with high blood pressure or abnormal blood pressure

4. Cleaning and maintenance

Turn off the product and unplug the power cord before cleaning. Wipe the surface of the product with a dry, soft cloth. Do not clean the product and accessories with harsh detergents such as benzene, thinners, and store in a cool, dry place when not in use. Also avoid exposure to direct sunlight or humid environment, otherwise it may cause color fading or product damage.

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators positively contributes to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005.

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Short Declaration of Conformity

The manufacturer declares under its sole responsibility that the Medivon Hand product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility, Directive 2014/35/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU as regards the list of restricted substances. The full text of the Declaration of Conformity is available at the following website: www.medivon.pl/deklaracje

Importer: Armare S.A.
ul Grzybowska 87
00-844 Warszawa, Poland
Made in China
Marka: Medivon
E-mail: serwis@medivon.pl
www.medivon.pl
SN: 02265904119280986

Producent: Xiamen Healthpal Electronic Co., Ltd.
NO.148 BINGLANGXILI, SIMINGDISTRICT. XIAMEN,FUJIAN, CHINA
E-mail: vilin@xmhealthcare.com

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Shortened warranty terms

This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860. We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: <http://medivon.pl/warranty>