BEAUTIFLY

Hydra Peel


PRO

series


Deep cleaning beauty instrument

 /beautifly.eu  @beautifly.eu  beautifly.eu


 /beautiflyeu  @beautiflyeu

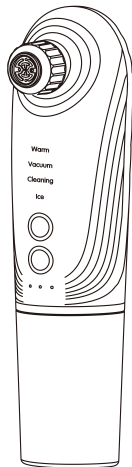
f1

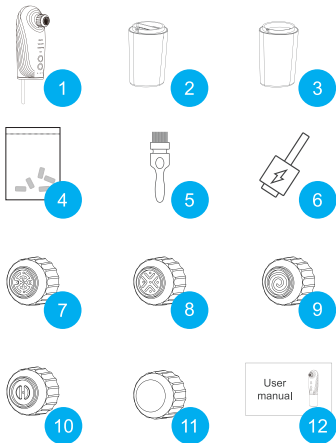
f2




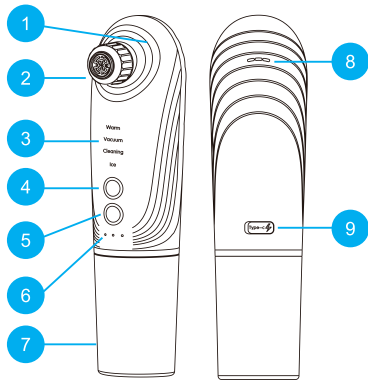
f3

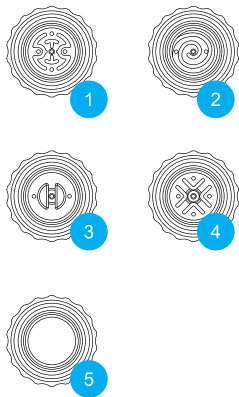f4

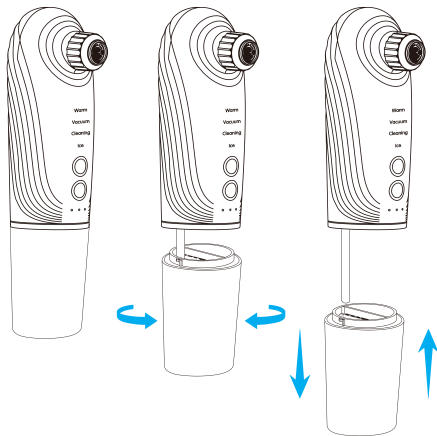f5

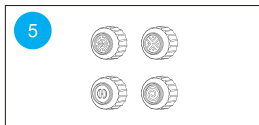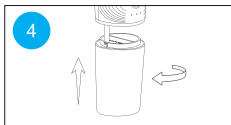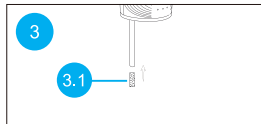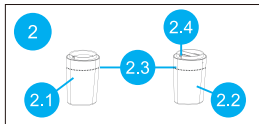f6

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www.beautifly.eu

www.youtube.com/beautiflyeu

Product size: 58x55x207mm
Product net weight: about 286g
Input voltage: 5V DC 1A
Output voltage: 3.7V
Battery capacity: 1500 mAh
Charging time: 3~4 hours
Function modes: 4 modes
Intensity levels: 3 levels
Tank capacity: 2 in 1 tank 40+40ml / all in 1 82ml
Implementation standards: GB 4706.1-2005 / GB 4706.15-2008
/ GB 4706.15-2008

Contains: [page 2, f1]

1. Main unit x1
2. 2 in 1 tank x1 (only for cooling in vacuum mode)
3. All-in-1 tank x1 (only for heating in vacuum/cleaning mode)
4. Filter cap x5
5. Small brush x1
6. Charging cable x1
7. Silicone head x1
8. Silicone cross head x1
9. Spiral head x1
10. Precision head x1
11. Self-cleaning cover x1
12. Instruction x1

Device description: [page 2, f2]

1. Air vent
2. Head
3. Operation mode indicator light
4. Toggle button (tilt indicator light)
5. Intensity change button (charging/low battery indicator)
6. Gear indicator (charging/low battery indicator)
7. Water tank

8. Speaker opening
9. Type C charging plug

Head description: [page 3, f3]

1. Silicone head (large)
suction: excessive (often)
Function: blackhead, dead skin, cutin, pore cleansing (universal type)
location: large areas of the face
area: alar, nasal groove
2. Spiral head (small)
suction: medium (commonly used)
Function: dead skin, cutin, pore cleansing (universal type)
location: small area
area: T zone, nasal sulcus
3. Precision head (small)
suction: strong (persistent dirt)
Function: blackhead, acne, cutin, pore cleansing (blackhead preferred)
location: small area
area: T zone, nasal sulcus
4. Silicone cross head (large)
suction: strong (persistent dirt)
Function: blackhead, acne, cutin, pore cleansing (blackhead preferred)
location: large areas of the face
area: cheeks, forehead
5. Self-cleaning nozzle
function: used to automatically clean the device

- Use point suction or gentle lifting and moving to clean parts that require cleaning. Don't stay in one part for too long.
- After cleansing, due to accelerated skin circulation, slight irritation or slight redness may occur, which is normal and disappears automatically after some time.
- This product is not suitable for plastic surgery, rhinoplasty or other body parts that have undergone cosmetic skin treatments.
- Do not disassemble, repair or modify the device yourself.
- Do not wash or immerse the device in water, you can wipe the device with a damp paper towel after use.
- After using the device, throw away the used water to prevent it from spoiling for a long time.
- Do not place the device in a place exposed to a strong magnetic field.

Function and description of the cleaning head

- Mode 1 - heat therapy: this mode stimulates the skin with red light, opens pores through the thermal circulation system, softens blackheads, acne, cosmetic residues, etc.; [use 2 in 1 tank]
- Mode 2 - vacuum adsorption: thanks to the technology of immediate suction of warm vacuum, it will soften blackheads, acne, grease, cosmetic residues, etc., and then clean them thoroughly.
- Mode 3 - hydrogen-oxygen penetration: after cleaning the pores of the skin, the hydrogen-oxygen water circulation system provides the skin with the moisture and oxygen it needs; [use all-in-1 tank]

- Mode 4 - cold firming: in this mode, blue light calms and tightens, increases cell activity, whitens and rejuvenates the skin, lightens discolorations, narrows pores in the circulatory system of cold therapy, lifts the skin, seals it with moisture and oxygen, etc.; [use 2 in 1 tank]

The above 4 functional modes form an overall closed-loop cleansing cycle system, which makes facial skin cleansing more convenient, thorough and effective, makes the skin full of elasticity and vitality, helps to deeply improve facial skin problems and delay aging skin.

Heat therapy (3 minutes) + vacuum adsorption (3 minutes) + hydrogen-oxygen penetration (4 minutes) + cold firming (3 minutes); Each mode waits for user's instructions (at this time, the next mode light flashes), user presses the mode button to make the device work normally. If the user does not issue commands within 1 minute, the system will automatically enter standby mode.

Automatic cleaning function of the device: after cleaning the facial skin, remove the suction tip of the device, clean it and put it aside, put the cleaning cap on the suction head, fill the tank with clean water, and switch on mode 3 (hydrogen-oxygen penetration), the device will start to clean itself.

How to replace the water tank: [page 3, f4] Rotate counterclockwise, then pull down to remove the water tank, insert up, then right to install the water tank.

Device On/Off/Mode Switch

- Press and hold the power button for 2 seconds to turn on or off the device, short press the power button to switch the function mode.

- Briefly press the intensity button to increase or decrease the strength. One light will illuminate for first gear, two lights for second gear and three lights for third gear.

- When the machine is tilted more than 45 degrees, the tilt indicator light will turn red and an alarm will sound, the machine will stop operating and resume operation when in the upright position.

- During operation, when the two button lights flash red and white alternately and the voice alert indicates that the battery level is low, the product will stop working and the device should be charged at this time.

- The intensity level button light flashes while charging and turns off when fully charged. the charging process takes about 3.5 hours.

Device operation diagram [see page 3, f5]

1. Use facial cleansers to remove makeup.

2. According to the needs of the water tank (All-In-1/2-In-1), the water cannot exceed the maximum water level line.

- 2.1 All-in-1 tank.

- 2.2 2 in 1 tank

- 2.3 Highest watermark

- 2.4 Add water from this side

3. Insert the filter cap into the water inlet pipe.

3.1 Filter caps

4. Attach the water tank to the bottom of the main unit with the logo facing you and screw it clockwise.

5. Attach the appropriate cleaning head. The size of the cleaning head is different according to your own needs.

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www.beautifly.eu

www.youtube.com/beautiflyeu

Product size: 58x55x207mm
Product net weight: about 286g
Input voltage: 5V DC 1A
Output voltage: 3.7V
Battery capacity: 1500mAh
Charging time: 3~4h
Function modes: 4 modes
Suction gear: 3 gears
Tank capacity: 2-in-1 tank 40+40ml / all-in-1 82ml
Implementation standards: GB 4706.1-2005 / GB 4706.15-2008
/ GB 4706.15-2008

Included: [page 2, f1]

1. Host x1
2. 2-in-1 tank x1 (only for cooling vacuum mode)
3. All-in-1 tank x1 (only for heating vacuum/cleaning mode)
4. Filter plug x5
5. Small brush x1
6. Charging cable x1
7. Cross silicone nozzle x1
8. Oval silicone nozzle x1
9. Spiral nozzle x1
10. One word nozzle x1
11. Self-cleaning cap x1
12. Manual x1

Device description: [page 2, f2]

1. Vent
2. Nozzle
3. Function mode light
4. Switch button (tilt indicator light)
5. Gear button (charging/low battery indicator light)
6. Gear indicator (charging/low battery indicator light)
7. Water tank

8. Speaker hole
9. Type-C charging plug

Head description:[page 3, f3]

1. Oval silicone nozzle (large)
suction: excessive (common)
function: blackhead, dead skin, cutin, pore cleaning
(all-round type)
location: large areas of the face
area: alar, nasal groove
2. Spiral nozzle (small)
suction: medium (commonly used)
function: dead skin, cutin, pore cleaning (all-round type)
location: small area

area: t zone, nasal groove
3. One word nozzle (small)
suction: strong (stubborn dirt)
function: blackhead, acne, cutin, pore cleaning (blackhead preferred)
location: small area
area: t zone, nasal groove
4. Cross silicone nozzle (large)
suction: strong (stubborn dirt)
function: blackhead, acne, cutin, pore cleaning (blackhead preferred)
location: large areas of the face
area: cheeks, forehead
5. Cleaning cap
function: used to connect the water circuit of the machine
and clean the whole machine in a circular manner

Precautions

- Please use point suction or gently lifting and sliding to clean the parts that need to be cleaned, and do not stay in one part for too long.
- After cleansing, due to the accelerated skin circulation, the skin will experience some slight irritation or slight redness, which is normal and will disappear automatically after a while.
- This product is not suitable for plastic surgery, rhinoplasty and other parts that have undergone skin cosmetic surgery.
- Do not disassemble, repair or modify this product by yourself.
- The design of this product does not support water washing. after use, you can wipe the body with a wet paper towel, and do not put it in water.
- After the machine is used, the used water must be poured out to avoid deterioration for a long time.
- Do not place it in an environment such as strong magnetism.

Cleaning Head Function And Description

- Mode 1 - warm heat export (warm): this mode wakes up the skin with red light, opens pores through the thermal circulation system, softens blackheads, acne, oil, cosmetic residue, etc.; [use 2-in-1 tank]
- Mode 2 - vacuum adsorption (vacuum): through vacuum hot instant suction technology, it will soften blackheads, acne, grease, cosmetic residues, etc., then thoroughly clean them.
- Mode 3 - water oxygen penetration (cleaning): after the pores

of the skin are cleaned, the hydrogen-oxygen water circulation system supplies the skin with the moisture and oxygen it needs; [use all-in-1 tank]

- Mode 4 - ice therapy firming (ice): in this mode, the blue light calms and tightens, enhances cell activity, whitens and rejuvenates the skin, fades marks, tightens pores through the ice therapy circulatory system, lifts the skin, locks in moisture and oxygen, etc.; [use 2-in-1 tank]

The above 4 functional modes construct a closed-loop overall cleansing cycle system, which makes facial skin cleansing more convenient, thorough and efficient, makes the skin full of elasticity and vitality, helps to deeply improve facial skin problems, and delay skin aging.

Warm heat export (3 minutes) + vacuum adsorption (3 minutes) + water oxygen penetration (4 minutes) + ice therapy firming (3 minutes); each mode waits for user instructions (at this time the next mode light is blinking), the user presses the mode button to it work normally. if the user does not give instructions within 1 minute, the system will automatically enter the standby mode.

Automatic cleaning function of the instrument: after the facial skin is cleaned, please take off the suction nozzle of the instrument, clean it and put it away, and install the cleaning cap on the suction head, add clean water to the integrated water tank, and set it to mode 3 (water oxygen penetration) after turning on the machine, the instrument can perform self-cleaning automatically;

How to replace the water tank: [page 3, f4] rotate counter clockwise, then pull down to remove the water tank, insert

upward, and then clockwise to install the water tank.

Instrument operation steps and status

- Long press the power button for 2 seconds to turn on or off, short press the power button to switch the function mode.
- Short press the gear button to increase or decrease gear, one light will be on for the first gear, two lights will be on for the second gear, and three lights will be on for the third gear.
- When the machine is tilted more than 45 degrees, the tilt indicator light will be red and alarm, the machine will stop working, and it will resume working when it is upright.
- During work, when the two button lights flash red and white alternately, and the voice prompts low power, the product will stop working, and the machine needs to be charged at this time.
- The gear button light flashes when charging, and the light turn off when fully charged. the charging process takes about 3.5 hours.

Instrument operation diagram [see page 3, f5]

1. Use facial cleansers to remove makeup.
2. According to the needs of the water tank (All-In-1/2-In-1), the water cannot exceed the maximum water level line.
 - 2.1 All-in-1 tank
 - 2.2 2-In-1 tank
 - 2.3 Highest water mark
 - 2.4 Add water to this side
3. Insert the filter plug into the water inlet pipe.
 - 3.1 Filter plugs

4. Insert the water tank into the bottom of the main unit, with the logo facing yourself, and screw it into place clockwise.

5. Replace with the required cleaning head. The size of the cleaning head is different according to your own needs. You can choose the appropriate cleaning head in different positions.

6. After cleansing, apply moisturizing products.

More information about the company and products [page 4, f6]:

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under his sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of legislation. Member States in the field of electromagnetic compatibility and Directive 2011/65/EC of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment and Commission Delegated Directive (EU) 2015/863 of 2011. On 31 March 2015, amends Annex II to Directive 2011/65/EU of the European Parliament and of the Council on the list of restricted substances.

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Shortened warranty terms

This Guarantee is granted by a Polish company, Beautifly Sp. z o. o. based in Warsaw. Plac Bankowy 2, 00-095 Warsaw, Poland. You can find us in the official Polish register of enterprises under the KRS number: 0000887306, NIP: 5252854218 and REGON: 388382021 We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: www.beautifly.eu/warranty-beautifly