


medivon.

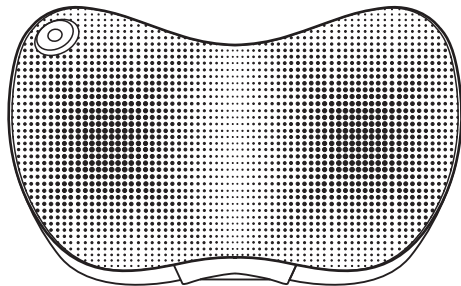User manual

Medivon Travel

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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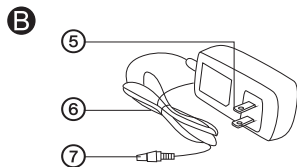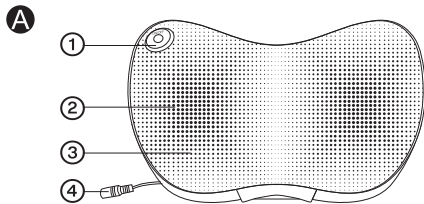Description of the device

1. Massager switch
2. Movable massage and heating heads
3. Material cover for the massage heads
4. Socket for power connection
5. Plug in
6. Power cable
7. Power supply terminal
8. Car adapter

Dangers

1. Do not use replacement parts that are not approved by the manufacturer.
2. Non-professional service personnel should not disassemble and replace parts of the leg massager without permission.
3. If the power cord is damaged, it must be replaced by appropriate professional personnel to avoid danger.
4. Do not use or store the product in a place with high humidity or high temperature.
5. After the massage is finished, unplug the power cord.
6. Do not use the product for more than 15 minutes at a time.
7. If any side effects occur, stop using the massager immediately.
8. If the temperature is too high or the product begins to smell bad or emit a burning smell, stop using it immediately and unplug the power.
9. Do not move the product by pulling the power cord to prevent the power cord from coming loose.
10. This product is not intended for persons who lack physical, sensory or mental capabilities or lack of experience and knowledge in this area (including children).
11. People who are not sensitive to heat should use the heating function carefully
- preferably under the supervision of a caregiver.

Device operation

The device is very easy to use - just press one button to start a relaxing massage of the selected body part.

Pressing the button again changes the direction of rotation of the massage heads. In these two modes, the warm-up function works by default. To turn off the heating, press the button a third time. The device will continue the massage but this time without the warm-up function.

Warnings

If you have any of the following conditions, consult your doctor before use:

- Pregnant women, children or elderly people
- People suffering from skin infections or having skin redness and burns
- who use a pacemaker or similar implanted medical device;
- Patients with heart disease or malignant tumors, etc.
- People who are undergoing treatment or experience discomfort
- People suffering from osteoporosis
- People with high blood pressure or abnormal blood pressure

Cleaning and care

Before cleaning, turn off the product and unplug the power cord. Wipe the surface of the product with a dry, soft cloth. Do not clean the product and accessories with corrosive detergents such as benzene, thinners. Store in a cool, dry place when not in use. Also avoid exposure to direct sunlight or humid environments, otherwise it may cause colors to fade or damage the product.

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Safety Warnings, Maintenance, and Storage Instructions

1. Before using the product, read the operating instructions, safety warnings, maintenance, and storage instructions, and follow the recommendations contained therein. Use the product only for its intended purpose.
2. Before using the product, remove any plastic wrapping and all packaging components and keep them out of the reach of children.
3. This product is not a toy and is not intended for use by children. This product should only be used by adults.
4. Do not use the product if it is not working or is malfunctioning, has been dropped, or has been damaged in any other way. Do not use the product if the power cord or plug is damaged. Do not use the product if it emits an unpleasant odor (e.g., burning), overheats, or emits smoke. In such cases, do not open or repair the product yourself. If the product is still covered by warranty or guarantee, contact the seller or guarantor. If the warranty or guarantee period has expired, or if the defect is not inherent to the product, return the product to an authorized service center for repair.
5. Do not use the product or its accessories near heat, water, or moisture (unless the product is specifically designed for use in water or a humid environment), and do not immerse it in water or other liquids (unless the product is specifically designed for such use).
6. This product is intended for use at room temperature between 10 and 40 degrees Celsius. Do not expose the product to direct sunlight.
7. After use, store the product in a dry place at room temperature between 10 and 40 degrees Celsius, out of the reach of children.
8. If the product becomes dirty, unplug it from the power source (if applicable), wipe it with a dry or damp cloth, and allow it to dry before storing or using it again. Do not use abrasive materials or corrosive cleaning agents, as this may damage the product.
9. If the product is intended for wall or ceiling mounting, ensure it is securely connected after installation. Improper installation may cause the device to fall and damage it, or cause injury to users or others.
10. Do not leave the product unattended when it is turned on. If the device is equipped with a battery, do not use the device while charging. The power source must match the product's required voltage. After use, unplug the product from the power source and allow it to cool before moving or storing it. Do not twist or coil the power cord, as this may damage it.
11. Never cover any ventilation openings on the product and leave at least one meter of free space around the product during operation to allow air to flow in and out of the product.
12. Product accessories that are not supplied by the manufacturer may cause the product to malfunction or be damaged.
13. Additional information regarding massage devices: This product should not be used by individuals who have undergone any medical procedure within the last three months. It is not recommended for individuals with pacemakers or implanted cardiac implants, those who have had a stroke, those suffering from cardiovascular disease, particularly hypertension or hypsomnia, or during pregnancy. If you have any skin complaints, consult a physician before using this product. Before using this product, wait at least one hour after your last meal and remove any jewelry, watches, belts, or other metal objects. Do not use this product while sleeping. If you experience dizziness, weakness, tingling, or pain while using this product, discontinue use immediately and consult a physician.

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Short Declaration of Conformity

The manufacturer declares under its sole responsibility that the product Medivon Travel is in conformity with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

The full text of the declaration of conformity is available at the following internet address: www.medivon.pl/deklaracje

Importer: Armare S.A.
ul Grzybowska 87,
00-844 Warszawa, Poland
Made in China
Marka: Medivon
E-mail: serwis@medivon.pl
www.medivon.pl
SN:01265904119281013

Producent: Xiamen Healthpal Electronic Co., Ltd.
NO.148 BINGLANGXILI, SIMINGDISTRICT.
XIAMEN.FUJIAN, CHINA
E-mail: vilin@xmhealthcare.com

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005.

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Shortened warranty terms

This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860. We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: <http://medivon.pl/warranty>