

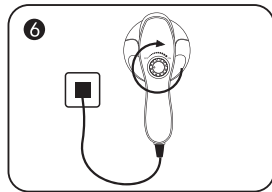
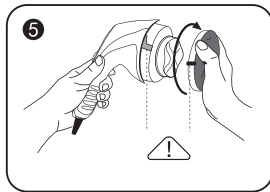
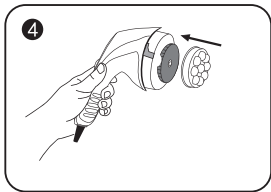
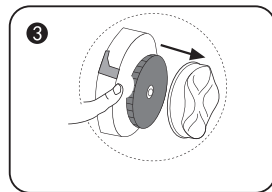
medivon®

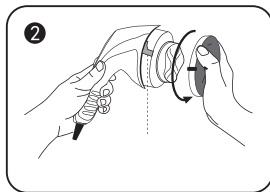User manual

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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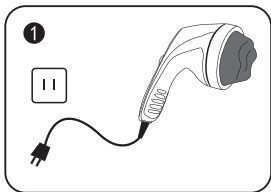Device manual

1. Before replacing the massage head, be sure to disconnect the power
2. Hold the device with both hands, and then turn the cover with the protective mesh inwards - unscrewing it.
3. Find the two gaps between the head and the massager. Then, using your fingers, remove the head by pressing in these places.
4. Select your favorite massage head and then press it to the massager as in the previous step.
5. Reattach the head cover, paying attention to the hooks of both elements
6. Connect the device to the power supply, then turn on the massager with the knob

To clean the massager, use only a soft cloth and neutral detergent.

No water or any liquid may enter the device.

Keep the massager in a dry place after use, avoid high temperatures, direct sunlight and gases that can cause corrosion. Always disconnect the device from the power supply when not in use.

Warnings and Safety Instructions

Do not use the device if you are pregnant or suffering from heart problems or any physical injury. If you are unsure, consult your physician before use. Wait one hour after your last meal before using it. Keep the massager away from children - it is not a toy!

Do not use the device on the same body area for more than 15 minutes in one massage session. Do not use the device in the bathroom or in damp places. If the device does not work properly, immediately disconnect the power supply and contact the service. Do not use the device when you feel unwell. Before replacing the head, make sure that the device has stopped working. If the power cord is damaged, do not attempt to start the device. If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or a similarly qualified person in order to avoid danger to life or health.

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005

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The manufacturer declares, under its sole responsibility, that the ME-800 device complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances. The full content of the declaration of conformity is available at www.medivon.pl/deklaracje

Importer: **Armare S.A.**

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Model: ME-800 | SN:09255904119283680 | Marka: Medivon | email: serwis@medivon.pl | Made in China

Producent: **Shenzhen 3C Linkworld Technology Co.,Ltd**

Room 40203, 4th Floor, Building 1, Huahan Science and Technology Industrial Park,
No. 16 Jinniu West Road, Pingshan District, Shenzhen City, China

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Shortened warranty terms

This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860. We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website:
<http://medivon.pl/warranty>