medivon.

Gun massager

Professional percussion therapy

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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1. Device operation

The device is very easy to use - just one press of a button is enough to start a relaxing full body massage that will allow you to relax after a hard day at work or effectively massage tense muscles after training.

After starting with the power button on the back of the device, you can change the massage modes cyclically, gradually increasing the frequency and power of the vibrations using the + and - buttons

To turn off the device, slide the power button back to the OFF position

Before use, it is recommended to fully charge the device with the charger. Only the charger included in the set should be used for this purpose. A fully charged battery will be indicated by the backlight. It is recommended not to fully discharge the device to extend battery life. It is forbidden to use the massager while charging the battery. It may cause damage or fire.

Note: we replace the massage heads only when the device is turned off. It is enough to remove them from the device by pulling (they do not have a screw thread).

2. Massage and settings

Start by choosing a massage head that is suitable for you and your needs. Make sure it is firmly attached and not damaged. Then start the device by pressing the POWER button.

The strength of the massage is regulated by the user. Depending on the force with which we press the working massager to the body, we will get different final effects. For the first massage, it is recommended to use the massager gently.

Using the M button, we start the heating mode that is adapted to the heating head (the head must be pressed to the very end of the mounting bracket). In this mode, use the + and - buttons to select the temperature of the head in the range of 35 to 60 degrees. The cooling head must be pre-cooled for about an hour in the freezer.

3. Dangers

1. Do not use spare parts that are not approved by the manufacturer.
2. Unprofessional service should not disassemble and replace parts of the leg massager without permission.
3. If the power cord is damaged, it must be replaced by appropriate professional personnel to avoid a hazard.
4. Do not use or store the product in a place with high humidity or high temperature.
5. After completing the massage, disconnect the power cord.
6. Do not use the product for more than 15 minutes at a time.
7. If any side effects occur, please stop using the massager immediately.
8. If the temperature is too high or the product starts to smell bad or smells burnt, stop using it immediately and disconnect the power.
9. Do not carry the product by pulling on the power cord to prevent the power cord from loosening.
10. This product is not intended for people who have no physical, sensory or mental abilities or have no experience and knowledge in this field (including children).
11. People insensitive to heat should use the heating function carefully - preferably under the supervision of a caregiver.

4. Warnings

If you have any of the following conditions, consult your healthcare practitioner before use:

- Pregnant women, children or the elderly
- People who suffer from skin infections or have reddened and burned skin
- who use a pacemaker or similar implanted medical device;
- Patients with heart disease or malignant tumors, etc.
- People who are treated or who experience discomfort
- People suffering from osteoporosis
- People with high blood pressure or abnormal blood pressure

5. Cleaning and maintenance

Turn off the product and unplug the power cord before cleaning. Wipe the surface of the product with a dry, soft cloth. Do not clean the product and accessories with harsh detergents such as benzene, thinners, and store in a cool, dry place when not in use. Also avoid exposure to direct sunlight or humid environment, otherwise it may cause color fading or product damage.

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005

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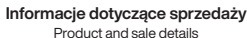
The manufacturer declares under its sole responsibility that the Gun Bandit product is in compliance with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011 of the European Parliament and of the Council /65/EU of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to the Directive of the European Parliament and of the Council 2011/65/EU as regards the list of restricted substances. The full content of the declaration of conformity is available at www.medivon.pl/deklaracje

Importer: Armare S.A. | ul Grzybowska 87, 00-844 Warszawa, Polska

Model: Gun Bandit | SN:12255904119281358 | Marka: Medivon | email: serwis@medivon.pl | Made in China

Producent: Dong Guan Great Ant Innovation R&D Limited

No. 31 Shanglang Road, Xiabian Village, Chang'an Town, Dongguan City, China | e-mail: joannawang@oneaudiodigital.com

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This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860.

We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website:
<http://medivon.pl/warranty>