
medivon.

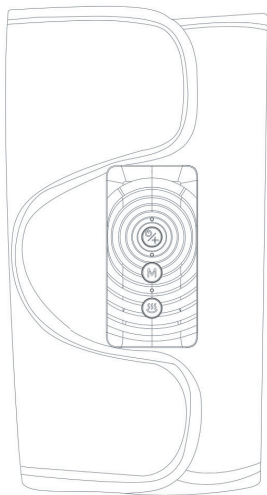User manual

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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Device Manual

Fully charge the device before using it for the first time and after a long period of non-use. After charging the battery, you can hold the ON/OFF button for at least 2 seconds to start the massager. You can freely change the operating modes with the "M" button and the frequency with the "+" button (pressing the button repeatedly changes the options 1-2-3).

After pressing the heating button, the device will start a warming massage.

When the battery level is low, the red diode will turn on, after charging it will turn green. The supply voltage is 5V/2A (USB-C), use this for charging. Do not use the device while it is charging.

After the massage is completed, the device can be turned off by pressing the button again.

The massager is ideal for people who need deep relaxation or improve blood circulation. With systematic use of the massage, you will reduce the symptoms of fatigue, swelling of the feet and calves, or chronic pains of unknown origin. The construction of the massager is very light and fits all calf sizes.

The strength of the massage should be adjusted to your own preferences. The device must not be used by persons with injuries or who have undergone surgery in the last 3 months. Do not use the device in high temperatures and damp rooms. Absolutely keep away from water. Also, people with hypertension, hypotension or severe strokes cannot use the device. People with skin problems should consult a doctor before use.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment. The device was placed on the market after 13.08.2005.

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The manufacturer declares, under its sole responsibility, that the Medivon Dione device complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances. The full content of the declaration of conformity is available at www.medivon.pl/deklaracje

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska
Model: Medivon Dione