
medivon.

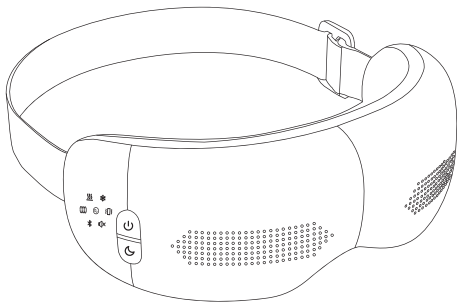User manual

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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Device Manual

Fully charge the device before using it for the first time and after a long period of non-use. After charging the battery, you can hold the ON/OFF button for at least 2 seconds to start the massager.

When the massager is turned on, pressing the power button cyclically changes the day massage modes, and the moon button cyclically switches between night modes.

The diodes on the device indicate which functions are currently active in a given mode: cooling, heating, vibration.

A quick double press of the moon button plays the next track, while holding it down mutes the music.

Recharge the device when the battery level is low. Do not fully discharge the battery. The supply voltage is 5V/1A. Do not use the device while it is charging.

After finishing the massage, the device can be turned off by pressing the ON/OFF button again.

To connect to the massager and play your own music, just search for the device in the list of visible bluetooth devices and connect. The music will play automatically when you select a song.

The strength of the massage should be adjusted to your own preferences. The device must not be used by persons with injuries or who have undergone surgery in the last 3 months. Do not use the device in high temperatures and damp rooms. Absolutely keep away from water. Also, people with hypertension, hypotension or severe strokes cannot use the device. People with skin problems should consult a doctor before use. Pregnant women cannot use the massager.

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005

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The manufacturer hereby declares that the Horizon Lynx device complies with Directive 2014/30 / EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65 / EU of the European Parliament and of the Council of 8 June 2011 on restriction of the use of certain hazardous substances in electrical and electronic equipment. The full text of the declaration of conformity is available at the internet address: www.medivon.pl/deklaracje

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska
Model: Horizon Lynx
Marka: Medivon
SN: 06255904119281501
e-mail: serwis@medivon.pl

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Shortened warranty terms

This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860.

We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: www.medivon.pl/warranty