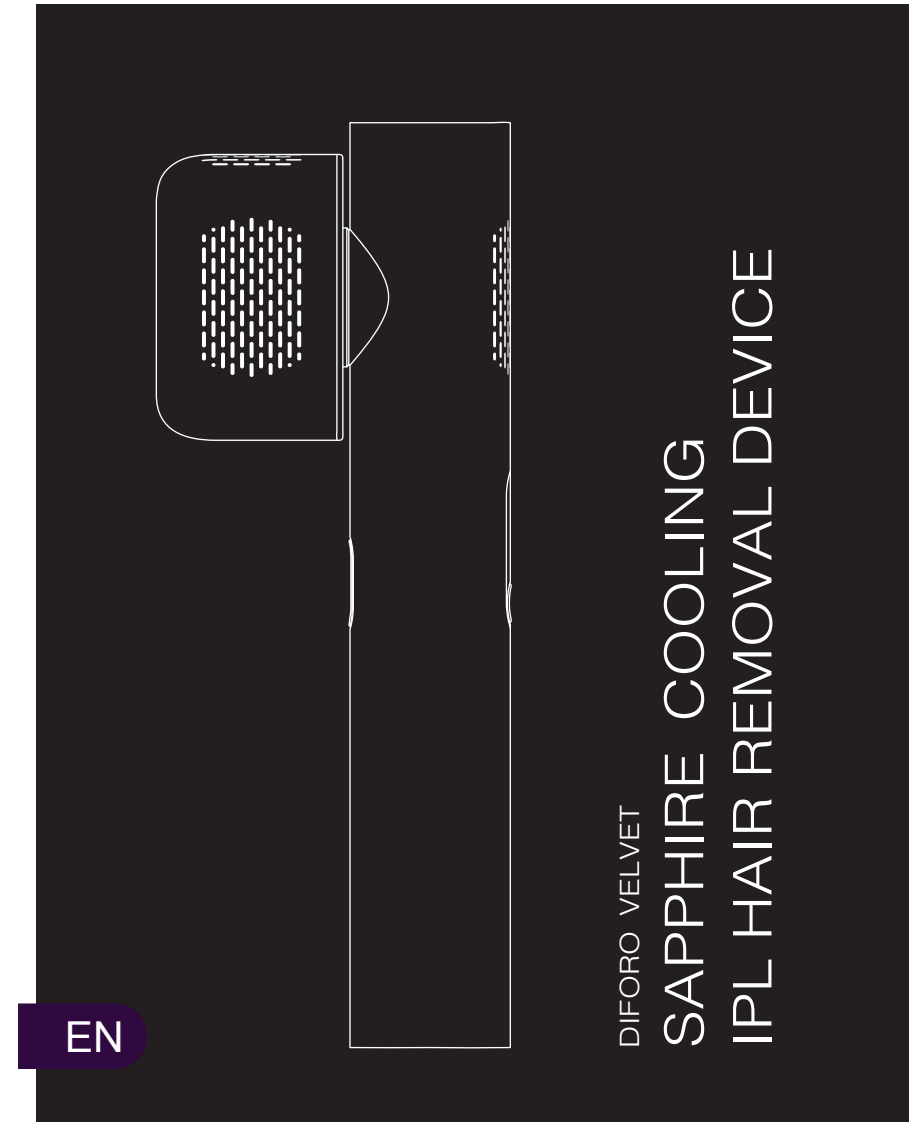
Product principle of operation


Using the principle of photothermal action, the IPL freezing point hair removal device generates high heat inside the hair follicle, thus degrading the structure of the hair follicle. Weakening the hair follicles and removing the hair in a painless manner. The hair removal device has a contact sapphire cooling head, which avoids the heat from spreading outside the hair follicle structure, while allowing the heat to diffuse from the epidermis. In such cases, the skin is not damaged, so there is no pain or discomfort. Blue light has an anti-inflammatory effect and revitalizes the skin, so in addition to the hair removal effect, we strengthen and regenerate the dermis.

Limitations on use

1. Women who are pregnant, breastfeeding and menstruating;
2. When the skin color is black or dark brown;
3. Patients with scars or skin diseases;
4. For very fair skin, always perform a skin test;
5. Those who have recently been exposed to strong sunlight or who cannot avoid sun exposure after depilation;
6. Patients with heart rate regulators, insulin pumps or other medical devices;
7. minors (during the developmental period, rising hormone levels).

Set contents

- IPL hair removal device
- Charger AC
- Protective goggles
- Manual and warranty

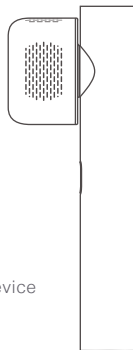Charger AC


Protective goggles

Manual and warranty

- IPL hair removal device

Product specifications

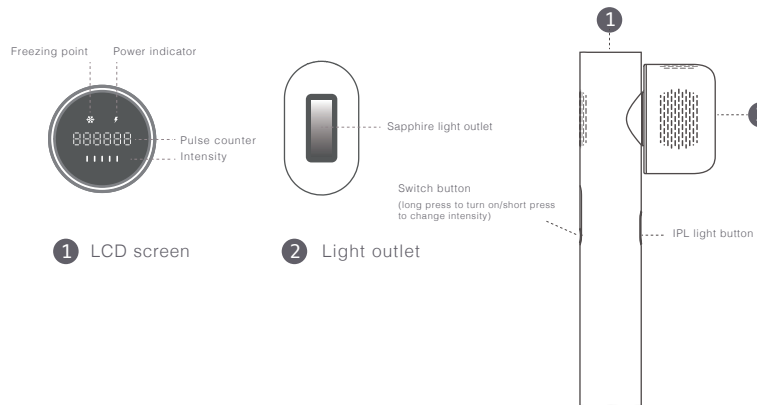Operating modes

Single mode (suitable for small areas such as fingers and hair on the lips):

Place the light outlet on the area of skin to be treated and completely align it with the skin. Briefly press the flash button to start, after the device flashes, move the light head to the next section of skin. Then press the flash button again and repeat the above steps in succession.

Tip: If the skin becomes too hot during the epilation process, just pause for a few seconds and let the light outlet cool down again before proceeding with hair removal.

Continuous mode (applies to large areas such as arms, legs, etc.):


Place the outlet of the lamp on the area of skin to be treated and completely align it with the skin. Press and hold the flash button and immediately start moving the device slowly and evenly across the skin to keep it in close contact with the skin. Flashes will be emitted at specific intervals. During this time, do not stop moving the device and move it forward and backward in the area that needs epilating, make sure the epilation area is evenly covered.

Tip: The skin speed can be adjusted according to the flash speed, to make the coverage more uniform, the flash frequency is faster at the 1st level, slightly slower at the 2nd level etc. The higher the level, the lower the output speed of the light.

Rotation

- The head of the device can be rotated for use, there is no need to twist the wrist and arm, and operation is more convenient. Note that the direction of rotation of the device can only be rotated clockwise (when the head is placed vertically, along the housing) and counterclockwise (when the head is placed horizontally). Do not rotate the guide head in the opposite direction, as this will cause damage to the head under external force and adversely affect its use.

- For hair on the lips, armpits, elbows and other difficult to access parts of the body, the head of the device can be rotated to the "ten" state, which is more convenient. can be rotated to the "ten" state, which is more convenient and more comfortable to use. more convenient to use.

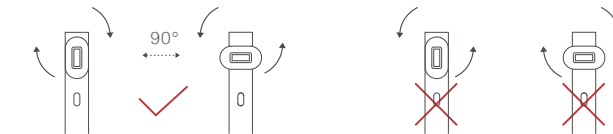Depilation cycles

First cycle: For 1-2 weeks, apply twice a week. Week 3-5, apply once a week;

Second cycle: Another 5-12 weeks, apply once every two weeks.

Third cycle: After this time period, when the hair grows 1-2 mm, one use.

- Sustained use can inhibit hair growth so that hair will not grow back during the growth cycle. The growing hair will become thinner and shorter and significantly personified.
- It should be used with caution, depending on the growth and location of the hair. Hair on the lips and arms is relatively soft, and the cycle of use of the epilator should be shorter. However, hair in the underarms and bikini area is relatively tougher, so the cycle and frequency of use can be increased accordingly. The thicker the hair, the greater the range of use should be.
- Hair removal devices should not be used too often, as this will affect depilation and cause damage to the skin (such as folliculitis). Hair removal time should be strictly controlled. If you experience allergy symptoms while using the device, contact a specialist or dermatologist.

Intensity in relation to skin color

Choose different levels depending on your skin and hair color. The higher the intensity level, the greater the energy and faster the effect. It is recommended to start with the first one and gradually increase to a tolerable, comfortable level. (Except for black skin)


Skin						
Hair						
Light blonde/beige	1-5	1-5	1-4	1-3	x	x
Brown	1-5	1-5	1-4	1-3	x	x
Dark brown	1-5	1-5	1-4	1-3	1-2	x
Black	1-5	1-5	1-4	1-3	1-2	x


Note: Skin shades with an "x" in the table are not suitable for this product. White, gray, red, light blonde hair is not suitable for this product. For different reasons, skin will react differently due to different exposure times.

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How to use

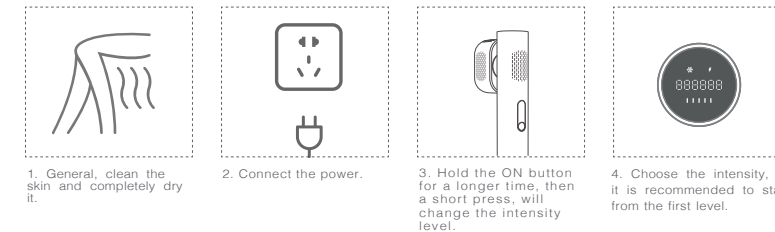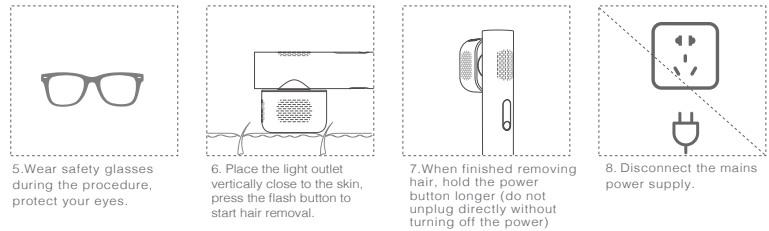Helpful tips

- When using the IPL device, the head of the lamp should fit tightly against the skin to make full contact with the hair follicle root, ensuring that the follicle absorbs more energy and you get a more effective hair removal effect. At the same time, the sapphire touches the skin to pre-cool the epidermis, reduce skin heat and better protect the skin during the treatment.
- When the light outlet is not completely in contact with the skin, the light will not be triggered to avoid accidental illumination and eye damage.
- Using the device in continuous mode may cause the temperature of the case to rise. If this effect is noticed, you can try to stop epilation for 30 seconds. The sapphire glass, during this time will cool down enough to make the treatment more comfortable again.

Caring after use

Precautions after hair removal

- Avoid contact with water for 6 hours after epilation, and use warm water to moisturize the epilation site.
- Do not use perfumes, lotions, masks or cosmetics after epilation.
- Exposure to UV radiation should be avoided for 24 hours after hair removal treatment. This includes artificial tanning in a tanning bed.
- If you experience discomfort after use, consult a specialist.

Cleaning and maintenance of the equipment

- To ensure the best results and longest life of the product, clean it before each epilation. Especially the sapphire filter and power head. Remove dust from the air inlets and outlets.
- The device must be turned off and unplugged before cleaning. The device must cool down.
- Store the device in a dry and dust-free place at room temperature.
- Wipe and clean the sapphire light outlet with a soft, clean and dry cloth.
- Never wash this product and accessories with water.
- Wipe the product with a clean or slightly damp tissue or cloth. If the cleaning effect is not satisfactory, you can wipe the light outlet with low concentration alcohol. The metal body must not be wiped with alcohol. Do not use abrasive harsh cleaners (toothpaste, etc.) or strong liquids to clean the product.
- When discarding the product, take it to a designated recycling center.

Precautions

Be sure to use the product in accordance with the regulations in the "User Manual". If you have any questions during use, please contact customer service online. Improper use may cause damage to the device.

- 1) Do not use in the bath or shower;
2. do not use this product near water-filled bathtubs, shower trays, sinks or other vessels containing water. Store the product and power supply in a dry place;
- 3) If the device is damaged, do not touch its internal parts to avoid electric shock;
- 4) Do not put any non-factory accessories in the device;
5. do not use this product around the eyes. Wearing safety glasses is an activity that protects you when triggering the light. If you lose your glasses, you can use other glasses that have a good sunscreen;
6. it is recommended to wear safety glasses every time you use this product. A flash of light that is not aimed directly at your eyes will not damage your vision, but a flashing light can significantly affect your comfort and cause vision problems;
7. you may need an assistant or partner when epilating your back. Get help when you are not completely sure that you are in full control of the device;

8. First of all, shave the depilation site. We recommend using a manual shaver. Excessively long hair absorbs heat on the surface of the skin and may cause redness and swelling, and burning hair will cause odor and discomfort;
9. If the skin in the irradiation zone is not dry or clean, it will affect the effect of use and may cause the lamp to fail to flare or flare accidentally;
10. We recommend that you exercise caution and do not give the device to children under 15 years of age. IPL hair removal is not recommended during puberty and increased hormonal activity;
11. Use different light intensities depending on your complexion, skin color and hair. Use the chart or online consultation;
12. If the power cord is damaged, it must be replaced by a specialist from the manufacturer's service department to avoid danger;
13. The power head and IPL light pulse window should be checked and cleaned regularly. A dirty device may be damaged during use;
14. Do not place the device near objects that have high electromagnetic fields;
15. Do not use the device when dust or other dirt has clogged the air inlets. Do not place any stickers or fabrics on the inlets. Obstructed inlets will cause overheating and equipment failure;
16. Do not point the light outlet in a direction other than the epilation zone and area;
17. Do not use light pulses through clothes and other objects to avoid burning, discoloration and damage;
18. Do not use against animals and plants;
19. Clean and protect the device thoroughly after use;
20. adjust the intensity of irradiation according to the condition of the skin. Remember to use the first level for the first use;
21. If the skin becomes slightly hot and red after use, this is a natural reaction;
22. The product must not be left unattended both on and off. Pack it up and protect it from children;
23. Do not use the device in rooms where the equilibrium between the environment and the working device will be more than 20 degrees. Such a difference may cause malfunctions. The proper working temperature is 0-45 °C;
24. Do not use near flammable and explosive materials (including sprayers and gas stations, etc.), which may cause fire or electric shock;

25. Do not expose the product to direct sunlight or ultraviolet light for a long time, which may seriously damage the product;

26. If the built-in fan does not work (no murmur is heard) after power on, do not use it to avoid failure. In this case, turn off the unit and contact the manufacturer's service department.

Do not use in particularly sensitive areas:

1. Head

2. Eyebrows

3. Ears

4. Neck and neck

5. Eyes and eyelids

6. Acne area
7. Nipples and surrounding areas

8. Inside the navel

9. Intramammary, genitals, mucous membranes, etc.

10. After a surgical procedure.

11. In the areas of tattoos, body painting.

Product parameters

Product name	Freezing Sapphire IPL Hair Removal	Intensity level	5 levels
Product model	VELVET	Energy efficiency	7.5~23J
Wavelength	550-1200nm	Power	48W
Dimensions	89*40*233 mm	Rated Inpu	24V-2A
Weight	95 g	Pulse duration	800-850µs

Warning

Always remember to unplug your IPL device from the electrical outlet after treatment. Take care of your safety and the safety of others.

Do not use this IPL device in bathtubs, showers, sinks, or other vessels that contain or come into contact with water.

Do not dispose of electrical appliances as unsegregated municipal waste. Contact your local municipality or manufacturer for information on available collection systems.

Warranty card

The manufacturer guarantees its products for a period of 24 months from the date of purchase. Accessories included with the original product are warranted for a period of 12 months from the date of purchase.

We recommend smooth operation of products and complaints.

- ① Carefully read and follow the instruction manual.
- ② Keep the instruction manual, warranty card and proof of purchase.
- ③ Use the original accessories supplied with the device.
- ④ If there are problems with the device, contact the manufacturer.
- ⑤ Do not make repairs or modifications yourself.
- ⑥ Use the product for its intended purpose.

During the warranty period, all repairs are performed free of charge.

Any defects or faults will be repaired within a maximum of 14 days, calculated from the time of delivery of the device to the manufacturer's authorized service center.

Limitations:


- The warranty will not be recognized in case of mechanical damage caused by inattention or fortuitous events.
- If the device is used contrary to its intended use and specifications.

DIFORO Brand
CUSTOMER SERVICE **+48 453 271 100**

Consultation hours : 08:00 ~ 16:00 (8 AM ~ 4 PM)
(Monday to Friday)

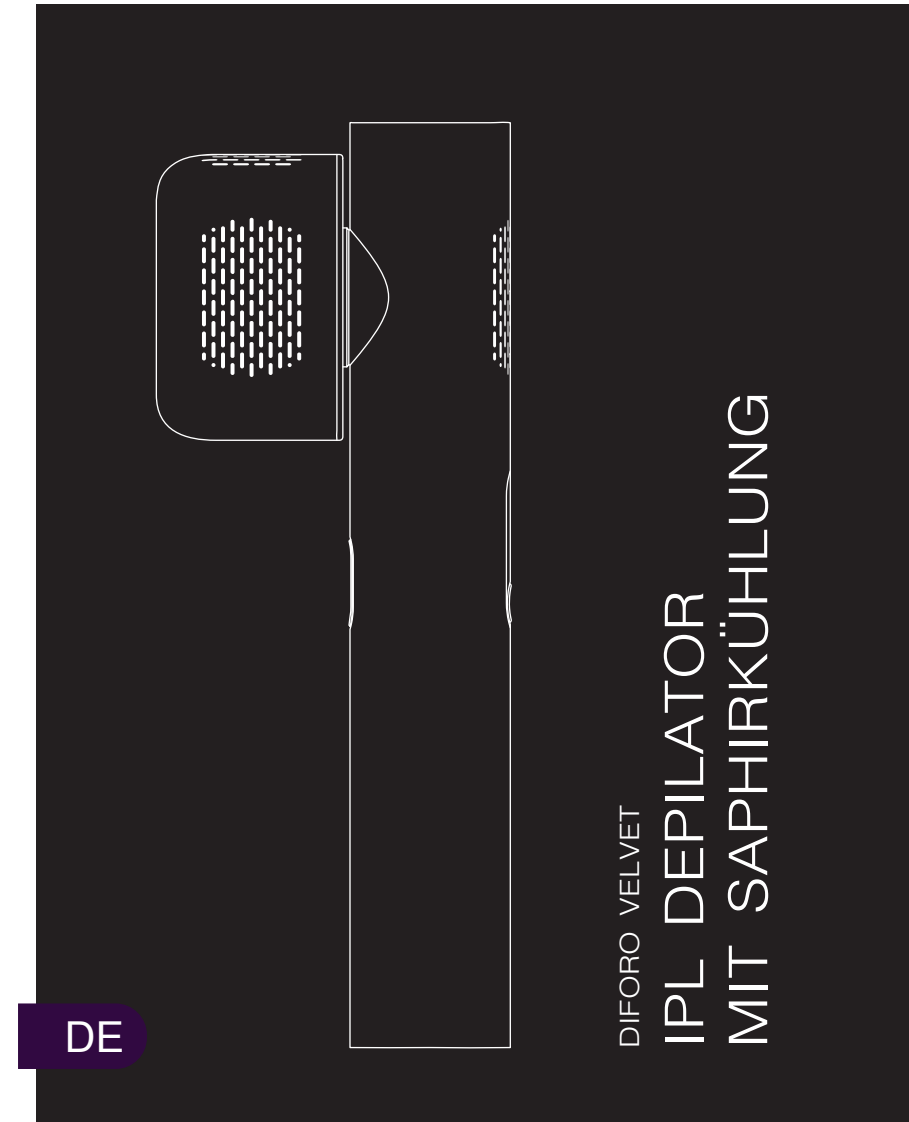PRODUCENT
DIFORO by VOBÉ sp. z o.o.
Rudna Mała 167A, 36-054 Mrowla, POLAND
www.diforo.pl get@diforo.pl
NIP: 5252799270, KRS 0000800189
Made for Diforo in China

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