
medivon.

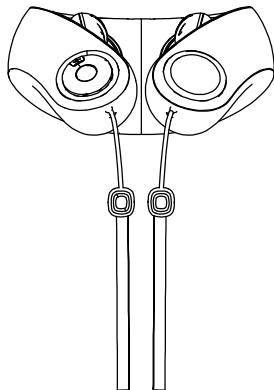Neck massager

Shiatsu therapy

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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The device is easy to use – simply press the button to start the massage. Hold the POWER button for 2 seconds to start the massage; holding it again activates the default massage mode. The adjustment button allows you to change the speed of the head, and another button allows you to adjust the heat level or turn it off. Before first use, charge the device via USB. Do not use the massager while charging – this may damage the battery. It is recommended not to completely discharge the battery to prolong its life.

Begin the massage in a comfortable position, and people with long hair should tie it back to avoid getting tangled in the head. The massage intensity depends on the pressure applied, so be careful, especially during first use. The massager is equipped with speakers – you can connect to a mobile device via Bluetooth and listen to music during the massage. A short press of the music button changes the track, a double press turns it off, and a long press plays the music.

Do not use parts not approved by the manufacturer. Replacement or repair should be performed by qualified personnel. A damaged power cord must be professionally replaced. Do not use the device in wet or high-temperature environments, or for longer than 15 minutes at a time. Unplug the power cord after use. In the event of a malfunction or overheating (burning odor, high temperature), turn the device off immediately. Do not pull the cord. This device is not intended for use by individuals with physical or cognitive disabilities without supervision. Individuals who are heat-sensitive should exercise particular caution when using the heating function.

Consult your doctor before use if you are pregnant, have skin infections or injuries, have implanted medical devices, have heart disease, cancer, osteoporosis, high blood pressure, or are undergoing treatment for back pain.

Unplug the device before cleaning. Wipe with a dry, soft cloth. Do not use harsh chemicals or solvents. Store in a cool, dry place, away from sunlight and moisture – otherwise, discoloration or damage to the casing may occur.

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005

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The manufacturer hereby declares that the Collar Grand device complies with Directive 2014/35/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65/EU of the European Parliament and of the Council of June 8, 2011 on restrictions on the use of certain hazardous substances in electrical and electronic equipment. The full content of the declaration of conformity is available at the following internet address:

Importer: Armare S.A.
ul Grzybowska 87,
00-844 Warszawa, Poland.
Made in China.
SN: 07255904119282997
Marka: Medivon
email: serwis@medivon.pl
Model: Medivon

Manufacturer: Ningbo Fulljoy Electronic Technology Co., Ltd
Address: No. 108, Huiyuan Road, Fenghua District, Ningbo, China
email: chengya3@fulljoy.com.cn