EMS FOOT MASSAGER

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
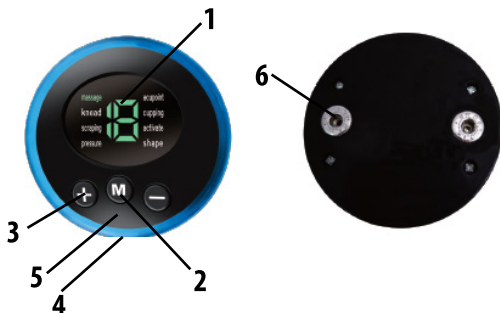
Mat and remote control

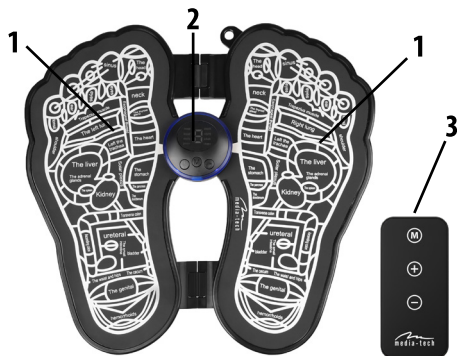1. The location of the user's feet so that they adhere to the largest possible surface
2. Mat power latches for controller placement
3. Remote control with +/- adjustment buttons, ON/OFF button and mode change button

Charging and installation

Before starting the massage, we recommend charging the controller battery to 100%. To do this, connect the socket [4] with the included USB/microUSB cable to the 5V/1-2A charger or USB port on the computer. Do not use the so-called fast chargers. **Do not use the device while charging**
The charging of the device is signaled by the LED [5], which will turn off after charging is complete.

Remember to **fully** charge the device after each use and charge it at least once every 3 months, even when the device is not in use, to keep the battery in good condition.

Remove the battery discharge protection from the remote control. Connect the controller to the charging latches so that a "click" is heard

Usage

- a. Moisten the soles of your feet with a damp towel for a better massage effect. The procedure should be performed sitting down so that the leg muscles are unburdened and relaxed. Do not put your full weight on the mat
- b. Place your feet on the designated places on the mat so that they adhere as much as possible
- c. Use the [+] button to turn on the device
- d. Use the mode button [M] to select one of the massage modes
- e. Use the +/- buttons to set one of the 19 levels of intensity.

The mat will turn off automatically when one foot is removed.

The massage involves the use of mild electrical impulses (EMS) to stimulate the muscles of the feet and calves. This helps to improve blood circulation, relieve fatigue and muscle tension, and stimulate the acupressure points on the feet. The mat has different massage zones, covering the sole, heels, toes and arches of the feet.

The first use may require getting used to a slight sensory discomfort, especially in people with sensitivity to the so-called. "ticklers".

We recommend that you select the massage intensity level "1" the first time and find the treatment modes that best suit you.

EMS treatment is recommended maximum once a day, up to 15 minutes of continuous massage. The device has a protection that automatically turns off the treatment after about 15 minutes.

IMPORTANT!

**The MT6527 is not a medical device and cannot replace recommended therapy in medical centers.
In case of doubt, it is recommended to consult a doctor.**

IMPORTANT!

The mat works based on EMS muscle electrostimulation technology with low current. People who have implanted electronic circulatory support and heart function, as well as problems with arrhythmia or after a heart attack, should consult a doctor about the safe use of this device.

Before using the device, medical consultations should be performed by people with diagnosed arterial hypertension and atherosclerosis.

FEATURES

- EMS electrostimulation foot and calf massage mat
- Improve blood circulation and regeneration of tired legs
- Very easy to use
- Durable, hard, foldable design for use anytime and anywhere
- 8 massage modes
- 19 intensity levels
- Control module with LED display
- Remote control for greater convenience of use
- Power supply: rechargeable 150mAh lithium-polymer battery for up to 160min on a full charge
- Battery charging: micro USB 5V cable ~1 hour
- Automatic shut-off when feet are removed from mat
- Dimensions: 31x29cm
- Weight: 0.24kg

Useful tips

Dry skin significantly reduces the effect of EMS pulses.

We recommend wetting the soles of your feet with a damp towel each time to improve contact.

Important!

For the mat to work, a closed circuit is necessary. Touch both feet (or hands) flat on the two footprints on the mat simultaneously. You should feel a tingling sensation in both hands and a slight contraction in your forearms.

Touching only one part will have no effect.

PL: Pełna treść deklaracji zgodności CE jest dostępna na stronie www.media-tech.eu
EN: The full text of the CE declaration of conformity is available at www.media-tech.eu.
HU: A CE megfelelőségi nyilatkozat teljes szövege elérhető a www.media-tech.eu oldalon.
CZ: Úplné znění prohlášení o shodě CE je k dispozici na www.media-tech.eu.
SK: Úplné znenie vyhlásenia o zhode CE je dostupné na www.media-tech.eu.
LT: Visas CE atitikties deklaracijos tekstas yra prieinamas adresu www.media-tech.eu.
RO: Textul complet al declarației de conformitate CE este disponibil pe www.media-tech.eu.
HR: Cijeli tekst CE izjave o sukladnosti dostupan je na www.media-tech.eu.
DE: Der vollständige Text der CE-Konformitätserklärung ist unter www.media-tech.eu verfügbar.
BG: Пълният текст на декларацията за съответствие CE е достъпен на www.media-tech.eu.
GR: Το πλήρες κείμενο της δήλωσης συμμόρφωσης CE είναι διαθέσιμο στη διεύθυνση www.media-tech.eu.

PRODUCED & EU IMPORTED

MEDIA-TECH
Polska Sp. z o.o.

ADDRESS: BRZEZINY 13B,
05-074 HALINÓW, POLAND

E-MAIL: SUPPORT@MEDIA-TECH.EU

Due to continuous development, product specifications and appearance are subject to change without prior notice. For technical support, please visit www.media-tech.eu.

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EN

The device is marked with a crossed-out wheeled bin symbol, in accordance with European Directive 2012/19/EC on waste electrical and electronic equipment. Products marked with this symbol should not be disposed of or thrown away with other household waste at the end of their useful life. The user is obliged to dispose of used electrical and electronic equipment by delivering it to a designated point where such hazardous waste is recycled. Collecting this type of waste in separate locations and recycling it properly contributes to the protection of natural resources. If the product contains batteries covered by Regulation (EU) 2023/1542 of the European Parliament and of the Council, it cannot be disposed of with normal household waste. Please refer to local regulations regarding the separate collection of batteries, as proper disposal will prevent potential negative effects on the environment and human health. Proper recycling of waste electrical and electronic equipment has a beneficial effect on human health and the environment. For information on where and how to dispose of waste electrical and electronic equipment in an environmentally safe manner, users should contact their local authority, waste collection point, or the point of sale where the equipment was purchased.

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EN - General safety rules

- Before first use and operation of the Media-tech device, read the user manual carefully and follow the precautions below
- The product is not intended for children aged 0-3 years
- The product may contain packaging elements (foils, blisters) that, if used incorrectly, may cause serious bodily harm or even death
- Do not immerse the product in water unless it is waterproof (marked at least IPX7)
- Do not use in a humid environment, unless it is waterproof at IPX4
- Unless otherwise marked, the product is intended for home use at room temperature
- Use only the original accessories supplied in the set. Non-original ones may cause damage or other serious events
- Before starting the device, make sure that all elements are properly installed and secured.
- For products powered by a battery or accumulator: use only batteries that meet the technical requirements of the device
- In the event of damage to the power cable or other electrical part, do not attempt to repair the device yourself.
- Connect the device only to a suitable power source with parameters consistent with the product specifications.
- Unauthorized disassembly of an electrical, electronic or battery-based product, modifications or attempted repair may result in electric shock and/or ignition and fire.
- If an electrical or battery-based product begins to emit smoke, a burning smell, or becomes abnormally hot during operation, immediately stop using it, disconnect it from the power source and put it in a safe place away from flammable materials.
- A device equipped with an O3 ozone generator emits a minimum concentration of this gas - lethal to bacteria and fungi - not dangerous to people and animals. However, we do not recommend direct inhalation directly above the device.
- A device equipped with a UV diode emits a minimum value of these rays - lethal to bacteria and fungi - not dangerous to people and animals. However, we do not recommend observing this light from a very close distance
- In the case of speakers and headphones, too high a volume setting can damage hearing
- If the product contains colored and pulsating lighting, it should be used with great caution in the vicinity of people with epilepsy
- None of the Media-tech products are medical products and may be used as part of medical therapy
- Hygiene products (irrigators, toothbrushes) should be used strictly for their intended purpose in the areas of the body designated for this purpose. Other use may result in serious bodily harm or disability
- The use of mechanical massage products should be subject to medical consultation in the case of people with cardiological, respiratory and bone diseases
- The use of EMS microcurrent massage products should be subject to medical consultation in the case of people with cardiological and skin diseases
- EMS massagers used contrary to the instructions or for too long (above the permissible duty cycles) may cause local, superficial damage to the epidermis
- Pregnant women should consult a doctor about the possibility of using all massage devices

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