medivon®

User manual

Leg massager

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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
1. Device operation

1. Check if the device is in good condition.
 2. Starting for 2 seconds to turn on the device, the product works automatically with: medium switch / high temperature operation, press a second time for 2 seconds to turn off the device.
 3. Adjust the massage training Z 1/2/3 and one of the 3 heating modes: H1/H2/H3
 4. Universal use, the use may cause traces of use on the feet, calves, legs and/or thighs, which may not be removed immediately. If you are concerned about using it, put on stockings or pants before using the massager.
- Note: Do not use the massager on clothes or pants fastened with a control zipper. This could be a product solution.
5. The product has memory: the device remembers the current operating state (air pressure, heating, vibration mode and gear) after starting the machines. This option remains in effect after the device is turned on.


Note: When using with an external device or lying down to apply the risk of falling

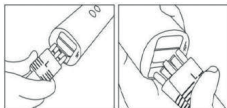How to connect the device?

Insert the connector marked with the same letters into the corresponding holes ("L" -> "L" and "R" -> "R")

Note: treat yourself so that the letter "L" is usually up and the letter "R" is down. If the device is not in this structure, air pressure cannot be generated.

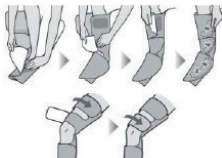Sequence of action:

1. Add the foot to the massage device, then stick the Velcro on the feet, calves (two Velcro), thighs and accessories.
2. There are 2 items on the thigh that may cause side effects.
3. It is easier to attach with Velcro after adjusting the Velcro straps on the foot, calf and thigh.

Note: Just check if you are sitting in the chair. During installation, use Velcro to prevent it from sticking to the materials.

2. Dangers and Warnings

1. Do not use replacement parts that are not approved by the manufacturer.
2. Non-professional service should not disassemble and replace parts of the device without permission.
3. If the power cord is damaged, it must be replaced by appropriate professional personnel to avoid danger.
4. Do not use or store the product in a place with high humidity or high temperature.
5. After the massage is finished, unplug the power cord.
6. Do not use the product for more than 25 minutes at a time.
7. If any side effects occur, stop using the massager immediately.
8. If the temperature is too high or the product begins to smell bad or emit a burning smell, stop using it immediately and unplug the power.
9. Do not move the product by pulling the power cord to prevent the power cord from coming loose.
10. This product is not intended for use by persons who lack physical, sensory or mental capabilities or lack of experience and knowledge in this area (including children).

If you have any of the following conditions, consult your doctor before use:

- Pregnant women, children or elderly people
- People suffering from skin infections or having skin redness and burns
- People who use a pacemaker or similar implanted medical device;
- Patients with heart disease or malignant tumors, etc.
- People suffering from osteoporosis
- People with high blood pressure or abnormal blood pressure
- People who feel discomfort when using the device

3. Cleaning and care

Before cleaning, turn off the product and unplug the power cord. Wipe the surface of the product with a dry, soft cloth. Do not clean the product and accessories with corrosive detergents such as benzene or thinners. Store in a cool, dry place when not in use. Also avoid exposure to direct sunlight or humid environments, otherwise it may cause colors to fade or damage the product.

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005

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The manufacturer declares under his sole responsibility that the Medivon Airy device complies with Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, as well as Directive 2011 of the European Parliament and of the Council /65/EU of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to the Directive of the European Parliament and 2011/65/EU as regards the list of restricted substances. The full text of the declaration of conformity is available at the internet address www.medivon.pl/deklaracje

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Shortened warranty terms

This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860.

We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for

a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: www.medivon.pl/warranty

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Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska