medivon®



Medivon Portia Arm

Medivon Portia Leg

Medivon Portia Center

User manual

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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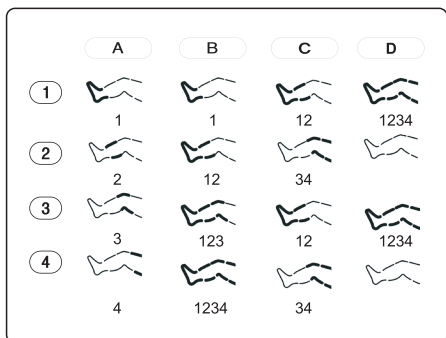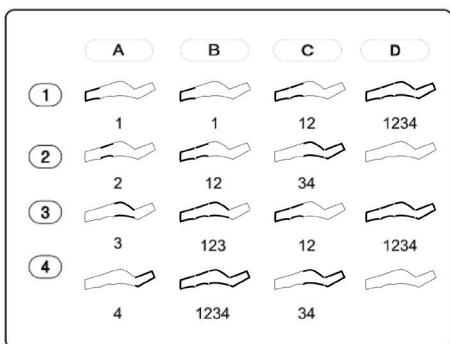Device manual

The device is ready for immediate use after purchase. Press and hold the ON/OFF button for approximately 2 seconds to switch the device on. The ON/OFF button is located on the main control unit (Medivon Portia Center).

Before starting the device, connect the pumping hoses to the control unit, paying attention to the small indentation on the connector, which indicates the correct connection direction. If only one pumping accessory is used, insert the provided plug to cover the second pumping port in order to prevent air leakage. The power input socket is located on the side of the main unit; use only the power adapter supplied with the set.

After switching on the device, repeatedly pressing the ON/OFF button cycles through the massage areas (1–2–3–4). Use the MODE button to change the operating modes (A–B–C–D). A diagram showing the mode combinations and their operating principles is provided in the illustration above. The TIMER button allows you to set the automatic shut-off time as desired. The INTENSITY button adjusts the massage strength and intensity in sequential steps. After completing the massage session, the device can be switched off by pressing and holding the ON/OFF button again.

The product is intended for home use only, at room temperature. Use the device strictly in accordance with the manufacturer's instructions. This product is not a medical device and does not replace medical advice or rehabilitation therapy.

Do not use the device if you are pregnant or if you have a pacemaker, electronic implants, or metal implants in the massage area. Consult a physician before use if you suffer from heart disease, epilepsy, thrombosis, osteoporosis, cancer, sensory disorders, or spinal conditions. Do not use the device on areas with wounds, burns, inflammation, swelling, varicose veins, or recent injuries.

Before using the product, remove any plastic protective films and all packaging components, and store them out of reach of children. Do not use the device while sleeping or while driving a vehicle. Do not operate the device near heat sources, water, or moisture. Discontinue use immediately if you experience pain, dizziness, numbness, or any other concerning symptoms. Do not exceed the recommended duration of a single massage session. Keep the device out of reach of children.

The product is not intended for use by children or by persons with reduced physical or cognitive abilities.

In case of contamination, disconnect the device from the power supply and wipe it with a damp cloth; do not immerse it in water; do not reconnect the device to the power supply until it is completely dry. Store in a dry place, away from heat sources. Do not modify the product or attempt to carry out repairs yourself.

The diagram of massage area combinations is shown on the first page of this manual. Proper operation of the accessory requires the original Medivon pumping station.

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Short Declaration of Conformity

The manufacturer declares under its sole responsibility that the Medivon Portia Center product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility, Directive 2014/35/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU as regards the list of restricted substances. The full text of the Declaration of Conformity is available at the following website: www.medivon.pl/deklaracje

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Producent: Zhejiang Xi lian Electronic Technology Co., Ltd.
No.2, Yihuan Road, Xinxin Village, Kunyang Town, Pingyang County, Wenzhou, 325499 Zhejiang, China.
E-mail: export@clavcare.com

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
Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005.

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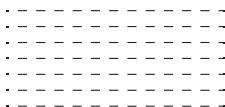Shortened warranty terms

This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860. We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements. In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them. Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: www.medivon.pl/warranty

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Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska | Marka: Medivon | E-mail: serwis@medivon.pl