



medivon.

User manual

Medivon Sooth

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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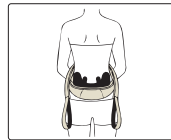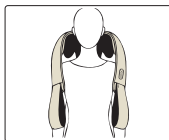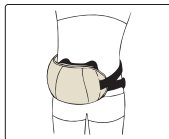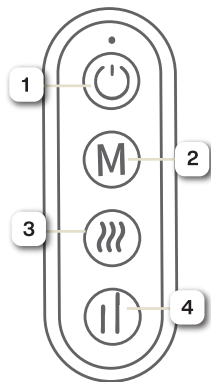Before use, it is recommended to **fully** charge the device's battery using the charger. For this purpose, **only** the charger included in the set should be used. It is recommended not to **completely** discharge the device to extend the battery **life**. It is forbidden to use the massager while charging the battery. This may cause damage or fire.

Controller operation:

1. Power button - press to turn on the massager, press again to turn off the massager.
2. Mode change button - click the button to **cyclically** switch between 2 modes: **vitality** and **relaxation**.
3. Heating button - press while the massager is working to turn the heating function on or off.
4. Intensity button - briefly click the button to switch between 2 massage intensity **levels**.

Examples of how to use the massager are shown in the **illustration** on page 1.

Safety instructions Before starting the massager, **please** read the instructions **carefully** to use the device **safely** and **efficiently**. Before first use, charge the product for about 2.5 hours to **fully** charge it. The product cannot be switched on or used **while** charging. Before using the massager, remove any **jewelry**, metal accessories, etc. to avoid discomfort during use. Wait one hour after the **last** meal before use.

Keep the massager out of the reach of children - this is not a toy! Use the massager for no more than 10 minutes per session. After this time, the device **will** automatically turn off. Take a break for a few minutes before resuming use. **Unplug** the adapter from the socket when not in use.

Do not use the device if you are pregnant, have heart problems or have any physical injuries. Do not use the device when you feel **unwell**. If you are unsure, **consult** a doctor before use.

Do not use or **leave** the massager in the bathroom or other places exposed to moisture. If the device does not work properly, **immediately** disconnect the power supply and contact the service center. If the power cord is damaged, do not attempt to use the device. The cable must be replaced by the manufacturer, its service agent or a **similarly** qualified person in order to avoid danger to **life** or health. Avoid contact with sharp edges or pointed objects that could cut or puncture the surface of the device. Store the device at room temperature, in a dry and shaded place, out of the reach of children.

Care and cleaning

Always **unplug** the device from the power supply before cleaning. Never touch the **plug** with wet hands. Wait for the device to **cool** down before **cleaning**.

Do not use **cloths** containing chemicals or alcohol or solvent-based **liquids**. Never immerse any part of the device in **liquid**.

If the device is particularly **dirty**, **follow** the instructions:

1. Dampen the cloth in water or a 3-5% mild detergent solution.
2. Wipe the dirty areas with a damp cloth.
3. Rinse the cloth in water and squeeze it out. Then wipe any detergent residue from the surface of the device (until it is dry).

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005.

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Short Declaration of Conformity

The manufacturer declares under its sole responsibility that the Medivon Sooth product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility, Directive 2014/35/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU as regards the list of restricted substances.

The full text of the Declaration of Conformity is available at the following website: www.medivon.pl/deklaracje

Importer: Armare S.A.
ul Grzybowska 87,
00-844 Warszawa, Poland
Made in China
Marka: Medivon
E-mail: serwis@medivon.pl
www.medivon.pl
SN:01265904119283635

Producent: Xiamen Healthpal Electronic Co., Ltd.
NO.148 BINGLANGXILI, SIMINGDISTRICT, XIAMEN,FUJIAN, CHINA
E-mail: vilin@xmhealthcare.com

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Shortened warranty terms

This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860. We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: <http://medivon.pl/warranty>