


medivon.

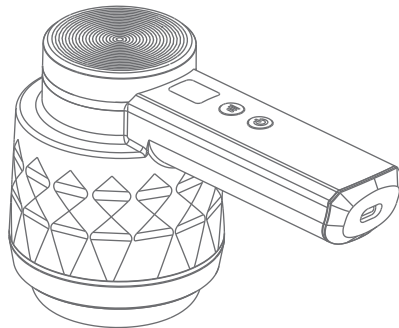Body Massager

Deep relaxation

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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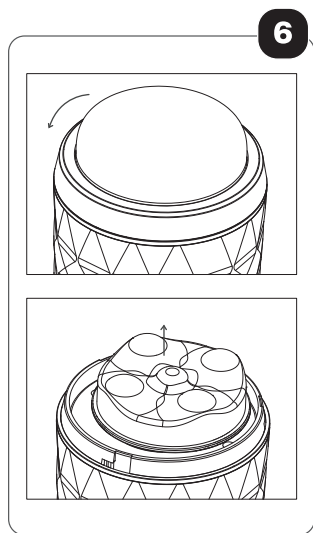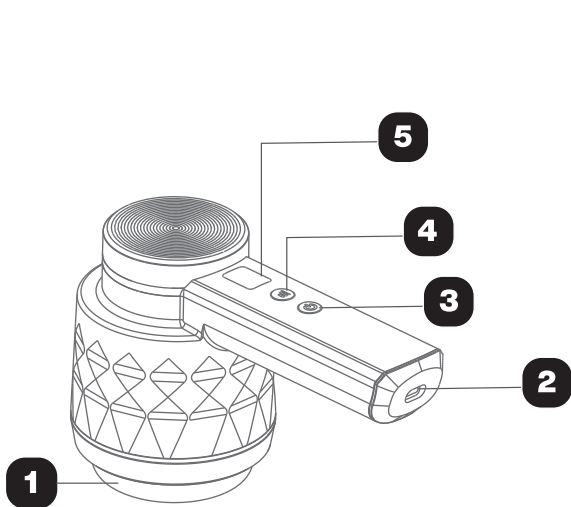1. Operation and application

Turn on the massager by pressing the power button (3) for 2 seconds; the device will start and enter working mode.

In working mode, a short press of the power button (3) allows you to cycle between working modes; There are five different modes available. A short press on the adjustment button (4) allows you to adjust the massage strength; 9 intensity levels are available.

To turn off the device, hold down the power button (3) for 2 seconds.

The massager is powered by a built-in lithium battery. If the battery is exhausted or low, please charge the device before next use.

Use the included cable for charging. Check that the power connection is correct and that the cable is properly connected to the device (2).

Voltage instability may occur during charging; Do not use the product while charging. The charging process takes approximately 4 hours. After charging, unplug the power.

Replacing massage heads:

As shown in the picture (6), turn counterclockwise to unscrew and remove the massage head mesh cover (1). After replacing the head, screw on the mesh cover clockwise. The heads can be replaced by pulling them upwards.

2. Dangers and warnings

People with medical implants, undergoing treatment, heart disease, high blood pressure, acute illness, pregnancy, osteoporosis, spine fractures, skin damage, fever above 38°C and people who are insensitive to heat should consult a doctor before use. Minors and people with intellectual disabilities may use the product only under supervision. The product is intended for use by one person, do not use it with others. After moving the product from a cold to a warm place, wait an hour before using it. Recommended use time is 10 minutes. Do not handle the product with wet hands and do not use it outdoors. Only use a power supply with the appropriate voltage. Do not smoke while using the product. Please use a qualified USB adapter for charging. Do not use in a humid environment. Do not use the product immediately after meals or simultaneously with other therapeutic devices. If there are any abnormalities during use, discontinue use immediately and consult a doctor. If the power cord is damaged, contact your distributor for replacement. When charging is complete or when not in use, unplug the plug from the outlet. Do not use the product while standing on it or damage the power cord. Do not use on an electric heating blanket or space heater. Do not disassemble the product unless you are an authorized service technician. In the event of a power outage, immediately turn off the product and unplug it from the power supply. Do not use the product while sleeping, intoxicated or feeling unwell. If there are any abnormalities (smoke, unpleasant odor, etc.), discontinue use immediately. Do not clean the product with gasoline, solvents or insecticides. In case of water spill, do not use the product. Do not use during a lightning storm. Before disposing of the product, please remove the battery. Keep these tips for safe and effective use of the product.

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005.

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This device was placed on the market after 13.08.2005.

The manufacturer hereby declares that the Medivon Lift device complies with Directive 2014/35/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65/EU of the European Parliament and of the Council of June 8, 2011 on restrictions on the use of certain hazardous substances in electrical and electronic equipment. The full content of the declaration of conformity is available at the following internet address: www.medivon.pl/deklaracje

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska


Model: Medivon Lift | Made in China | Marka: Medivon

email: serwis@medivon.pl | SN:11255904119283789

Producent: Xiamen Everwinn Electronic Technology CO. LTD

5th Floor, Building 3, No. 289, Wengjiao Road, Haicang District, Xiamen, Fujian, China

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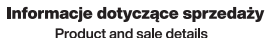This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860. We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website:
<http://medivon.pl/warranty>