
medivon.

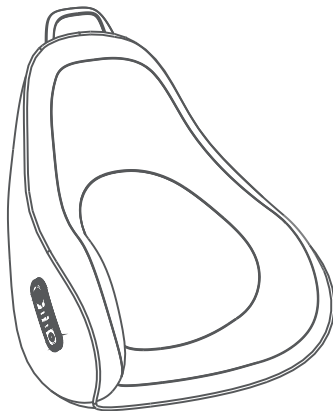Shiatsu Back Massage

Deep relaxation

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1. Operation and application

The device is wireless, but it is recommended to fully charge the battery using a USB-C cable before first use.

To turn on the device, hold down the POWER button (5) for about 2 seconds. The percussion massage mode and warm-up mode will be activated first; a short press of the same button changes the direction of the massage and turns off the percussive massage, another cyclical press of the button turns off the warm-up mode, and a long press for two seconds turns off the device.

The intensity button (3) pressed successively adjusts the intensity levels for the head massage.

While charging (4), the device stops working and the power indicator flashes green; When the battery is fully charged, the indicator is solid green. After turning on, the device works for 15 minutes by default, then turns off automatically. A long press on the power button restarts the massage. Possible massage areas (6) are: neck, abdomen, calves, waist, back, thighs. An example application is presented in the graphic.

The device has an automatic timer function, but it is recommended that each use should not exceed 30 minutes.

Prolonged disuse of the battery may shorten its life. It is recommended to charge the battery at least once or twice a month.

2. Dangers and warnings

Pregnant and menstruating women, people with serious heart disease, hemorrhagic diseases, sepsis, fever should not use this product. Children and people with intellectual or physical disabilities should use the device under supervision. People sensitive to heat must be careful due to the heating function of the device. Do not stand or put excessive pressure on the device to avoid damage. Do not use without protective cover in place and on skin with wounds or inflammation. After the massage, drink 500 ml of water or hot tea to support your metabolism. Do not use the device if the housing or cover is damaged. Do not attempt to repair or disassemble the device yourself. Avoid putting your fingers under the rotating massage elements. Only qualified personnel may perform repairs and modifications to the device. Unplug the appliance from the power supply before performing maintenance and do not use wet hands to handle the plug to avoid electric shock. Unplug the device when not in use. In the event of a power outage, immediately unplug the device. Do not use the device with any power source other than the recommended one to avoid fire or damage. When disconnecting the cable, grasp the plug, not the cable. If the cable is damaged, replacement should be carried out by qualified personnel. Do not use the device in damp or very dusty places to avoid electric shock or malfunction. Do not place the device near heat sources or in places exposed to direct sunlight. Avoid using the device during a lightning storm.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment.

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This device was placed on the market after 13.08.2005.

The manufacturer hereby declares that the Medivon Zona device complies with Directive 2014/35/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65/EU of the European Parliament and of the Council of June 8, 2011 on restrictions on the use of certain hazardous substances in electrical and electronic equipment. The full content of the declaration of conformity is available at the following internet address: www.medivon.pl/deklaracje

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska

Model: Medivon Zona