medivon.

User manual

Compression massager

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
1. Device operation

1. Check if the device is in good condition.
2. Attach the massager to your lower leg and make sure the power button is facing to the right. Then adjust to proper tightness. (Note: the massager can be used under clothes).
3. Choose the settings that are right for you.
4. Please note that the device may leave traces of use on the feet, calves, knees and/or thighs, which may not disappear immediately. If you are concerned about traces of use, put on stockings or pants before using the massager.

The product works primarily through a multi-stage air bag that is inflated and deflated multiple times in sequence, mimicking the deep tissue compression style technique used by a real massage therapist. The device improves local blood circulation, relaxes and increases muscle vitality.

The device has 6 airbag modes, 3 intensities, 3 heating levels and 4 vibration modes that can be adjusted and used individually or on both legs according to the following models:

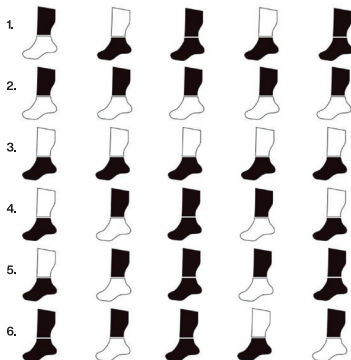2. Dangers and Warnings

1. Do not use replacement parts that are not approved by the manufacturer.
2. Non-professional service should not disassemble and replace parts of the device without permission.
3. If the power cord is damaged, it must be replaced by appropriate professional personnel to avoid danger.
4. Do not use or store the product in a place with high humidity or high temperature.
5. After the massage is finished, unplug the power cord.
6. Do not use the product for more than 15 minutes at a time.
7. If any side effects occur, stop using the massager immediately.
8. If the temperature is too high or the product begins to smell bad or emit a burning smell, stop using it immediately and unplug the power.
9. Do not move the product by pulling the power cord to prevent the power cord from coming loose.
10. This product is not intended for use by persons who lack physical, sensory or mental capabilities or lack of experience and knowledge in this area (including children).

If you have any of the following conditions, consult your doctor before use:

- Pregnant women, children or elderly people
- People suffering from skin infections or having skin redness and burns
- People who use a pacemaker or similar implanted medical device;
- Patients with heart disease or malignant tumors, etc.
- People suffering from osteoporosis
- People with high blood pressure or abnormal blood pressure
- People who feel discomfort when using the device

3. Cleaning and care

Before cleaning, turn off the product and unplug the power cord. Wipe the surface of the product with a dry, soft cloth. Do not clean the product and accessories with corrosive detergents such as benzene or thinners. Store in a cool, dry place when not in use. Also avoid exposure to direct sunlight or humid environments, otherwise it may cause colors to fade or damage the product.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment. The device was placed on the market after 13.08.2005.

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The manufacturer declares under his sole responsibility that the Medivon Misto device complies with Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, as well as Directive 2011 of the European Parliament and of the Council /65/EU of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to the Directive of the European Parliament and 2011/65/EU as regards the list of restricted substances. The full text of the declaration of conformity is available at the internet address www.medivon.pl/deklaracje

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska
Model: Medivon Misto