


medivon.

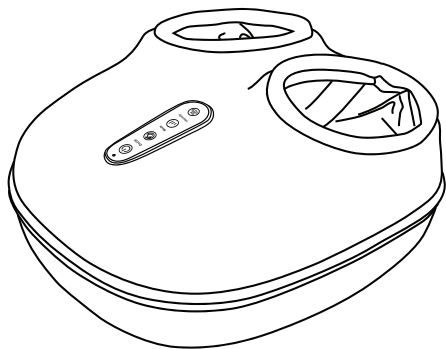User manual

Medivon Artemis

PL | EN | DE | FR | ES | PT | IT | CZ | SK
SL | EE | LT | LV | HU | RO | RU | NB | FI
EL | HR | SV | DA | NL | TR | BG | AR

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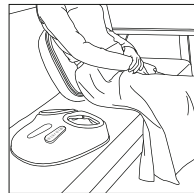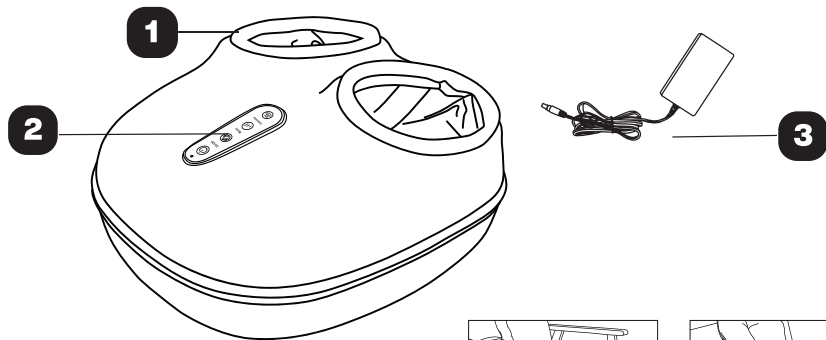1. Operation and application

1 - foot holes, 2 - remote control, 3 - adapter

The device has a rolling, kneading and heating massage function. The massager is controlled by a detachable remote control. There are three automatic massage programs with an acoustic signal. After connecting to the power supply, the device goes into standby mode.

Press the button to turn on the massager. A ready tone will appear and the device will enter mode 1 (medium massage intensity, heating on). If you want to turn off the device before the time is up, press the on/off button again; otherwise, the massager will automatically turn off after 15 minutes. After turning off the device, a long sound will appear, informing you that the massager has switched to standby.

Press the button to change one of the three available massage modes.

Press the button to change one of the three available massage intensities (low, medium, high)

Press the button to turn the heating function on or off.

There is a zipper in the center of the device. Unfold it to remove the top of the massager. The lower part of the device can then serve as a massage cushion for the back, arms and calves.

2. Dangers and warnings

People with medical implants, undergoing treatment, heart disease, high blood pressure, acute illness, pregnancy, osteoporosis, skin damage, fever above 38°C and people who are insensitive to heat should consult a doctor before use. Minors and people with intellectual disabilities may use the product only under supervision. After moving the product from a cold to a warm place, wait one hour before using it. Recommended use time is 15 minutes. Do not handle the product with wet hands and do not use it outdoors. Only use a power supply with the appropriate voltage. Do not smoke while using the product. Please use a qualified USB adapter for charging. Do not use in a humid environment. Do not use the product immediately after meals or simultaneously with other therapeutic devices. If there are any abnormalities during use, discontinue use immediately and consult a doctor. If the power cord is damaged, contact your distributor for replacement. When charging is complete or when not in use, unplug the plug from the outlet. Do not damage the power cord. Do not use on an electric heating blanket or space heater. Do not disassemble the product. In the event of a power outage, immediately turn off the product and unplug it from the power supply. Do not use the product while sleeping, intoxicated or feeling unwell. If there are any abnormalities (smoke, unpleasant odor, etc.), discontinue use immediately. Do not clean the product with gasoline, solvents or insecticides. In case of water spill, do not use the product. Do not use during a lightning storm. Before disposing of the product, please remove the battery. Keep these tips for safe and effective use of the product.

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005.

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Short Declaration of Conformity

The manufacturer declares under its sole responsibility that the Medivon Artemis product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility, Directive 2014/35/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU as regards the list of restricted substances.

The full text of the Declaration of Conformity is available at the following website: www.medivon.pl/deklaracje

Importer: Armare S.A.
ul Grzybowska 87,
00-844 Warszawa, Poland
Made in China
Marka: Medivon
E-mail: serwis@medivon.pl
www.medivon.pl
SN:01265904119283888

Producent: Xiamen High Top Electronic Technology Co. Ltd.
F2-5, No.118, Siming Park, Tong'an Industrial Zone,
Xiamen City, Fujian Province, China.
E-mail: sales6@xmhighstop.com

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Warranty Card

Pieczętka / Stamp

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Product and sale details

Seven horizontal dashed lines arranged vertically, providing space for handwritten product and sale details.

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Shortened warranty terms

This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860. We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: <http://medivon.pl/warranty>