
medisana®

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Explanation of symbols

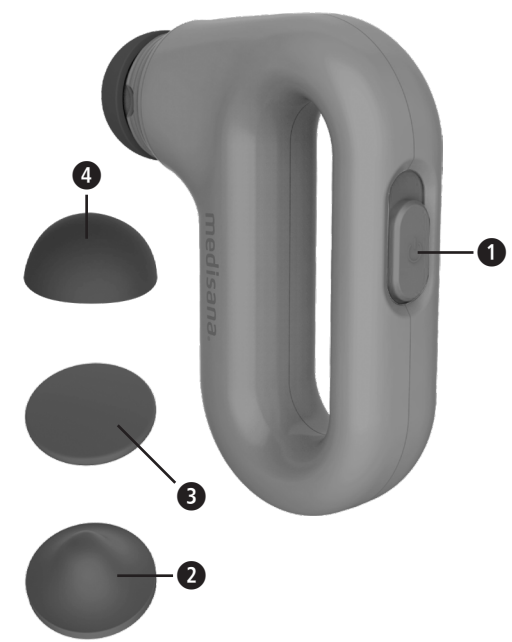IMPORTANT
Failure to follow these instructions may result in serious injury or damage to the device.

WARNING
These instructions must be followed to prevent possible injury to the user.

CAUTION
These instructions must be followed to prevent possible damage to the device.

NOTE
These instructions provide useful additional information regarding installation and operation.

LOT

LOT number

Manufacturer

Date of manufacture

Recycling symbols/codes:
These are used to provide information about the material and its proper use and recycling.

WARNING

Please ensure that the polythene packing is kept away from the reach of children! Risk of suffocation!

GB

Safety Instructions

Please read these instructions carefully before using the device, in particular the safety instructions. Keep the instructions in a safe place for future reference. If you pass the device on to third parties, be sure to include these instructions for use.

- The device is intended for private use only. The appliance must not be used for commercial or medical purposes.
- If you have any health concerns, consult your doctor before use.
- Use the unit indoors only! Do not use the device in damp rooms (e.g. when bathing or showering).
- Use the device only for its intended purpose according to the instructions for use. In case of misuse, the warranty becomes void.
- If you feel pain or find the massage uncomfortable, stop the application and consult your doctor.
- You must not use the device if you suffer from one or more of the following illnesses or complaints: Circulatory disorders, varicose veins, open wounds, bruises, skin tears, phlebitis or thrombosis.
- Do not use the appliance if you are pregnant. Special care is necessary when the appliance is used on, by or near children, sick and helpless persons.
- Do not use the device near the eyes or other sensitive parts of the body.
- Do not use this device to support or replace medical applications. Chronic conditions and symptoms could worsen.
- Do not use the device while you are sleeping or in bed.
- Do not use the device before falling asleep. The massage has a stimulating effect.
- Do not use the device while driving a vehicle or operating machinery.
- This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the product in a safe way by a person responsible for their safety.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children must not play with the appliance.
- Cleaning and user maintenance must not be carried out by children without supervision.
- Do not allow the unit to come into contact with pointed or sharp objects.
- Do not cover the unit when it is switched on. Never use it under blankets or pillows. There is a risk of fire, electric shock and injury.
- Do not immerse the unit in water or other liquids.
- If liquid does enter the unit, switch it off immediately and remove the batteries. Wait until the unit is completely dry.
- Always switch off the appliance with the power switch ❶ after use.
- Do not use the appliance if any damage is visible, if it is not working properly, if it has been dropped or immersed in water. To avoid hazards, send the unit to the service centre for repair.
- Do not insert any objects into the openings of the unit.
- You may only clean the unit. In case of malfunctions, please do not repair the unit yourself, as this will automatically void all warranty claims. Contact your specialist dealer and have the repairs carried out by the authorised service.
- Only use accessories intended for the unit and supplied by the manufacturer.

SAFETY NOTES FOR BATTERIES

- Do not disassemble batteries!
- Clean the battery and device contacts if necessary before putting in the batteries!
- Remove discharged batteries from the device immediately!
- Increased risk of leakage, avoid contact with skin, eyes and mucous membranes! If battery acid comes in contact with any of this parts, rinse the affected area with copious amounts of fresh water and seek medical attention immediately!
- If a battery has been swallowed seek medical attention immediately!
- Replace all of the batteries simultaneously!
- Only replace with batteries of the same type, never use different types of batteries together or used batteries with new ones!
- Insert the batteries correctly, observing the polarity!
- Keep the battery compartment well sealed!
- Remove the batteries from the device if it is not going to be used for an extended period!
- Keep batteries out of children's reach!
- Do not attempt to recharge these batteries! There is a danger of explosion!
- Do not short circuit! There is a danger of explosion!
- Do not throw into a fire! There is a danger of explosion!
- Keep unused batteries in their packaging away from metal objects in order to prevent short circuiting!
- Do not throw used batteries into the household refuse; put them in a hazardous waste container or take them to a battery collection point, at the shop where they were purchased!

Scope of delivery

Please first check whether the scope of delivery is complete and that the item is not damaged in any way. If in doubt, do not operate the unit and contact your dealer or service centre.

The scope of delivery includes:

- 1 medisana Mini Massager **HM 330**
- 3 Massage heads
- 1 Quick Reference Guide
- 2 batteries (type AAA, LR03) 1.5V

Packaging is recyclable or can be returned to the raw material cycle. Please dispose of packaging material that is no longer required properly. If you notice any transport damage when unpacking, please contact your dealer immediately.

Device and controls

- ❶ Power switch
- ❷ Massage head "Trigger point"
- ❸ Massage head "Flat"
- ❹ Massage head "Round"

Insert / change batteries

Before you can use your device, you must insert the enclosed batteries (2 x 1.5V batteries, type AAA/LR03). To do this, open the battery compartment cover and insert the batteries with the correct polarity (as marked in the battery compartment). Close the battery compartment again.

Replace the batteries if the unit no longer works after being switched on or if the massage suddenly stops..

Fitting / changing the massage head

3 different massage heads are included - you can exchange them as you like depending on your personal preference. Please check the connection of the components (appliance and massage head attachment) before the appliance is put into operation.

Massage attachment "Trigger point" ❷: Suitable for massaging deep tissue and trigger points, e.g. on the neck, palms or soles of the feet.

Massage attachment "Flat" ❸: Massage of the entire body musculature.

Massage attachment "Round" ❹: Particularly suitable for chest, buttock, leg and back muscles.

Basic operation

- Slide the power switch ❶ up one level to switch on the appliance. The massage starts at low level. To perform the massage at high level, slide power switch ❶ up to second level.
- Place the massage attachment with light pressure on the area to be massaged and slowly move the appliance over the target area. Massage each area only for a short time and then switch to another area or switch off the appliance. NEVER massage the same area for more than 3 minutes.
- Slide power switch ❶ all the way down to switch off the unit.

Cleaning and care

- Before cleaning the item, remove the batteries and let the unit cool down for at least 15 minutes.
- Clean the unit with a soft cloth slightly moistened with water. Never use harsh cleaning agents, petrol, thinner or alcohol.
- Then dry the unit with a dry, clean cloth.
- Do not expose the unit to direct sunlight, protect it from dirt and moisture.
- Make sure the unit is completely dry before using it again.
- Store the unit in a clean and dry place.

Technical specifications

Name and model:	medisana Mini Massager HM 330
Power supply:	3V=== , 2 x 1.5V batteries (type AAA, LR03)
Operating conditions:	Use only in dry indoor areas
Storage conditions:	clean, cool and dry
Dimensions L x W x H:	approx. 133 x 111.5 x 44 mm
Weight:	approx. 137 g (without batteries)
Article number:	88261
EAN number:	40 15588 88261 6

Disposal

This product must not be disposed of together with domestic waste. All users are obliged to hand in all electrical or electronic devices, regardless of whether or not they contain toxic substances, at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority or your dealer for information about disposal.

In accordance with our policy of continual product improvement, we reserve the right to make technical and visual changes without notice.

The current version of this instruction manual can be found under www.medisana.com

Warranty and repair terms

Your statutory warranty rights are not restricted by our guarantee below. Please contact your dealer or the service centre in case of a claim under the warranty. If you have to return the unit, please enclose a copy of your receipt and state what the defect is. The following warranty terms apply:

- The warranty period for **medisana** products is 3 years from date of purchase. In case of a warranty claim, the date of purchase has to be proven by means of the sales receipt or invoice.
- Defects in material or workmanship will be removed free of charge within the warranty period.
- Repairs under warranty do not extend the warranty period either for the unit or for the replacement parts.
- The following is excluded under the warranty:
 - All damage which has arisen due to improper treatment, e.g. non-observance of the user instructions.
 - All damage which is due to repairs or tampering by the customer or unauthorised third parties.
 - Damage which has arisen during transport from the manufacturer to the consumer or during transport to the service centre.
 - Accessories which are subject to normal wear and tear.
- Liability for direct or indirect consequential losses caused by the unit are excluded even if the damage to the unit is accepted as a warranty claim.

Medisana, Carl-Schurz-Str. 2, 41460 NEUSS, GERMANY

Information about service can be found here: <https://www.medisana.com/servicepartners>

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