
PL EN DE FR ES PT IT CZ SK SL EE LT LV
NB FI HU RO AR RU EL HR SV DA NL TR BG

BEAUTIFLY

ABS Tens Sculptor

FitRecovery

series

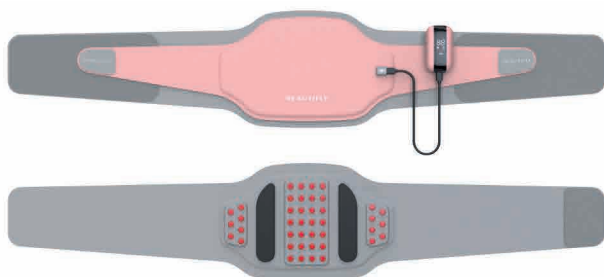ABS Tens Sculptor Belt

User manual

/beautifuly.global

/beautifuly_global


/@beautifuly.global


/@beautifuly.global


f1

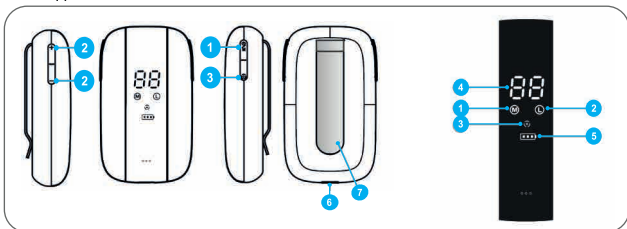f2

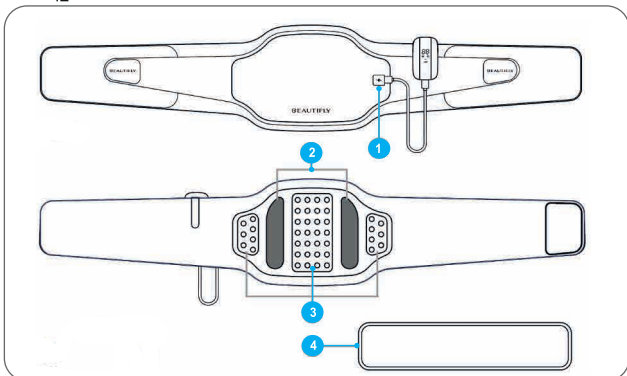f3

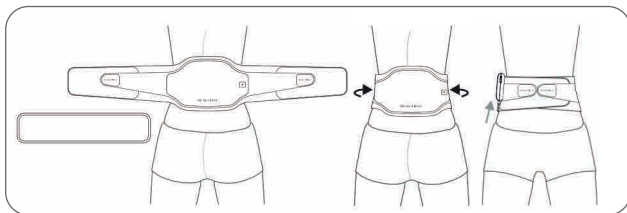f4

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- www.beautify.eu
- www.youtube.com/beautifyeu

Frequency:	60-100Hz
Battery Capacity:	2500mAh
Input current:	4.1mA
Charging time:	240 minutes
Intensity levels:	15
Massage modes:	7
Strap Size:	1100*215mm
Extension belt:	about 400*100mm

Security:

Before using the product, read this user manual and follow the recommendations contained therein when using the device.

Device description [f1]:

1. "M" - enable/disable and select EMS mode.
2. "L" "+/-" increase/decrease the intensity level.
3. Enable/disable RLT mode.
4. Display.
5. Battery capacity.
6. Type-C charging port.
7. Hanger.

[f2]

1. Type C input.
2. EMS electrodes.
3. Red and infrared light therapy.
4. Extension strap.

How to use?

1. Before first use or after a long period of non-use when the battery is discharged, please charge the device according to the charging instructions to ensure proper operation.
2. Before use, it is recommended to cleanse the skin and apply a non-greasy emulsion or water to moisturize it.
3. Depending on your waist circumference, decide whether you need to use an extension belt [f2-4]. Place the belt directly on your skin and fit it snugly, positioning the electrodes centrally so that you can comfortably select options on the display screen [f3]. The belt can be used on the back and stomach.
4. Connect the belt with the remote control using the included USB Type-C cable.
5. Long press [f1-1] to turn on the device.
6. Select settings according to your preferences:
 - Hold the [f1-1] button for about 2 seconds to turn on: EMS (Electromuscular Stimulation), which causes muscle contractions, bringing relaxation and rest, stimulates the release of endorphins, improves circulation, supports pain relief and healing. Working time: 15 minutes. A short press of the [f1-1] button allows you to cycle through the available 7 modes.
 - Short press the "+" and "-" buttons [f1-2] respectively increase or decrease the intensity level "L" from 1 to 15.
 - Hold the [f1-3] button for approximately 2 seconds to activate: RLT (Red and Infrared Light Therapy), which promotes cellular regeneration, relieves pain and reduces inflammation by penetrating deep into the tissues. Working time: 15 minutes.
7. Long press [f1-1] to turn off the device.

NOTE: If the connection between the device and the belt is incorrect or the EMS electrodes are not in contact with the skin, the EMS intensity level will automatically return to level 1. The device will beep twice and automatically turn off after approximately 30 seconds. As this is a product intended for contact with the skin, we suggest cleaning the belt thoroughly before using it by another person.

Landing:

The product does not have a dedicated adapter, so use only a 5 V/2 A power supply. Remember to always turn off the device before charging. Insert the USB Type-C cable into the charging port. While charging, the white light will flash, and when the product is charged, all 3 lights will be on constantly.

Note: No lights and 3 beeps mean the belt is unloaded. To ensure long-term performance of the device, we recommend regular charging every two months.

Care and storage:

Before cleaning, turn off the product and unplug the power cord. Wipe the surface of the product with a soft, damp cloth and allow it to dry naturally. Do not clean the product and accessories with corrosive detergents. Store in a cool, dry place, away from children and pets. Avoid exposure to direct sunlight or humid environment, otherwise it may cause color fading or damage.

Warnings and dangers:

If you feel any discomfort, stop using the device immediately. This is not a medical device and is not intended for medical purposes. Do not use it if you have: a pacemaker, defibrillator, artificial lung/heart, migraines, malignancy or suffer from skin diseases or infections such as neurodermatitis, purpura and acne. The product is not recommended for use by pregnant women, people recovering after surgery, suffering from sensory nerve loss, cognitive impairment or having an open wound on the skin surface.

The product is not a toy. The device is not recommended for use by children and adolescents under 18 years of age. The product is not intended for people who do not have physical, sensory or mental abilities or do not have experience and knowledge in this area (including children).

The device is not waterproof and should not be washed or cleaned under running water.

When using, do not press on the internal organs to avoid discomfort and/or injury. Please do not use it for a long time to avoid low temperature burns when using the heating function.

If any side effects occur, stop using the device immediately.

More information about the company and products [f4]:

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[EN] Rules for handling electrical and electronic equipment: Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution.

Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment.

Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m² devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment.

The product was introduced to the market after 13 August 2005.

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[EN] Summary declaration of conformity: The manufacturer hereby declares under his sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of legislation, Member States in the field of electromagnetic compatibility and Directive 2011/65/EC of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment and Commission Delegated Directive (EU) 2015/863 of 2011. On 31 March 2015, amends Annex II to Directive 2011/65/EU of the European Parliament and of the Council on the list of restricted substances.

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Model: ABS Tens Sculptor FitRecovery

Made in China.

SN: 03255903003702375

Marka: Beautifly

e-mail: serwis@beautifly.eu