

medivon.

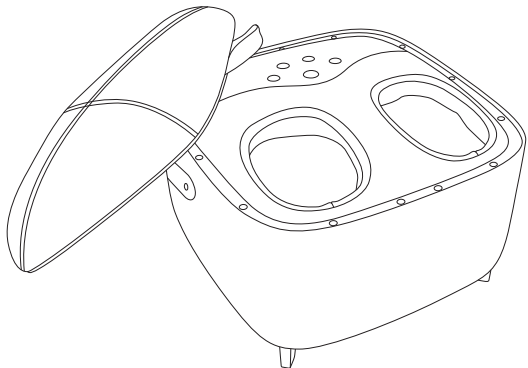User manual

PL | EN | DE | FR | ES | PT | IT | CZ | SK |
SL | EE | LT | LV | NB | FI | HU | RO | AR
RU | EL | HR | SV | DA | NL | TR | BG

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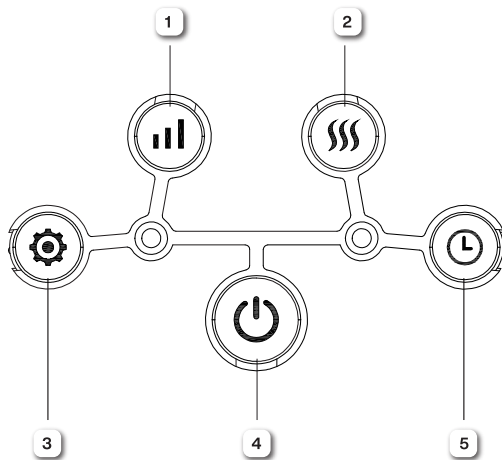Configuration and operation

Place the product on a stable surface and connect it to a power source using the included power adapter. The device is ready to use - it can be used as a piece of furniture with storage for small items or, after unfolding, as a classic foot massager.

To start the device after unfolding, put your feet into the massager and then press the POWER button (4). The default demo mode will start. Using the Shiatsu massage button (3), we can cycle between massage modes. The next button (2) allows you to turn the foot heating function on or off. Cyclic pressing of the button (1) controls the intensity of pressure and air massage.

The default massage time (5) is set to 15 minutes - blue. After pressing the button, you can reduce the duration of the massage to 10 minutes - green color, or to 5 minutes - red color.

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Waste electrical and electronic equipment as well as batteries and accumulators cannot be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage may result in damage to health or environmental pollution. Proper household handling of waste electrical and electronic equipment and batteries and accumulators positively contributes to the reuse and recovery, including recycling, of waste equipment. Waste electrical and electronic equipment as well as batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and performed the same functions as the equipment sold. When delivering equipment intended for households to the buyer, the distributor is obliged to collect used equipment from households free of charge at the place of delivery of the equipment, provided that the used equipment is of the same type and performed the same functions as the delivered equipment. A distributor running a retail unit with a sales area within the meaning of Art. 2 point 19 of the Act of March 27, 2003 on spatial planning and development of an area of at least 400 m² devoted to the sale of equipment intended for households, is obliged to accept waste equipment from households, free of charge, in this unit or in its immediate vicinity, none of the external dimensions exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market on August 13, 2005

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The manufacturer hereby declares that the Medivon Atlantis device complies with Directive 2014/35 / EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to the making available on the market of electrical equipment designed for use within certain voltage limits, Directive 2014/30 / EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65 / EU of the European Parliament and of the Council of 8 June 2011 on restriction of the use of certain hazardous substances in electrical and electronic equipment. The full text of the declaration of conformity is available at the internet address: www.medivon.pl/deklaracje

Importer: Armare S.A. | ul Grzybowska 87, 00-844 Warszawa, Polska | Model: Medivon Atlantis | SN:12255904119284571 | e-mail: serwis@medivon.pl | Made in China.
Producent: Shenzhen Forrest Health Technology Co., Ltd | 201, Building 1, No. 30, QinFu Road, LiuYue South Community, HengGang Street, LongGang District, Shenzhen, China
email: jason@forrestmassager.com

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Shortened warranty terms

This Guarantee is provided by Armare S.A., a company governed by Polish law, with its registered office in Warsaw, ul. Grzybowska 87, 00-844, Poland. You can find us in the official Polish register of companies under the KRS number: 0000730590, NIP: 5213833614 and REGON: 380721860. We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts: 24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: <http://medivon.pl/warranty>