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
PL

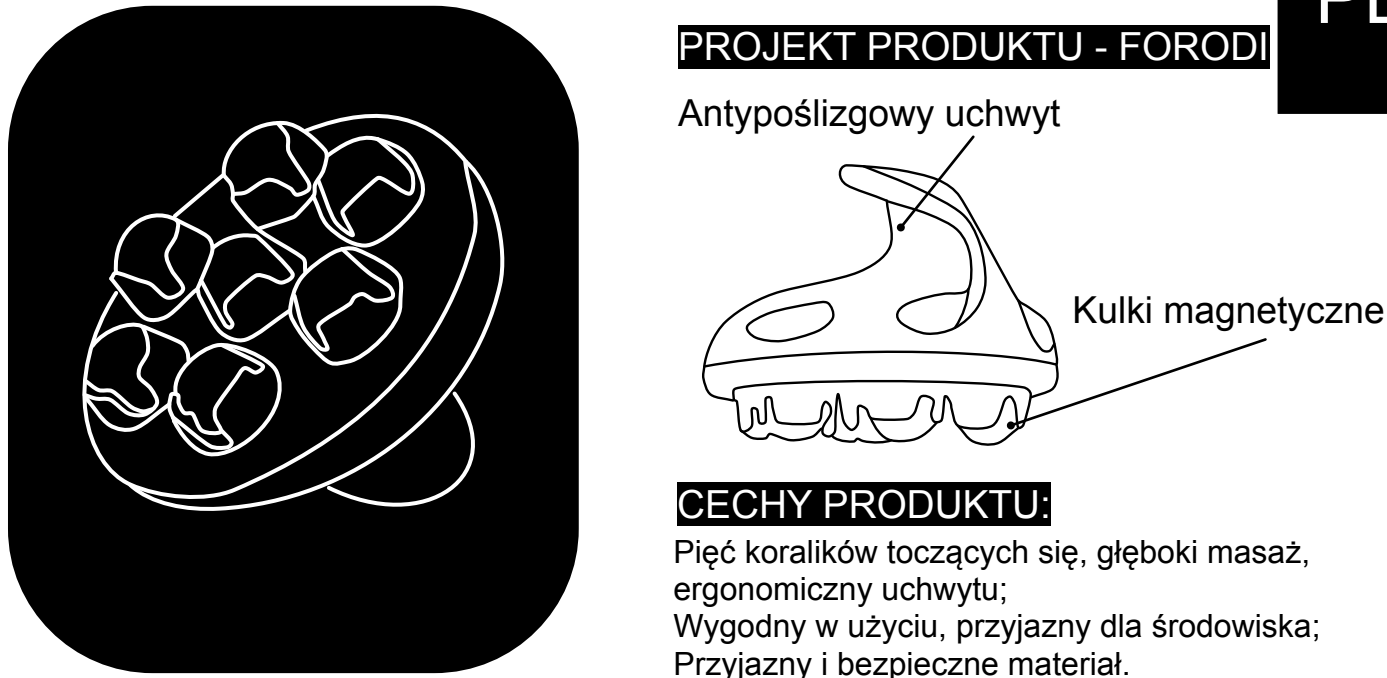PROJEKT PRODUKTU - FORODI

CECHY PRODUKTU:

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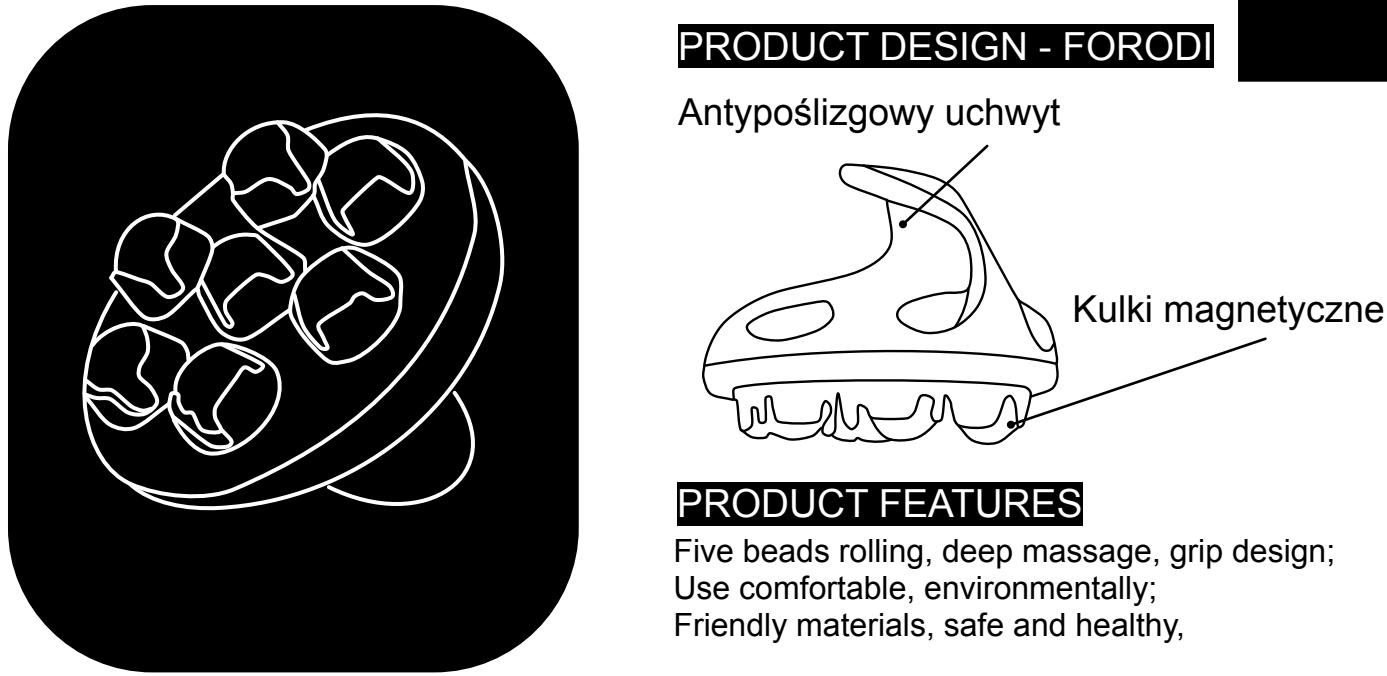
PRODUCT DESIGN - FORODI

PRODUCT FEATURES

Five beads rolling, deep massage, grip design; Use comfortable, environmentally; Friendly materials, safe and healthy,

- NOTE:**
- Do not use in pregnant woman and those with skin inflammation.
 - Please place this item out of reach for children.

CLEANING:

Method 1: Magnetic ball cleaning:
Find a thick soft towel, mix mild and neutral detergent first, soak the towel, and then put the massager on the towel to roll back and forward to clean the ball.
Method 2: Deep cleaning:
1. First use a wooden rod or plastic rod to open the steel ball. Do not use an iron rod or a knife to avoid damaging the buckle
2. After removing the five steel balls, soak them in neutral detergent for a while, and then clean the stains in steel balls and clips with and unused toothbrush.
Note: It is normal phenomenon that there will be black friction marks in buckle after use, and there will be metal and plastic friction!

Warranty:

The manufacturer guarantees its manufactured products for a period of 24 months from the date of purchase. Accessories included with the original product are warranted for a period of 12 months from the date of purchase.
For the efficient operation of products and complaints, we recommend.
Study and follow the operating instructions carefully.
1. Keep the manual, warranty card, proof of purchase.
2. In case of problems with the device, contact the Manufacturer.
3. Do not make repairs or modifications yourself.
4. Use the product according to its intended use.
5. During the warranty period, all repairs are performed free of charge.

Limitations:

The warranty will not be recognized, in case of mechanical damage due to inattention or fortuitous events.
If the device is used contrary to its intended use and specifications.

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
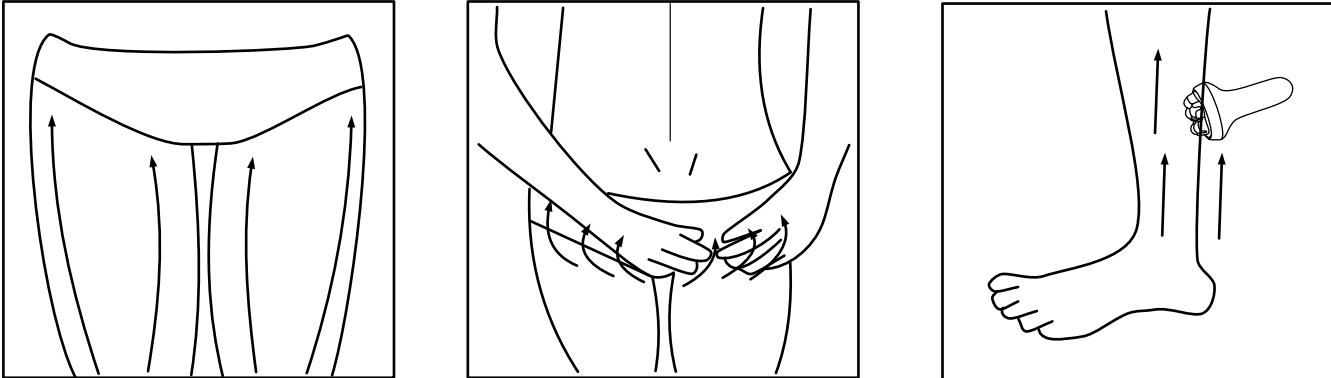
USAGE METHOD

Usage: Apply 2-3ml massage oil to corresponding body parts, and the massage with meridian brush for about 10 minutes.

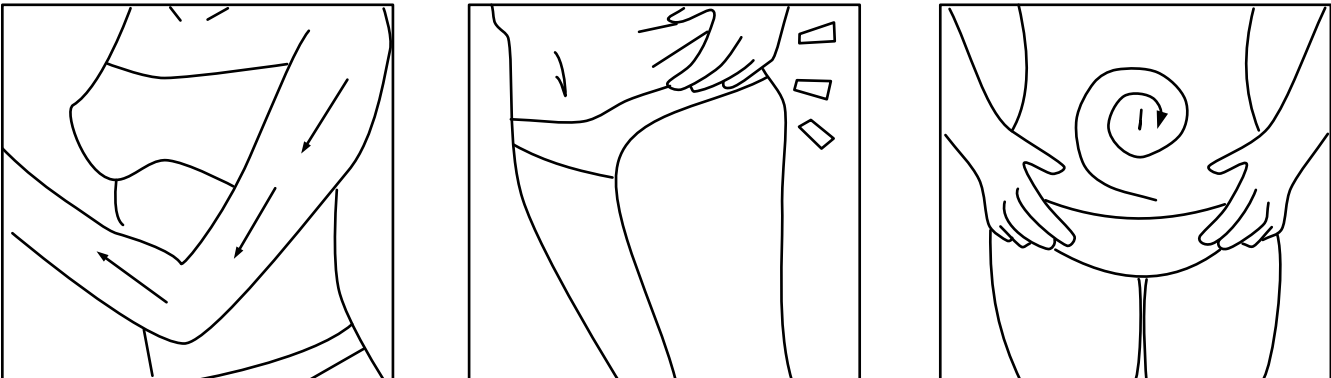•Arm Massage


Start at the upper end of the arm,Massage from top down, the inside, Finally massage lateral.

•Lumbar Massage

Pinch the excrescences on both side of the waist. After appropriate stimulation, slides from bottom top with product.

•Abdominal Massage

Massage was rotated in the abdomen in clockwise direction.

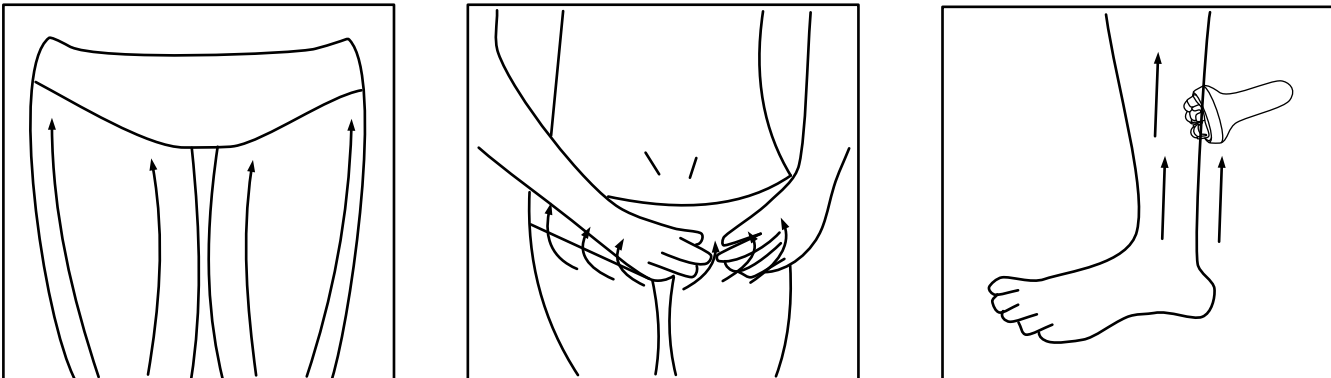•Thigh Massage

Helical fashion massage oil is used and massage is given from knee down to up with product.

•Buttock Massage

Starting at the hip base, massage the hips from the bottom up, especially where it is easy to hoard fat.

•Calf Massage

From the ankle, massage the calf from down to up from muscularity.

The following persons should not use:

- Patients with microvascular dilation caused by long-term use of steroids or impaired liver skills;
- People with allergies, atopic dermatitis, skin diseases, sensitive skin or sensitive skin;
- Those who have problems such as pigmentation and inflammation due to chemical fibers, or who have or have had dermatitis due to cosmetics;
- Those who have recently received chemical peels, and those who have had skin problems after using beauty equipment;
- * Do not use other parts except designated parts.
- * Those who are receiving doctor's treatment should consult a doctor before using.

FORODI Brand

MODEL: F5-Roller
Consultation hours: 08:00 ~ 16:00 (8 AM ~ 4 PM)
CUSTOMER SERVICE: +48 453 271 100
Made in China

PRODUCENT: FORODI by
VOBE Sp. z o.o., Rudna Mała 167A,
36-054 Mrowla, Poland
NIP: 5252799270
KRS 0000800189