IRYGATOR DO JAMY USTNEJ

CLEO

ORAL IRRIGATOR

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PL / ENG

SPIS TREŚCI

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RoHS

INDEX

I.	Warnings	20
II.	Technical data	22
III.	Product and accessories description	23
IV.	How to use	26
V.	Usage tips	29
VI.	Clean and maintenance	30
VII.	Troubleshooting	31
VIII.	Warranty conditions	32


I. WARNINGS

1. For the first-time usage, please select Child mode, in order to avoid the gums bleeding due to the lack of adaption to usage.
2. For patients with periodontal disease or those who have had oral surgery in the last 2 months, please consult your doctor before use.
3. Use this product only as indicated in these instructions or as recommended by your dental professional.
4. It must only be supplied at safety extra low voltage.
5. Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Please contact your local distributor or manufacturer's service department.
6. If product is opened or disassembled for any reason, warranty is void.
7. Keep charging unit and cord away from heated surfaces.
8. Do not use this product outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
9. Oral irrigators shall be used only with cold or warm water, or specific solutions as defined in the manufacturer's instructions.
10. Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
11. Do not direct water into the nose or ear. The potentially deadly amoeba *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into these areas.
12. Magnetic field: Pacemakers and similar implanted medical devices may be affected by magnet.
13. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
14. Children shall not play with the appliance.
15. Closely instruct and supervise children and individuals with special needs in the proper use of this product.
16. Save this user manual.


II. TECHNICAL DATA

Model	AOW09-W / AOW09-B
Voltage	li-ion battery 1800 mAh
Charging time	4-8 h
Working tmie	20+ days
Noise	≤ 72 dB
Water pressure	2-10 bars
Flow rate	200 ml/min
Tank capacity	350 ml
Waterproof	IPX7
Dimensions	82 x 82 x 318 mm
Net weight	368 g

III. PRODUCT AND ACCESSORIES DESCRIPTION

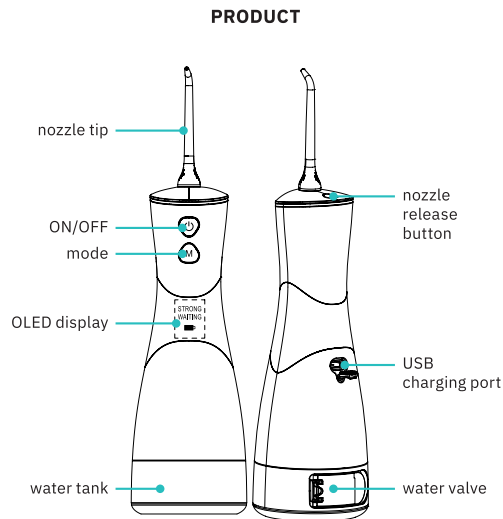ACCESSORIES

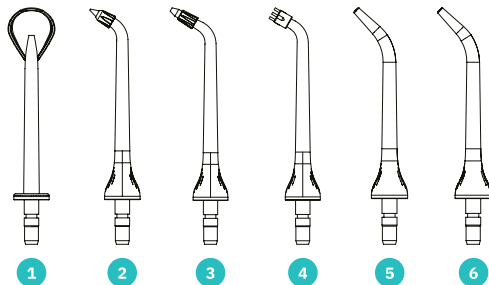1. Tongue scraping

Effectively reduces tongue coatings by removing bacteria and sulfur compounds responsible for bad breath.

Regular use improves oral hygiene and provides a long-lasting feeling of freshness.

2. Periodontal pocket nozzle

Enables precise irrigation of periodontal pockets, supporting the treatment of gum inflammation.

It is dedicated to hard-to-reach areas and enhances daily care.

3. Orthodontic nozzle

Facilitates cleaning of orthodontic appliances, effectively removing food debris and dental plaque.

Thanks to its functionality, it improves hygiene for individuals wearing braces.

4. Dental Plaque nozzle


Ensures effective cleaning around these structures and helps maintain hygiene in hard-to-reach prosthetic areas.

5. Standard nozzle

Designed for daily cleaning of interdental spaces and gum pockets, effectively removes dental plaque.

It is versatile and suitable for most individuals.

+ USB charging cable

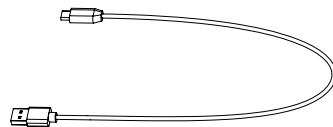IV. HOW TO USE

1. PREPARE BEFORE STARTING

Very important! For the first time to use, please fully charge the device (during charging, the product **cannot** be used).



2. INSERT NOZZLE

Insert the nozzle into the hole until you hear the *click* sound.

3. POUR WATER INTO THE WATER TANK

Ensure that during pouring water into the water tank, the device is switched off.

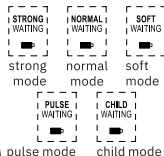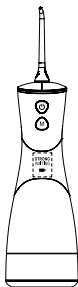
4. SELECT THE MODE

Press the **ON/OFF** button to power on, then select your suitable mode.

Warning: we recommend to select child mode for start to avoid gums bleeding.


5. POWER ON

Put the nozzle into your mouth and press the **ON/OFF** button two times.

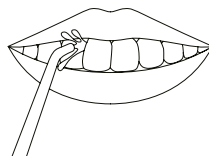6. CHANGE NOZZLE


Power off and remove the tip, press the nozzle release button, and then pull up the nozzle.

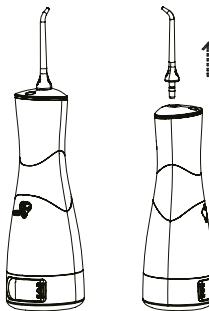- **ON/OFF** button functions:
 - first short press is to start up,
 - second short press is to the product operation,
 - third short press is to power off the product.
- **Mode button** is the working mode selection button. Short press to switch the 5 working modes:
STRONG-NORMAL-SOFT-PULSE-CHILD
- **Product operation description:**
 1. Press **ON/OFF** Button to start up, and the **WELCOME** screen will be displayed.
 2. After that, the product enters the working mode selection screen, short press **Mode button** to select the demand working mode.
 3. After selecting the working mode, short press **ON/OFF** Button again and the product will work.
 4. The product will automatically power off after 2 minutes of work. If you need the product to shut-down during work, just press **ON/OFF** Button again to enter **BYE!** screen
- The product can be charged with adapters of SV1A or above with TYPE-C USB cable.

V. USAGE TIPS


Tilt down to place the nozzle in the mouth and point the nozzle at the tooth, then turn on the device.

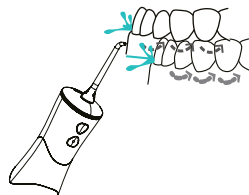Make water jet at 90° to the gum, and close your lips slightly to avoid water spray.

For best results, start with the molars (back-teeth) and clean towards the front teeth.


Slide the tip along the gum line and keep flossing a while at inter dental place, till everywhere are cleaned.

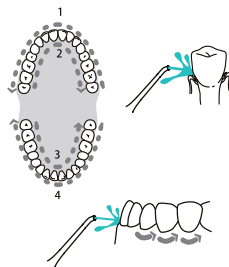VI. CLEAN AND MAINTENANCE

WATER TANK

1. Turn off the device, and open the back tank lid with the injection hole facing down.
2. Turn on the device after clearing the residual liquid in the tank, to drain the residual liquid in the device until there is no more water to drain.
3. Turn off the device.

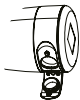Remove and clean water tank if necessary:

1. Keep your left hand to hold the product body, and your right hand to hold the water tank and take off it.
2. You can use brush to clean water tank if necessary.
(**Note:** brushes are not included)

BATTERY MAINTENANCE

- For the first use, the device must be fully charged approximately 6 hours.
- If the device is stored for a long time, please fully charge it in advance.
- In order to extend battery lifetime, please recharge it in time when power is insufficient.

VII. TROUBLESHOOTING

Problems	Reasons	Solutions
Leakage of water between jet tip and screw part	The seal ring on the fixed socket for the nozzle is damaged	Contact the after-sales
	The nozzle does not install correctly	Re-install correctly
Water tank leaks water	Seal ring is damaged	Replace the seal ring, call after-sales
Inadequate pressure	Nozzles damaged	Replace with new nozzle
	Low voltage	Charging
	Low battery	Charging
Water is not coming out	Nozzle blocking	Suck the Jet tip by mouth
	Water pipe blocking	Inject the water from the nozzle Shake after fill the water
Short working time after full charging	Insufficient charging time	Charging
	Circuit fault	Contact after-sales

VIII. WARRANTY CONDITIONS

1. The warranty is granted for a period of 24 months from the date of sale. Defects discovered during this time will be removed free of charge in Poland by the service facility within 14 business days from the date of delivery to the service facility. The condition for granting the warranty is that the user read the user manual and use the equipment in accordance with this manual.

2. Warranty repair is understood as the performance of specialized activities by the service facility, appropriate to remove the defect covered by the warranty.

3. If you find a malfunctioning device, the fastest and most convenient way will be to fill out the online form available at: webber.com.pl/reklamacja.html

4. If any incorrect operation of the device is found, please contact the Manufacturer's Service Department: via the website: www.webber.com.pl/reklamacje or by phone: 61 1020 4922

5. The complaining party has the right to replace the equipment with a new one if:

- a) during the warranty period, the service facility will make 4 repairs and the equipment will show defects that prevent it from being used as intended,
- b) the authorized service facility states in writing that it is impossible to remove the defect,
- c) the repair is not completed within the deadline specified in point. 1 or on another date agreed in writing with the complaining party.

6. The costs of repair or replacement, including in particular the costs of postage, transportation, labor and materials, are borne by the service facility.

7. The warranty does not exclude, limit or suspend the buyer's rights arising from the provisions of the Act of May 30, 2014. on consumer rights and other generally applicable legal provisions regulating the sale and liability of the manufacturer/seller for the goods.

A faulty device with a detailed description of the fault should be sent to the service address: **ul. Grunwaldzka 229, 60-179 Poznań, PL; tel. 61 1020 492**

Disposal of used electrical and electronic equipment (applies to European Union countries and other European countries with separate waste collection systems).

This symbol on the product or its packaging means that the product should not be treated as household waste. It should be handed over to the appropriate point that collects and recycles electrical and electronic equipment. Correct disposal of the product will prevent potential negative consequences for the environment and human health that could otherwise be caused by inappropriate disposal of the product.

Recycling materials helps preserve natural resources. For detailed information about recycling this product, please contact your local authority, waste disposal service or the store where you purchased the product.

Manufacturer:

METREX
ul. Grunwaldzka 229
60-179 Poznań, PL

www.webber.com.pl

RoHS

Webber 
med

METREX

ul. Grunwaldzka 229
60-179 Poznań, PL

www.webber.com.pl