concept

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ACKNOWLEDGEMENTS

Thank you for purchasing a Concept product. We wish you much joy with your new appliance every day you use it.

Please read the instruction manual carefully before initial use. To refer to this manual any time you need to, we recommend you to keep it in a safe place. And please pass it to any future owner of the appliance.

TECHNICAL SPECIFICATIONS

Voltage	5 V DC
Input	1,5 W
Functions	Clean, Soft, Massage
Oscillation	35 000 strokes/min
Charging	USB (DC) charging cable
Charging time	4 hrs.
Battery	Li-Ion 600 mAh
Water resistance	IPX7

IMPORTANT SAFETY PRECAUTIONS

- Do not use the device in any other way than as described in this manual.
- If you have had dental surgery in the last two months, consult your dentist about using the device.
- During the first days of using the device, you may notice a slight bleeding from the gums. If the gums bleed for more than 14 days, stop using the device and consult your dentist.
- Do not use brush heads other than those recommended by the manufacturer.
- If you have a pacemaker or other device installed, consult your doctor about using this toothbrush.
- The device is intended only for cleaning teeth, gums and tongue. Do not use it for other purposes.
- If you experience discomfort or pain, stop using the device and seek medical attention.
- Remove all the covering and marketing materials from the device before the first use.
- Make sure the mains voltage corresponds to the values on the rating plate of the device.

- This product is not a toy. Children should be supervised at all times when the appliance is switched on to ensure that they are using the toothbrush correctly and that they do not play with the toothbrush.
- If you put the device down during use, always switch it off.
- Before connecting or disconnecting the device from the mains socket, make sure that the device is switched off.
- When disconnecting the device from the wall outlet, never pull the power cord, grasp the plug end and disconnect it by removing it from the receptacle.
- Do not use the device in the bathtub.
- Do not leave the device in a place from which it can fall into a sink or a bathtub.
- Do not share the brush head with other people.
- Do not place metal objects (coins, staples, etc.) on the charger.
- Persons with reduced mobility, reduced sensory perception, mental disadvantage or persons not familiar with the operation must use the device only under the supervision of responsible person, who is familiar with the instructions.
- Pay extra attention if the device is used near children.
- Do not use the device as a toy.
- **Before cleaning and after use of the device, switch the device off and disconnect it from the mains.**
- Do not let the power cord touch hot surfaces.
- Do not use the device in an environment where there are explosive gases and flammable substances (solvents, varnishes, adhesives etc.).
- Do not use coarse and chemically aggressive substances to clean the device.
- Do not use the device if it is not working properly, has been dropped, damaged or immersed in liquid. Have it checked and repaired by an authorized service centre.
- Do not use the device with a damaged cord or plug and have it repaired immediately by an authorized service centre.
- Do not wrap the power cord around the device body.
- Do not use the device outdoors.

- The device is only suitable for home use and is not intended for commercial use.
- **Do not immerse the power cord, plug or device in water or any other liquid.**
- Do not use any accessories other than those recommended by the manufacturer.
- Do not repair the device yourself. Contact an authorized service centre.
- This device can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the device in a safe way and understand the hazards involved. Cleaning and maintenance performed by the user must not be carried out by children unless they are 8 years old and under supervision. Children under the age of 8 must be kept away from the device and its power supply. Do not let children play with the device.

Damages from failure to follow the manufacturer's instructions are not covered by the warranty.

CAUTION:


This device is waterproof. It is suitable for cleaning under running water. Do not use or clean the device while it is connected to the power cord.

PRODUCT DESCRIPTION

- 1 Replaceable brush head
- 2 On/Off/Mode change button
- 3 Handle
- 4 Charging connector
- 5 Everyday Clean mode indicator
- 6 Soft Clean mode indicator
- 7 Massage Clean mode indicator

Package contains

- 1x Toothbrush
- 4x Heads (2 with small head, 2 with large head)
- 1x USB (DC) Charging Cable

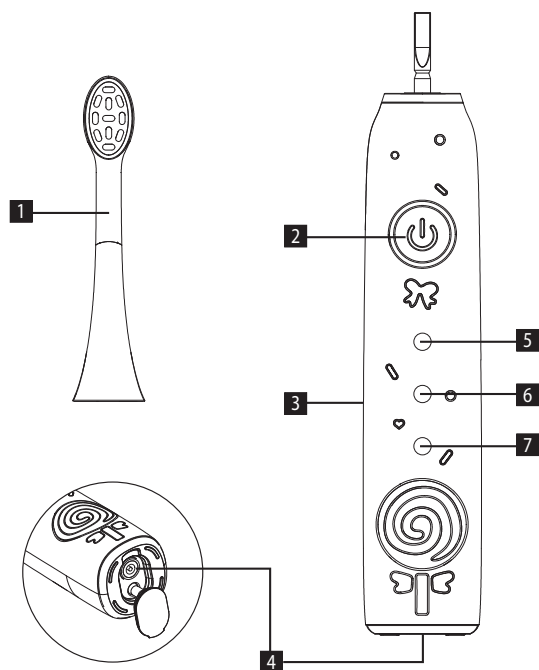OPERATING INSTRUCTIONS

Charging

- Fully charge the battery before using the device for the first time.
- If the device has not been used for more than 6 months, charge it for more than 8 hours.
- When using for the first time or when the battery has not been charged for a long time, the indicator may not light up for the first 5 minutes.
- When charging, the cleaning mode indicators will light up sequentially from the bottom to the top.
- When fully charged, the charging will turn off by itself.
- Charge only when the battery is low.
- When fully charged, the indicator light is white.
- To make better use of the battery, we recommend discharging the battery once a year.

Smart timer


- The recommended brushing time is at least 2 minutes. The device will automatically turn off after two minutes.
- Every 30 seconds, the device will pause to remind you to move to the next part of your mouth.

Cleaning modes


Clean - Standard choice for cleaning, suitable for everyday cleaning. Removes 100% of plaque.

Soft - Gentle cleaning of teeth and gums


Massage - Cleaning and massaging gums

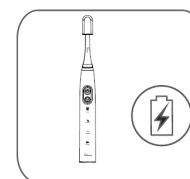
Fully charge the battery before using the device for the first time - the indicator lights up white.

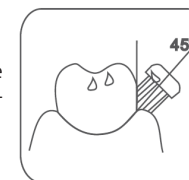
To clean your gums effectively, perform the cleaning at an angle of 45°.

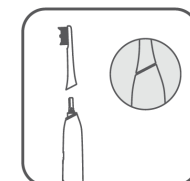
Place the brush head on the brush body.

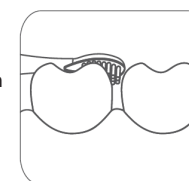Apply the brush lightly to avoid damaging the gums.

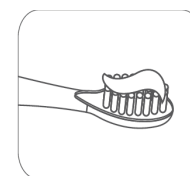Rinse the brush head, apply the paste, put the brush head in your mouth and then turn on the brush.

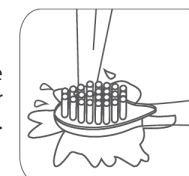
Rinse the brush head with water and dry.

TROUBLESHOOTING

Problem	Possible cause	Remedy
Weak vibrations	Too much pressure on teeth.	Reduce pressure on teeth.
	Improperly fitted brush head.	Fit the brush head again.
Loose brush head or abnormal sound	The brush head is damaged.	Replace the brush head.
Shorter operating time with a charged battery	Charging error.	Charge the toothbrush correctly.
	Empty battery.	Charge the brush to full charge.
The toothbrush does not work	Power has run out.	Charge the battery.

CLEANING AND MAINTENANCE

Brush head - After use, rinse with water and dry.
 Brush body - Rinse with water and dry. Do not use water warmer than 40 °C or any cleaning agents.
 Replacing the brush head - we recommend replacing the brush head every 3 months.
 Store the device in a dry and dust-free place.
 Do not expose to high temperatures.
Warning: The device contains a rechargeable lithium battery. Recycle accordingly.

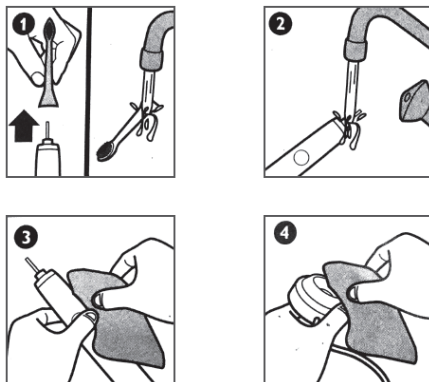SERVICE

Extensive maintenance or repair, which requires access to internal parts of the product, has to be performed by a qualified specialist or the service centre.

ENVIRONMENTAL CONCERNS

- Packaging materials and old appliances should be recycled.
- Packaging materials may be disposed of as sorted waste.
- Dispose the plastic bags made of polyethylene (PE) of as sorted waste.

Recycling of the appliance at the end of its lifespan:

This appliance is labelled in accordance with European Directive 2012/19/EU concerning used electrical and electronic appliances (Waste Electrical and Electronic Equipment - WEEE). The symbol on the product or on its packaging indicates that this product may not be treated as household waste. Instead, it should be taken to the appropriate collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. The disposal of such material has to be done according to the recycling regulations. For more detailed information about recycling of this appliance see the authorized local office or service for household waste disposal or the shop, where you purchased the appliance.

The product meets all the necessary requirements of the EU directives applicable to the specific product.

Changes in text, design, and technical specifications are subject to change without prior notice and we reserve the right to change them.

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