Nanoleaf LED Light Therapy Face Mask

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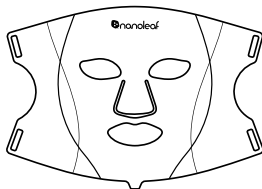
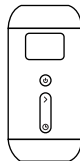
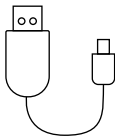
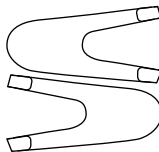
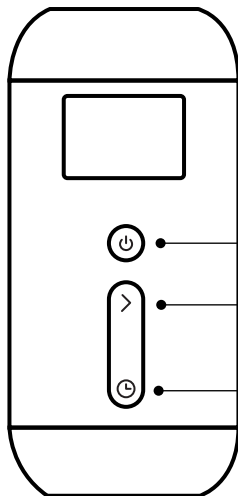
Power Button

Hold for 2s to Turn On / Off

By default, your device is set to turn off after 10 minutes. Press the power button once more when on to pause the timer.

Mode Selection


Press once to switch between Red, Green, Blue, Purple, Cyan, Yellow, & White light therapy.

Each light mode is combined with Near Infrared Light (NIR) to enhance the treatment effects. The screen will display the corresponding light mode. Press twice quickly to cycle between 5 different levels of intensity. The screen will display the corresponding light intensity.

Timer

Press once to add 1 minute to your timer. Long press to add more minutes quickly up to a max of 30 minutes. Your device will automatically turn off when your timer is up.

Contrôles

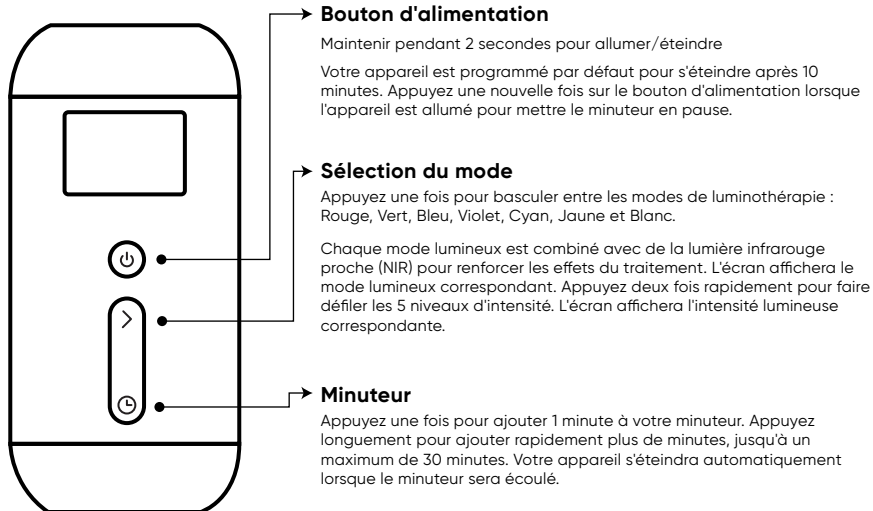Light Therapy Modes

Each light mode is combined with Near Infrared Light (NIR) that has 850nm penetration to enhance the treatment effects.

● Red Light (640nm + 850nm)

Mode R

Increases cell activity and promotes collagen formation to improve skin elasticity, and reduce the appearance of wrinkles.

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● Green Light (525nm + 850nm)

Mode G

Enhances the stimulation of cell energy to help even out skin tone, reduce redness and hyperpigmentation.

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● Blue Light (460nm + 850nm)

Mode B

Antibacterial and anti-inflammatory to help reduce acne breakouts and repair the skin.

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● Cyan Light (460nm + 525nm + 850nm)

Mode C

Soothes and calms stressed out skin. Promotes healing and treats inflammation.

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Light Therapy Modes

Each light mode is combined with Near Infrared Light (NIR) that has 850nm penetration to enhance the treatment effects.

Yellow Light (640nm + 525nm + 850nm)

Mode Y

Promotes tissue healing to help treat rosacea, reduce redness and improve overall skin tone.

White Light (640nm + 525nm + 460nm + 850nm)

Mode L

Penetrates the skin deeper to accelerate skin metabolism and enhance nutrient absorption.

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Purple Light (640nm + 460nm + 850nm)

Mode P

Multi-frequency treatment that targets acne and acne scars.

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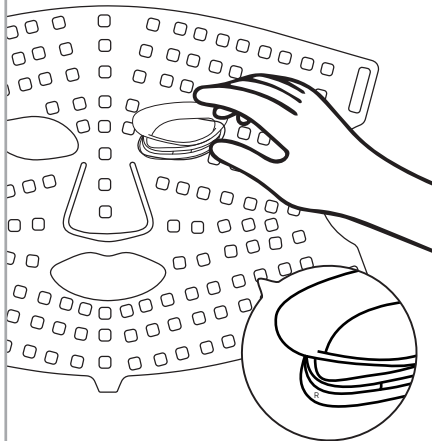
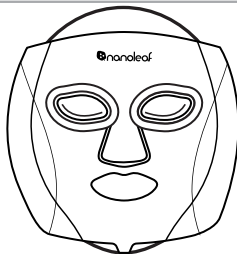
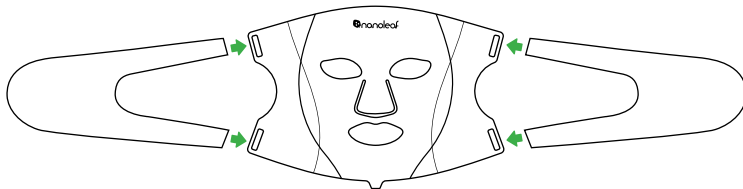
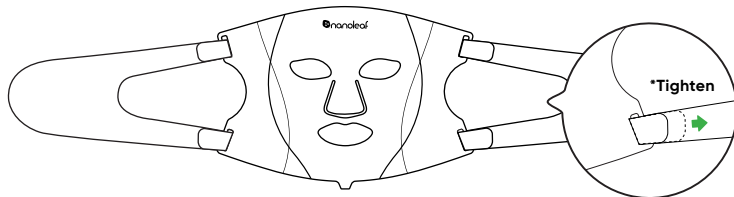
Things You Need to Know

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

WARNING

This is not a medical device and cannot be used for professional treatment or as a substitute for professional medical advice. The actual effect of use varies from person to person.

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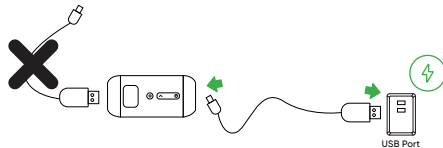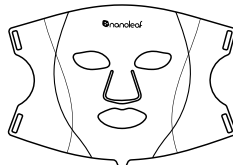Do not use the device while it is charging.

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If the device becomes excessively hot, emits smoke, or produces an unusual odor during charging, please immediately disconnect the charger and contact customer service.

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
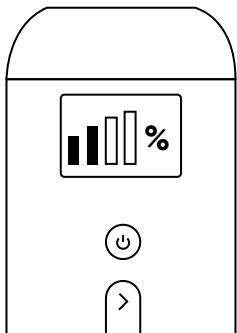
Portable

This product is rechargeable and portable and does not need to be connected to a wall power source while in use.

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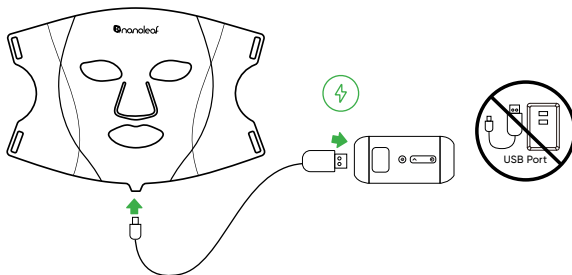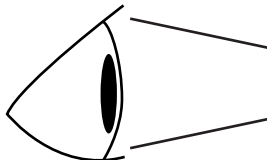
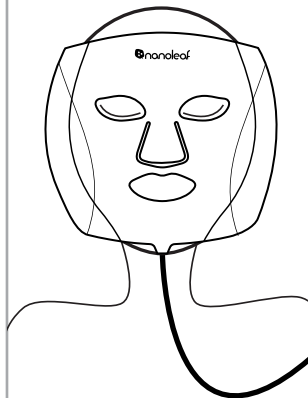
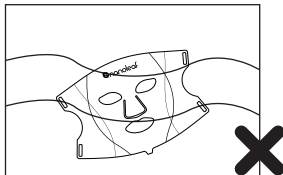
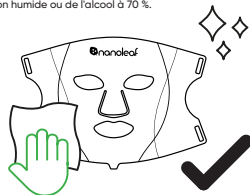
EN WARNINGS This product contains an infrared spectrum. Do not use this device for more than 30 minutes at a time. Pregnant women, children, individuals with epilepsy, thyroid patients or those on thyroid medication should not use this product. Do not use with photoallergies, or while taking drugs that can cause photoallergies or photosensitivity. Avoid wounds, and areas with skin trauma. Individuals that utilize special medical implants, heart rate monitors or have special medical conditions should consult with a doctor before use. Please note that improper use of LED face masks can potentially lead to skin darkening or other issues. Individuals that utilize special medical implants, heart rate monitors or have special medical conditions, facial wounds, or photosensitivity should consult with a doctor before use. Do not dispose of battery into fire or dismantle it under any condition. This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

WARNINGS: SAFE USE OF BATTERIES 1. Do not attempt to disassemble or modify the battery. Doing so can lead to leakage, fire, or explosion, and will void the warranty. Proper disposal and recycling methods should be followed for used batteries. Always refer to the manufacturer's guidelines for battery care and maintenance. 2. Do not leave your device's battery unused for extended periods. For optimal battery health, charge the battery at least once every three months. 3. To maintain battery life, ensure the device is charged within a reasonable time after use. An excessive discharge time can adversely affect battery performance. 4. Avoid long-term storage of the device with a depleted battery. Charge the battery to at least 50% before storing for long periods to prevent permanent damage.

CE INFORMATION NanoGrid Limited (Nanoleaf®) hereby declares that this device is in compliance with the essential requirements and other relevant provisions of Directive 2014/53/EU. The declaration of conformity is available at <https://nanoleaf.me/eudocs>. **UKCA DECLARATION** NanoGrid Limited (Nanoleaf®) hereby declares that this device is in compliance with the product requirements and other relevant provisions of the Radio Equipment Regulations 2017. The declaration of conformity is available at <https://nanoleaf.me/en-GB/products/declarations/>. **CLEANING AND MAINTENANCE** Do not clean the device with liquid or sprays when it is connected to power. Clean with a soft dry cloth, and do not use harsh cleaning agents.

WARRANTY INFORMATION This product has a one year warranty from the date of purchase. Damages caused by time and use is not within the scope of the service guarantee. The warranty scope of this product does not include failures caused by any factors other than the defects of the product itself.

Disposal of your old product: Your product is designed and manufactured with high-quality materials and components, which can be recycled and reused. When this crossed-out wheeled bin symbol is attached to a product, it means that the product is covered by European Directive 2012/19/EU. Please find out about your local collection system for electrical and electronic products. Please act according to your local rules and do not dispose of your old products with your normal household waste. The correct disposal of your old product will help to prevent potential negative consequences for the environment and human health.

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Nanoleaf Europe SAS