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EN USER MANUAL

MODEL: BL-101

Thank you for purchasing Multi-Therapy Device BL-101. Established in 1939 **TECH - MED**[®] TECHNIKA MEDYOZNA has always been with you for over 75 years – we use our experience to continually improve our products, which incorporate all the latest technologies in diagnostics and are constantly subjected to quality control. Your choice is the best possible one. We wish you good health. We are always at your disposal.

SAFETY INSTRUCTIONS

To ensure proper operation of this device, please read this instruction manual carefully before using the device and keep it for further reference.

WARNING

The following persons should not use this device: pregnant women or while nursing; children; persons who have a pacemaker or other electronic medical equipment; persons with allergic reactions; persons with skin diseases (skin blemishes, burns, etc.); persons who epilation or plastic surgery was performed on the areas of use: on the eyelids and throat; while wearing makeup.

- This device has different heads suitable for use on the face (forehead, cheeks and chin) and scalp. The device should not be used on the eyelids.
- If you are sensitive to light, prone to rashes or allergic reactions, do not use the light therapy mode.
- Do not use the light therapy mode if the treated area has infections, egzema, burns, folliculitis, open lacerated wounds, herpes, wounds or lesions, and hematomas.
- People with diabetes, lupus erythematosus, porphyria or congestive heart disease, bleeding disorders, history of immunosuppressive diseases (including HIV infection and AIDS) should not use this device.
- Do not use this device while undergoing radiation therapy or chemotherapy.
- Do not use this device if you are taking any of the following medications:
 - If your skin is currently being treated or has recently been treated with alpha hydroxy acids (AHAs), alpha hydroxy acids (AHAs), beta hydroxy acids (BHAs), topical isotretinoin and azelaic acid.
 - If you have taken any form of Accutane A or Roaccutane A in the past six months. Such treatment makes the skin susceptible to injury, trauma and irritation.

- If you are taking photosensitizers or any medications, read the patient information before using the device. Do not use the device if the information leaflet mentions that the drug may cause photoallergic and phototoxic reactions, or avoid exposure to sunlight while taking the drug.
- If you are taking anticoagulants, including high doses of aspirin, use the device one week after stopping treatment.
- Do not repair, disassemble (except for replacement heads) or modify the device, as improper handling may damage it. Do not use the device for purposes other than those indicated in this manual.
- Do not look directly at the light when using the light therapy head, as this may cause discomfort. Use safety glasses to protect your eyes.

DEVICE DESCRIPTION

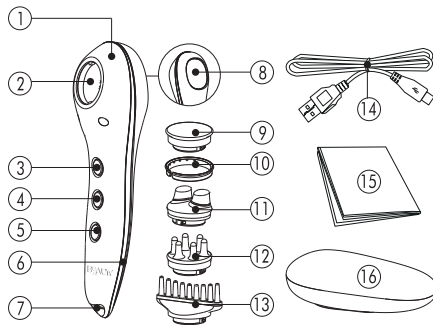
Figure 1

- | | |
|---|------------------------------|
| 1. Body of the device | 9. cosmetic head ION+/- |
| 2. connection point between the head and the device | 10. cotton ring |
| 3. power button | 11. microcurrent facial head |
| 4 ION+ /- MC (microcurrent) button | 12. scalp massage head |
| 5. LED button | 13. microcurrent scalp head |
| 6. conductive rod | 14. USB-C cable |
| 7. charging port | 15. cleaning cloth |
| 8. light therapy head | 16. storage box |

PREPARING THE DEVICE FOR USE

Charging the device

The device should be charged according to its charging system. Charging the device takes about 2.5 hours. When the battery is fully charged, the device can be used for about 1.5 hours.

The light indicator flashes red while charging and glows blue when the device's battery is charged.

If the battery charge is low, the microcurrent indicator will flash 3 times and the device will turn off automatically.

Figure 2

Removal/assembly of heads

Disassembly: grasp the body of the device with one hand, and with the other hand turn the head to the left and gently pull to remove it. (Figure 3)

Figure 3

Installation: grasp the body of the device with one hand, connect the head to the body of the device, and rotate clockwise until the position is locked (Figure 4)

Figure 4

BUTTONS

Power button:

You can turn on the device even if the head is not installed.

Press the power button, after a beep, the indicator will turn white, and the device will go into standby mode. If the device is not used, it will automatically turn off after 30 seconds.

Ion+/- / MC (microcurrent) button:

If the corresponding head is properly installed, press the ION +/- / MC button to activate the corresponding function. The device automatically recognizes the function of the installed head. If the accessory head is changed, the ION +/- / MC button will display the corresponding mode.

- After installing the ION+/- cosmetic head, press the ION +/- / MC button, the indicator will turn green after beeping, and the device will start the Ion+ (cleansing) mode. Press the ION +/- / MC button again, the indicator will turn yellow after the beep and the unit will enter Ion- (humidification) mode. Both modes change cyclically.
- After installing the scalp microcurrent head or facial microcurrent head, press the ION +/- / MC button. When the indicator beeps, it will turn green, and the device will start the low-intensity microcurrent mode. Please press the ION +/- / MC button again, after the sound is emitted, the indicator will turn yellow, and the device will switch to the high-intensity microcurrent mode.

LED Button:

- The LED light therapy head is built into the device. Press the LED button to activate the light therapy mode.
- Turn on the device and briefly press the LED button. When you hear a sound, the indicator will turn red and the red light mode will start. Press the LED button again, after hearing a sound, the indicator will turn yellow and the device will switch to flashing yellow light therapy mode.

TIME OF TREATMENTS

Face care:

ION+ face mode Cleaning 6 min	->	ION-face mode- 6 min	->	MC (microcurrent) face mode 6 min	->	Red light face 6 min	->	Yellow light face 6 min
--	----	--------------------------------	----	--	----	--------------------------------	----	-----------------------------------

Head care

ION+ head mode Cleaning 6 min	->	ION- Head Mode- 6 min	->	MC head mode (microcurrent) 6 min
--	----	---------------------------------	----	--

The device is designed for daily use.

USAGE

To ensure proper operation of the device, read the instructions carefully before use.

Daily facial care

Use the device on clean skin, making sure to remove makeup and dirt from the skin before use.

Do not look directly at the light when using the light therapy head, as this can cause discomfort to the eyes.

Ion + mode

Install the facial care head, place a cotton swab on the surface of the head and press the cotton ring down over the swab. Saturate the cotton swab with tonic (Figure 5)

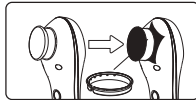
Figure 5

Press the power button and the ION+/-/MC button successively to start the purification mode. This mode has an audible vibration.

Grasp the device with both sides of your hand so that it touches the conductive bars. With the ION+/- head touching the skin, gently move it in a circular motion from the center to the sides of the face (Figure 6)

Figure 6

ION mode -

After using the ION + cleansing mode, please remove the cotton pad and wash your face with warm water. Switch to ION mode - moisturizing.

Apply lotion/serum to the face. Hold the device with both sides of your hand. Place the head against the face and gently move it in circular motions from the center to the sides of the face to facilitate the absorption of skin care products. (Figure 7)

Figure 7

MC (microcurrent) mode.

- After using ION- mode, install the microcurrent head to the face.
- Press the power button and the ion+/- / MC button successively. The device will start the low intensity microcurrent and sound vibration mode. Press the ION+/- / MC button again to start the high-intensity microcurrent mode.
- Move the head surface from the lower jaw to the cheekbone or from the corner of the mouth to the ear. Select the intensity of the microcurrent mode, after a short beep, move the accessory head surface as shown in Figure 8. Repeat the procedure 2 times.

Figure 8

Around the mouth: Select the intensity of the microcurrent, when the sound is heard, apply the surface of the head to the skin and move it as shown in Figure 9. Repeat the procedure 2 times.

Cheeks: To use the head on the cheeks, move the needle from the corner of the mouth to the ear. When you hear a sound, move the massage head slightly as shown in Figure 10. Repeat the procedure 2 times. (illustration 10)

Around the eyes: Massage according to the arrows in illustration 11 and repeat the procedure 2 times.

Figure 9

Figure 10

Figure 11

- Forehead: Close your eyes and slowly move the surface of the head along the forehead to the hairline, as shown in Figure 12. Repeat the procedure 2 times.

Figure 12

Red/yellow light therapy

- Before use, check if the skin is sensitive to light: Try the light therapy mode on the forearm for 3 minutes. If the area has become red after 3 minutes of use and the redness persists for more than 2 hours, it means the skin is sensitive to light. Stop using the light therapy mode.
- Do not look directly at the light therapy source to avoid discomfort to the eyes.
 - Press the power button and the LED button successively to turn on the red light mode. Press the LED button again to enter the yellow light therapy mode (intermittent mode).
 - Place the lamp head 2-3 cm away from the skin to evenly illuminate the skin surface. The light therapy head can also contact the skin when using this mode.

Use the device for 2 minutes on the left side of the face, 2 minutes on the right side of the face and 2 minutes on the forehead. Every 6 seconds, the device beeps to remind you to change the treatment area. (Figure 13)

Figure 13

Daily head care

ION + mode

- Mount the scalp massage head.
- Press the power button and the ION +/- MC button to start the ION + cleansing mode. This mode includes sound vibration.
- Comb through the scalp to cleanse the hair as recommended in Figure 14 in the order of 1-3.

Figure 14-1

Figure 14-2

Figure 14-3

ION mode -

- Press the ION +/- MC button again to switch to the ionic humidification mode. This mode includes sound vibration.
- Apply the preferred cosmetic to the scalp.
- Brush the scalp to increase the absorption of the applied product, apply as recommended in Figure 15 (order 1-3).

Figure 15-1

Figure 15-2

Figure 15-3

Microcurrent mode

- Place the microcurrent head on the scalp, press the power button and ION +/- MC button sequentially to start the microcurrent mode. The device will beep every 6 seconds to remind you to change the treatment area.
- Apply a scalp care cosmetic.

Brush the scalp and massage it in the order of 1-3 as recommended in Figure 16.

Figure 16-1

Figure 16-2

Figure 16-3

CLEANING AND MAINTENANCE

- After using the device, clean the heads and body of the device with a clean, slightly damp cloth. Do not clean directly under water.
- Do not use steel wool, abrasive cleaners or caustic liquids to clean the device.

STORAGE

This device is equipped with a lithium battery. Do not place or store the device in a hot or humid environment for a long time. When not using the device for a long period of time, store it after it is fully charged. Do not throw the device into a fire or dispose of it in general household waste. Dispose of the device at an official recycling center.

DEVICE SPECIFICATIONS

Name:	Multi-Therapy Device
Model:	BL-101
Input:	5 V $\equiv$ 0,5 A
Battery type:	Li-ion 3,7 V 500 mAh
Executive Standard no:	Q/XM1002-2020 Q/XM1021-2020

TROUBLESHOOTING

Problem	Possible cause	Solution
The device is not working properly.	The head is not properly installed.	Reinstall the head in accordance with According to the instructions.
	Battery exhausted	Charge the device according to the instructions.
The microcurrent is weak.	Skin too dry.	Use a toner or serum.
Other problems.		Please consult the manufacturer.

Note: If the device still does not work properly after applying the above solutions, stop using the device and contact the manufacturer.

DISPOSAL

The built-in battery contains substances that can pollute the environment. When disposing of this device, do not dispose of it in general household waste. Dispose of this device at an official recycling point to contribute to environmental protection.

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