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EN USER MANUAL

MODEL: BL-104

Thank you for purchasing the BL-104 Sonic Facial Cleanser.

has always been with you for more than 75 years - we use our experience to continuously improve our products, which incorporate all the latest technologies in diagnostics and are constantly subjected to quality control. Your choice is the best possible one.

We wish you much good health. We are always at your service.

SAFETY INFORMATION

To ensure the proper operation of this device, please read this manual carefully before using the device and keep it for future reference.

WARNING

- The device should not be used by: pregnant or breastfeeding women; people who are allergic to silicone; people with facial nerve disorders; people who have metal, silicone or plastic implants on the face due to fracture or plastic surgery; on areas with open wounds, chapped or swollen skin; if there are serious skin diseases or any other facial skin diseases.
- Check the device before use. Do not use if the device is damaged.

- If you experience discomfort while using the device, stop using it immediately and consult a doctor.
- Do not use the device while it is charging.
- Do not attempt to modify, disassemble or repair the device yourself.

NOTE

- Do not share the same device with other people.
- Keep the device out of the reach of children.
- Clean the device regularly to prolong its life.

This unit is equipped with a rechargeable lithium battery. Do not place or store the device in a hot or humid environment for a long time.

DEVICE DESCRIPTION

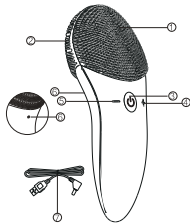
Figure 1

01 Cleansing silicone brush surface
02 Exfoliating silicone brush surface
03 Power button
04 Button +

05 Button -
06 Charging port
07 Charging cable

PREPARING THE DEVICE FOR OPERATION

Charging the device takes about 2.5 hours. When the battery is fully charged, the device can be used for about 3 hours. When the battery charge is low, the light indicator glows red. During charging, the light indicator flashes red. When the battery is fully charged, the light indicator will light up blue.

- **Turn on:** Briefly press the power button to turn on the device. Depending on your skin type and preference, briefly press the - button or the + button to adjust the intensity of cleansing.
- **Turning off:** The device is equipped with an automatic timer that automatically turns off after 90 seconds. To turn off the device before 90 seconds, briefly press the power button.

METHOD OF USE

- **Pore cleansing**
 - Moisten the face with water, apply the gel/cream cleanser to the palm of one hand, hold the device with the other hand and rub the cleansing surface on the gel/cream to create foam.
 - Apply the foam to the face, gently moving the cleansing surface in a circular motion from the inside to the outside of the face. Recommended cleansing area, order and time of application: T-zone (forehead, nose)

20 seconds, left and right cheeks 15 seconds each (30 seconds total), chin 10 seconds, neck 30 seconds (Figure 2).

Figure 2

- **Exfoliation of the epidermis**

- Moisten the face with water, apply the gel/cream cleanser to the palm of one hand, hold the device with the other hand and rub the exfoliating surface on the gel/cream to create foam.
- Apply the foam to the face and use the exfoliating surface to exfoliate the skin. Gently move the exfoliating surface in a circular motion from the inside to the outside of the face. Recommended cleansing area, order and time of application: T-zone (forehead, nose) 20 seconds, left and right cheeks 15 seconds each (30 seconds total), chin 10 seconds, neck 30 seconds (Figure 3).

For sensitive skin, it is recommended to use a silicone exfoliating brush once a week. For normal, combination or oily skin, use the silicone exfoliating brush according to the skin's needs.

- After use, turn off the device, rinse your face with water and apply moisturizer.

Figure 3

CLEANING AND MAINTENANCE

- Make sure the device is turned off before cleaning.
- Do not immerse the device in water. Wash the surface of the device under running water and dry with a dry cloth.
- Do not use lubricants, detergents, polishes or air fresheners.
- Do not use hot water above 40°C or caustic liquids to clean silicone surfaces.

- Do not use clay-based, silicone-based or gritty cleaners and exfoliators, as they may damage the surface of the brush head.

STORAGE

This device is equipped with a rechargeable lithium battery. Do not place or store the device in a hot or humid environment for an extended period of time. If you do not intend to use the device for a long period of time, please store it fully charged.

TECHNICAL SPECIFICATION

Name: Sonic facial cleanser

Model: BL-104

Input: 5.0 V $\overline{\text{=}}$ 0.5 A

Battery type: Li-ion 3.7 V 500 mAh

Performance standard: Q/XM 1001-2020

TROUBLESHOOTING

Problem	Possible cause	Solution
When the power button is pressed, the device does not work properly	Battery (almost) depleted	Charge the device.
The charging indicator does not light up.	Charging cable incorrectly inserted.	Insert the charging cable properly into the charging port.
Other problems.		The manufacturer should be consulted.
Note: If the device still does not work properly after applying the above solutions, stop using the device and contact the manufacturer.		

DISPOSAL

The built-in battery contains substances that can pollute the environment. When disposing of this device, do not dispose of it in general household waste. Dispose of this device at an official recycling center to contribute to environmental protection.

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