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EN USER MANUAL

MODEL: BL-107

Thank you for purchasing the BL-107 Thermolifting for the eye area.

has always been with you for more than 75 years - we use our experience to continuously improve our products, which incorporate all the latest technologies in diagnostics and are constantly subjected to quality control. Your choice is the best possible one. We wish you much good health. We are always at your service.

SAFETY INFORMATION

To ensure proper operation, please read this manual carefully before using the device and keep it for future reference.

WARNING

- The device is not a toy, keep it away from children.
- Check the device before using it. Do not use the device if it is damaged.
- Keep the USB cable and charging port dry.


- Do not use the device while it is charging.

NOTE

- The device should not be used by: pregnant, breastfeeding or menstruating women; allergic people; people with facial nerve disorders; people with serious facial skin diseases; people with eyelash implants.
- Do not use the device while wearing contact lenses.
- The device is not intended for use by persons (including children) with limited physical, sensory or mental abilities, or persons without experience and knowledge, unless they are supervised or instructed in the use of the device by a person responsible for their safety.
- Stop using the device immediately and consult a doctor if any discomfort is felt.
- When using red light therapy, close your eyes and do not look directly at the light source.
- Keep the device out of the reach of children.

- Do not share the same device with other people.
- Clean the device and accessories regularly to prolong their life.
- Do not attempt to disassemble, repair or modify the device yourself to avoid accidents.

DEVICE DESCRIPTION

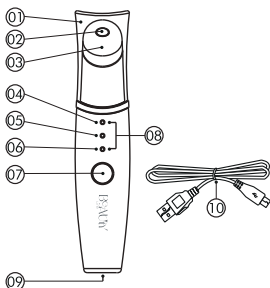- 1 Protective cap
- 2 Light source
- 3 Eye massage head
- 4 High temperature indicator
- 5 Sound indicator

- 6 Red light indicator
- 7 Power/mode button
- 8 Charging indicator
- 9 USB type C charging port
- 10 USB type C cable

BEFORE USE


The battery of the device is discharged before use. Charge the device with the USB cable before using it. The device has a low battery alarm system. When the battery level is low, the three indicators will flash 3 times and the device will turn off automatically.

Charging the device:

- Make sure the device is turned off.
- Connect the USB cable to the USB charging port on the bottom of the device, and then plug the power adapter into a wall outlet (power adapter not included). (Figure 2)
- During charging, the three indicators light up periodically. When the battery is fully charged, the three indicators light up solid. Charging the device takes about 2 hours.

- When the battery is fully charged, the device can be used for about 30 minutes. Use 10 times, 3 minutes each time.

Note: Do not attempt to replace the device's battery.

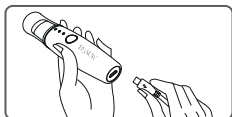USER MANUAL

To ensure proper operation of this device, please read the instructions carefully before use.

1. Operation of the device

Press and hold the power/mode button to turn on the device and hot mode. Press the power/mode button again to switch between modes. The device has 7 modes: hot mode - sound mode - red light mode - hot mode + sound mode - hot mode + red light mode - sound mode + red light mode - hot mode + sound mode + red light mode (Figure 3)

In any mode, press and hold the power/mode button to turn off the unit. The smart timer automatically turns off the device after 5 minutes.

2. Usage

It is not recommended to use the device on the eyelids to avoid injury.

- Apply a layer of eye cream/essence to the skin around the eyes (if using the mode with light therapy, do not apply any substances to the skin).
- Press and hold the power/mode button to turn on the device. Select the mode to be used according to your skin's needs (press the power/mode button briefly to switch between modes).
- Gently press the massage head against the skin, slowly move it and apply slight pressure from the inner to the outer corner of the eye. (Figure 4)

- Recommended time and area of use: use the device on the left eye for 1.5 minutes and on the right eye for 1.5 minutes. For better results, use daily.
- After using the device, turn it off and follow your normal skin care regimen.
- After each use, clean the body of the device and the massaging head with a clean, slightly damp cloth. Store the device in a dry, well-ventilated place to prolong its life.

CLEANING AND MAINTENANCE

- Before and after every single use, clean the housing of the device to ensure its best performance.

- Make sure the device is turned off before cleaning.
- Clean the device with a clean, slightly damp cloth. Do not immerse the device in water.
- Do not use alcohol, steel wool, abrasive cleaners, chemicals or caustic liquids to clean the device.
- Do not use lubricants, detergents, polishes or air fresheners to clean the device.

STORAGE

The device is equipped with a lithium battery. Do not place or store the device in a hot or humid environment for a long time. If the device will not be used for an extended period of time, store it after it is fully charged. Do not throw the device into a fire or dispose of it in general household waste. Take the device to an official recycling center for disposal.

PRODUCT SPECIFICATIONS

- Name:
Thermolifting for the eye area
- Model: BL-107
- Input: 5V $\equiv$ 0.5A
- Battery type:
Li-ion 3.7V 300 mAh
- Performance standard:
Q/XM 1015-2020
Q/XM 1021-2020

TROUBLESHOOTING

Problem	Possible cause	Solution
The device is not working.	Battery (almost) depleted.	Charge the device according to the instructions.

The massage head heats up slowly.	Room temperature below 15°C.	To make the massage head warm up faster, use the device in a room with a temperature above 15°C.
Other problems.		Please consult the manufacturer.
Note: If the device still does not work properly after applying the above solutions, stop using the device and contact the manufacturer.		

DISPOSAL

The built-in battery contains substances that can pollute the environment. When disposing of this device, do not dispose of it in general household waste. Take this device to an official recycling center to contribute to environmental protection.

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