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beautyline.com.pl

BEAUTY
Line

MULTIFUNCTIONAL MASSAGER
FOR THE EYE AREA

EN USER MANUAL

MODEL: BL-105

Thank you for purchasing the BL-105 Multifunctional massager for the eye area.

has always been with you for more than 75 years - we use our experience to continuously improve our products, which incorporate all the latest technologies in diagnostics and are constantly subjected to quality control. Your choice is the best possible one.

We wish you much good health. We are always at your service.

SAFETY INFORMATION

To ensure the proper operation of this device, please read this manual carefully before using the device and keep it for future reference.

WARNING

- This device is not a toy. Please keep it out of the reach of children.
- Check the device before using it. Do not use the device if it is damaged.
- Keep the charging cable and charging port dry.
- Do not use the device while charging.

NOTE

- The device should not be used by: allergic people, people with facial nerve disorders, people with serious facial skin diseases, etc.
- Do not use the device while wearing contact lenses.
- This device is not intended for use by persons (including children) with limited physical, sensory or mental abilities or lack of experience and knowledge, unless they are supervised or instructed in the use of the device by a person responsible for their safety.
- Immediately discontinue use and consult a physician if any discomfort is felt.
- Keep the device out of the reach of children.
- Do not share the same device with others.
- Clean the device and accessories regularly to prolong its life.
- Do not attempt to disassemble, repair or modify the device yourself to avoid accidents.
- Do not throw the device into fire or dispose of it in general household waste. Dispose of this device at an official recycling center.

DEVICE DESCRIPTION

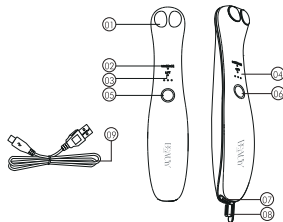
Figure 1

- 01 Eye massage head
- 02 Massage mode indicator
- 03 LIFT mode indicator
- 04 Intensity level indicator light
- 05 Charging indicator

- 06 Power/mode button
- 07 Type C charging port
- 08 Rubber plug
- 09 USB type C charging cable

PREPARATION INSTRUCTIONS

- Charge the device's battery before use. Charge the device with a type C charging cable before first use.
- The device has a low battery reminder function. When the battery level is low, the charging indicator flashes 10 times and 5 low battery beeps are emitted. Please charge the device before the battery runs out.

Charging the device:

- Make sure the device is turned off.
- Open the rubber stopper to access the charging port. Connect the charging cable to the device's charging port. Then plug the adapter into a wall outlet (adapter not included). (Figure 2)
- During charging, the charging indicator flashes. When fully charged, the charging indicator lights up with a steady light. Charging time is approximately 3 hours.
- A full battery allows about 60 minutes of continuous use. It can be used 7-8 times, 8 minutes each time.

Caution: Do not attempt to replace the device's battery.

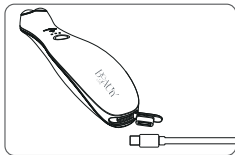
Figure 2

USAGE

To ensure proper operation of this device, please read the instructions carefully before use.

- Operation of the device

- 1) Mode instructions

- MASSAGE mode combines a 42°C hot compress with constant temperature + sound vibration to help soothe the skin.

This mode can be used with or without a skin care product. When using this mode with a skin care cosmetic, choose one that suits your skin's needs.

- LIFT mode combines EMS microcurrent + hot compression at a constant temperature of 42°C to help lift and tighten the skin and smooth fine lines around the eyes.

Use this mode with an eye essence, eye cream or gel to enhance performance and conductivity. It is not recommended to use this mode on dry skin.

2) Button instructions

- 1) On/Off Switch: Press and hold the power/mode button to turn the unit on/off.
- 2) Mode switch: Press the power/mode button briefly to switch between modes and intensity levels. The intensity level indicator on the device is "MASSAGE level 1" - "MASSAGE level 2" - MASSAGE level 3".

"LIFT level 1" - "LIFT level 2" - "LIFT level 3". Switch between 6 modes in a cycle. The MASSAGE mode has three intensity levels that control the intensity of the vibration, while the LIFT mode has three levels that control the intensity of the microcurrent. You can adjust the intensity level from low to high. (Illustration 3)

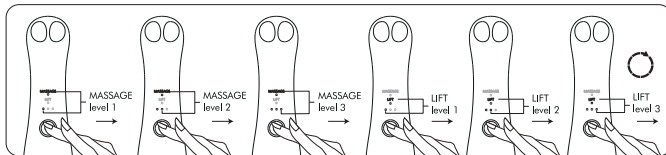
Figure 3

- 3) Quick mode switch: When the device is turned on, the MASSAGE mode automatically starts by default. Press the power/mode button twice in short intervals to switch modes directly.

Time reminder

If no button is pressed while the power is on, the device will automatically turn off after 6 minutes. Do not use the device if there is inflammation of the eyes or if they are swollen.

Eye massage

- Press and hold the power/mode button to turn on the device and the MASSAGE mode. Select the intensity level according to your skin needs and preferences. Gently press the massage head against the skin, slowly move or apply light pressure to the eye area. We recommend using the MASSAGE mode for 1-2 minutes per eye. (Figure 4).

Figure 4

Wrinkles around the eyes

- Before using the LIFT mode, use the MASSAGE mode to stimulate the eye area and increase the absorption of skin care products. Apply a layer of eye cream/essence to the eye area.

Switch to LIFT mode and select the level of microcurrent intensity according to your preference and skin needs. Gently press the massage head against the skin and slowly move it from the inner corner of the eye toward the temple. Then move the massage head along the outer corner of the eye, repeating the lifting movements. We recommend using LIFT mode for 1-2 minutes on each eye. (Figure 5)

Figure 5

* After using the device, turn it OFF and follow your skin care routine.

CLEANING AND MAINTENANCE

Clean the device housing before and after each use to ensure the best performance of the device.

- Before cleaning, make sure the device is turned off and the rubber plug of the charging port is tightly closed.
- Clean the device with a clean, slightly damp cloth. Do not immerse the device in water.
- Do not use alcohol, steel wool, abrasive cleaners, chemicals or caustic liquids to clean the device.
- Do not use lubricants, detergents, polishes or air fresheners to clean the device.

STORAGE

This device is equipped with a rechargeable lithium battery. Do not place or store the device in humid or hot places for a long time. If the device will not be used for a long time, charge it before storing.

DEVICE SPECIFICATIONS

Name: Multifunctional massager for the eye area

Model: BL-105

Input: 5.0 V $\equiv$ 0.5 A

Battery type: Lithium-polymer 3.7 V 800 mAh

Performance standard: Q/XM 1002-2020; Q/XM 1015-2020

TROUBLESHOOTING

Problem	Possible cause	Solution
The device is not working	Battery (almost) depleted	Charge the device according to the instructions.
The current is felt to be insufficiently strong.	Low conductivity.	Apply conductive essence/gel.
Other problems.	It requires a specific analysis.	The manufacturer should be consulted.
Note: If the device still does not work properly after applying the above solutions, stop using the device and contact the manufacturer.		

DISPOSAL

The built-in battery contains substances that can pollute the environment. When disposing of this device, do not dispose of it in general household waste. Dispose of this device at an official recycling center to contribute to environmental protection.

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