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EN USER MANUAL

MODEL: BL-108

Thank you for purchasing BL-108 Hot & Cool BL-108 Facial Massager

has always been with you for over 75 years – we use our experience to continually improve our products, which incorporate all the latest technologies in diagnostics and are constantly subjected to quality control. Your choice is the best possible one. We wish you good health. We are always at your disposal.

SAFETY INSTRUCTIONS

To ensure proper operation, please read this manual carefully before using the device and keep it for future reference.

WARNINGS

- The device is not a toy and must not be used by children.
- Always check the device before using it. Do not use it if the device is damaged or broken.
- If the device is used in a dry environment, its handling may cause static electricity. This is a normal phenomenon and does not affect the use of the device.

- Keep the USB cable and charging connector of the device dry.
- The device is not waterproof, do not wash it with water.
- Do not use the device while charging.

PRECAUTIONS

- Do not use the device in the following cases: being pregnant, during a menstrual cycle, or after childbirth; having a pacemaker or other electronic medical equipment; suffering from allergies, skin diseases or inflammations, etc.
- Do not use this device in the following areas: throat, areas after plastic surgery or depilation.
- Do not use the device with makeup on.
- The device is not intended for use by persons (including children) with limited physical, sensory or physical, sensory or mental abilities, or lack of experience and knowledge, unless they are under supervision or have been instructed to use the device by a person responsible for their safety.
- Immediately discontinue use of the device and consult a physician if any discomfort is felt while using the device.
- Keep the device out of the reach of children.
- Clean the surface of the device and accessories regularly to prolong their life.

- Do not share the device with others.
- Do not repair, disassemble or modify the device yourself.
- Remove jewelry/accessories (such as necklaces and earrings) and tie up hair to avoid entanglement with the device.
- If the device is not used frequently, charge it before storing.

DEVICE DESCRIPTION


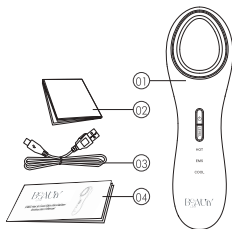
Figure 1

Set composition

01 Main unit
02 Cleaning wipe

03 USB type C charging cable
04 User Manual

Device components

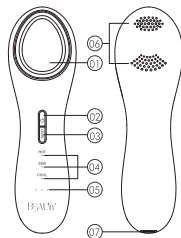- 01 Facial massage head
- 02 Power/mode switch button
- 03 Intensity level button

- 04 Mode indicator
- 05 Intensity level indicator
- 06 Ventilation holes
- 07 Type-C charging port

Figure 2

PREPARING THE DEVICE FOR OPERATION

Before using the device for the first time, the battery is discharged. Charge the device with the USB cable. Before using the device, remove the blue protective film from the surface of the massage head with a clean, slightly dampened cloth.

- Charging the device
 - Make sure the device is turned off.
 - Insert one end of the Type-C charging cable into the charging socket on the bottom of the device, and the other end into the adapter (adapter not included). Plug the adapter into a wall outlet to start charging.
 - During charging, the intensity level indicator lights up periodically. When fully charged, the indicator lights up with a steady light.
 - Charging the device takes about 3 hours.
 - A full battery allows about 40 minutes of continuous use.
 - The device cannot be used while charging.

- Turning on, off and switching modes

- Press and hold the power button for 1 second to turn on the device. When the device is on, press and hold the power button for 1 second to turn it off. (Figure 3)

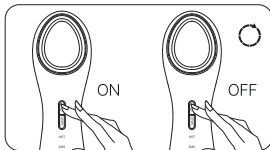
Figure 3

- The device has three modes of operation: HOT, EMS and COOL. The device activates HOT mode by default when you turn it on for the first time or when the battery runs out. Press the power/mode button to switch to EMS mode, then press it again to switch to COOL (cold massage) mode. (Figure 4)
- The device has an intensity level memory function. It automatically activates the mode used before the device was turned off. However, if the battery runs out, the device activates the default mode setting.

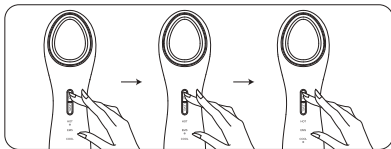
Figure 4

- Adjusting the intensity level
 - HOT and EMS modes have three intensity levels. To change the intensity level, press the intensity level button. You can adjust the vibration intensity in HOT mode and the current intensity in EMS mode. (Figure 5)
 - The temperature setting in COOL mode cannot be adjusted.

CAUTION: Do not attempt to replace the device's battery.

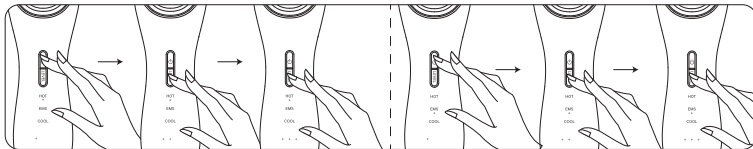
Figure 5

USE OF THE DEVICE

To ensure proper operation of the device, read the instructions carefully before use.

- Instructions for use
 - Clean the skin and remove makeup before using the device.
 - The device automatically shuts off after 6 minutes of use.

- When the device is turned off, the device's fan will only turn off after 1 minute. This is normal and does not affect the use of the device.

- **HOT mode**

Turn on the device and select HOT mode to start heat and vibration therapy. Adjust the intensity level to your skin's needs and preferences. This mode helps open pores, promotes absorption of skin care products and increases hydration levels. You can use this mode with a toner/essence or your favorite facial mask to increase nutrient absorption. It is recommended to use the device in this mode for 2 minutes on the left cheek, 2 minutes on the right cheek and 2 minutes on the forehead.

- **EMS mode**

Switch to EMS mode and adjust the current intensity level to suit your needs and skin tolerance. When using the device for the first time, start at intensity level 1. When using the EMS mode, apply a water-soluble gel to help conduct the EMS current. It is recommended to use the device in this mode for 2 minutes on the left cheek, 2 minutes on the right cheek and 2 minutes on the forehead.

TIPS: EMS massage helps define facial contours and smooth fine lines on the chin, cheeks and forehead. For better results, the device should be used daily or as needed by the skin.

- **COOL mode**

Switch to COOL mode to minimize the appearance of pores, lock in nutrients and reduce facial puffiness. Use this mode after applying lotion or moisturizer. It is recommended to use the device in this mode for 2 minutes on the left cheek, 2 minutes on the right cheek and 2 minutes on the forehead.

- When using EMS mode, move the surface of the massage head in an upward motion, as shown in the illustration. (Illustration 6)

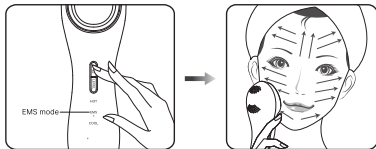
Figure 6

- When using HOT or COOL mode, lightly apply the massage head to the face or eye contour and gently massage as shown. (Illustration 7)

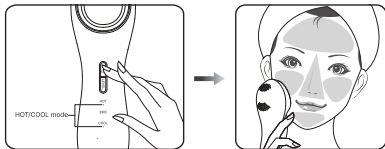
Figure 7

CLEANING AND MAINTENANCE

Clean the device with a clean, slightly dampened cloth before and after use.

- Always make sure that the device is OFF before cleaning.
- Do not immerse or clean the unit directly with water or other liquids. Clean with a clean, slightly dampened cloth.
- Do not use steel wool, abrasive cleaners, chemicals or caustic detergents to clean the unit.

STORAGE

The device has a built-in rechargeable lithium battery. Do not place or store the device in hot and humid places. If the device is not used often, charge the battery before storing it.

DEVICE SPECIFICATIONS

- Name: Hot & Cool Facial massager
- Model: BL-108
- Input: 5 V $\Rightarrow$ 1 A
- Battery type: lithium-polymer 3.7 V 1050 mAh battery.
- Performance standard: Q/XM 1020-2020

TROUBLESHOOTING

Problem	Possible cause	Solution
The device is not working	Battery (almost) depleted	Charge the device according to the instructions.
Note: If the device still does not work properly after applying the above solutions, stop using the device and contact the manufacturer.		

DISPOSAL

Lithium batteries contain substances that can pollute the environment. Do not dispose of this device with household waste. Instead, take it to an official recycling center to help protect the environment.

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