Feelo.

EASY Irygator do zębów

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Producent:
SHOP DENT Sp. z o.o.
Ul. Dziupli 12, 02-454
Warsaw, Poland
www.feelo.pl
Made in P.R.C

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Set elements

- 1 irrigator handle
- 1 water tank (300 ml)
- 2 standard nozzles
- 1 USB - DC V5 charging cable


Features and functions of the irrigator

- **Wide pressure range** – from 0.69 to 6.2 bar.
- **3 operating modes** – Clean (standard mode for daily oral irrigation), Soft (gentle cleaning for people with sensitive gums), Pulse (mode for simultaneous interdental cleaning and gum massage).
- **Efficient 1400 mAh lithium-ion battery** – one full charge lasts around 30 days in Soft mode, assuming the irrigator is used twice a day by one person.
- **Gradual pressure build-up** – the water stream gradually increases in strength over the first 3 seconds, making it easier to position the nozzle correctly.
- **Large 300 ml water tank** – allows for longer irrigation without frequent refills.
- **2-minute timer** – helps monitor irrigation time by automatically turning the device off after 2 minutes.
- **Slim and ergonomic handle** – fits comfortably in the hand for easy use.

Description of the irrigator

- a – Nozzle
- b – On/Off button
- c – Mode button
- d – Mode indicator
- e – Charging indicator

- f – Suction tube
g – Water tank
h – Nozzle release button

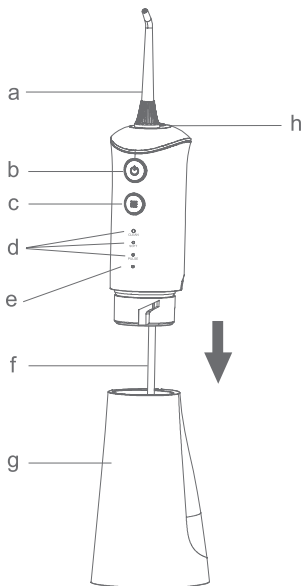Note

There may be residues of distilled water in the new unit (in the tank, suction tube, etc.), which was used for the performance test of the unit. It should not cause any concern when using the irrigator.

Charging the device (fig. 1)

1. Remove the cover on the back of the irrigator that protects the charging port, and make sure there is no water residue in the socket. Then firmly insert the DC 5V plug of the power cable into the irrigator's charging port.
2. Connect the other end of the cable (USB) to a USB charger and plug it into a power source.
3. The „CH“ indicator will blink while the device is charging. Once fully charged, the indicator will turn off.
4. After charging is complete, replace the cover on the charging port.
5. If the „CH“ indicator blinks during or after using the irrigator, it means the battery needs to be recharged. A full charge takes approximately 4 hours.

Important Note: Make sure you are charging with the appropriate voltage (5.0V) and current (1.0A). The power adapter must have output parameters that match the specifications listed on your irrigator's label. Use of a power supply with inconsistent parameters can lead to malfunction of the device. Failure to comply with this principles may result in damage to the device as well as void of the warranty.








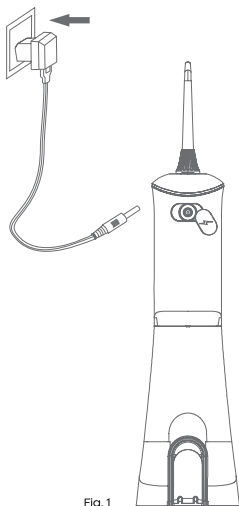
Fig. 1

Irrigator use

1. Attach the nozzle to the body of the device. (Fig. 2)
2. Open the tank cover and fill it with water at a temperature not exceeding 40°C, then close the cover. (Fig. 3)
3. Alternatively, you can unscrew the water tank, fill it with water, and then reattach it to the device. (Fig. 4)
4. Press the „Mode” button to select the cleaning mode. You can also switch modes while using the irrigator. (Fig. 5)
 - a. Clean – for standard cleaning
 - b. Soft – for gentle cleaning of sensitive gums

- c. Pulse – for interdental cleaning and gum massage
- 5. Hold the irrigator in an upright position with the nozzle tip pointed toward your teeth.
- 6. Slightly open your mouth to allow water to flow out freely.
- 7. Press the „On/Off” button to start irrigation. (Fig. 6)
 - a. Direct the water jet flow at a 90° angle to your teeth and gums.
 - b. Move the water jet slowly along the gumline and between your teeth.
 - c. Thoroughly clean interdental spaces as well as areas around braces, crowns, and bridges.
- 8. You can replace the nozzle if necessary. Turn off the irrigator first, then press the nozzle release button to remove and replace it. (Fig. 7)
- 9. After finishing irrigation, press the „On/Off” button again to turn off the device.

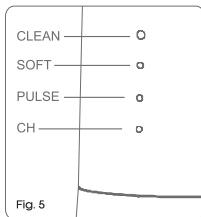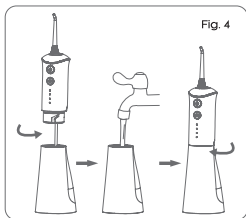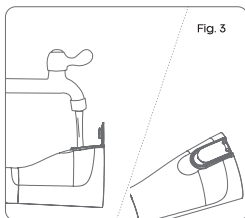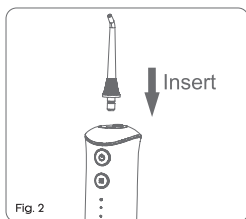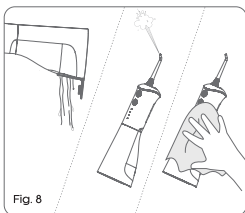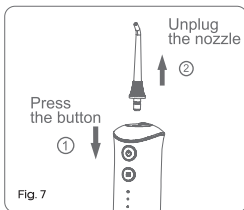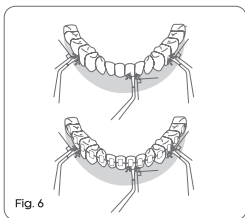After the irrigation is completed (fig 8.)

1. Open the tank cover and pour out any remaining water.
2. Briefly press the „On/Off” button to run the device and remove any residual water.
3. Turn off the irrigator, then close the tank cover.
4. Wipe the device with a dry cloth.

Note

To prevent bacterial growth, check for residual water in the unit. If the appliance is not going to be used for a long period of time, make sure it is properly dried after irrigation.

Cleaning and maintenance

Use only water or neutral cleaning agents to clean the irrigator. Do not use corrosive or abrasive substances (e.g., vinegar or citric acid solution, baking soda solution), as they may damage the device. Do not wash the irrigator with water hotter than 40°C. Keep the device away from high temperatures and avoid direct sunlight. The irrigator should be descaled once a month using only cleaners specifically designed for this purpose.

Main unit

1. Remove the dirt by wiping the device with a cloth dampened with soapy water, then dry it with a clean cloth.
2. Do not immerse the main part of the device in water.
3. Dry the device after cleaning.

Nozzle

1. Remove the nozzle from the irrigator and rinse it under running water after cleaning the interdental spaces.
2. Replace the nozzle every 6 months.

Water Suction Tube

1. Rinse with water and wipe with a soft cloth.
2. Do not bend, stretch, or twist the tube.

Water Tank





1. Rinse with water and wipe with a soft cloth.
2. If you are not going to use the device for an extended period, make sure the tank is dry before storing.

Disposal

The rechargeable battery inside your appliance is easy to remove for recycling (Fig.9). At the end of the product's life and prior to disposal, please remove the rechargeable battery from the handle by following the instructions below. During below operation, be sure

to protect your eyes, hands, fingers, and the surface on which you work.

1. Remove the water tank and unscrew the two screws beside of the suction tube. (Fig.10)
 2. Remove two rubber feet on the bottom the handle. (Fig.11)
 3. Insert a screwdriver into the slot located at the bottom of the handle and turn counterclockwise to release the bottom cap. (Fig.11)
 4. Open the bottom cover and put the inner out of the unit, then unscrew all the screws of the inner. (Fig.12)
 5. Cut off the battery line and pull out the batteries from the inner.
- If you have recycling questions, please contact your local waste management office.

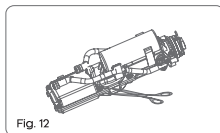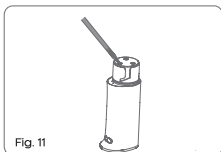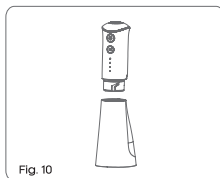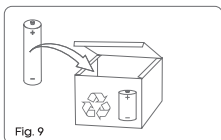Troubleshooting

Problem	Possible Cause	Solution
The irrigator does not work.	The irrigator has just been purchased or has not been used for more than 3 months.	Charge the device for at least 4 hours.
The irrigator only works for a few minutes, even just after it has been charged.	• The battery has reached the end of its service life. • The charging time was too short.	• Contact the Feelo service centre. • Charge the device for at least 4 hours.
The water pressure is very low.	• The nozzle is clogged. Hard water promotes limescale build-up, which can clog the tip. • The nozzle is de-formed.	Replace the nozzle with a new one.
Water does not flow out of the tip.	• The water tank is empty. • The device is tilted too much during use. • The irrigator has just been purchased or has been left unused for a long time.	• Fill the tank with water. • Hold the irrigator in an upright position. • Pull out the nozzle then flush nozzle hole on the body.

Specification

Model	Feelo Easy Irrigator do higieny jamy ustnej FE.WF.010.WHT.1
Input voltage	5V -- 1A
Rated input power	4.81±1.85W
Water pressure range	0,69-6,2 bar
Pulsation frequency	1400 pulses per minute

Water tank capacity	300 ml
Operation time	Approximately 30 days in Soft mode, assuming twice-daily use by one person.
Charging method	USB-DC 5V cable charging
Charging time	Approximately 4 hours
Battery type	Lithium-ion
Battery capacity	1400 mAh
Noise level	Up to 75 dB
Battery-low alert	Flashing charging indicator
Water resistance class	IPX7

User safety

1. Children and individuals who are unable to operate the device on their own should not use it. Keep the device out of reach of children.
2. Never attempt to modify, disassemble, or repair the device yourself.
3. Do not place or store the product in locations where it could fall into a bathtub or sink.
4. Do not reach for the device if it has fallen into water. Immediately disconnect it from the power source.
5. Do not use the product while bathing.
6. If you notice smoke, excessive heat, or unusual smells, turn off the device immediately and contact customer service.

Warnings! Read before use

- For domestic use only. 
- The irrigator can be used by children aged 8 years and over. Children and persons with limited physical, sensory or mental abilities or lacking experience and knowledge should use the device under the supervision of an adult supervisor, or should be adequately trained in the correct and safe use of the device. Children should not play with the device. Cleaning and maintenance should not be carried out by children.
- Use the device only with the charging cable provided in the set.
- Instructions on battery removal can be found in the „Disposal“ section on page 23 .
- Used equipment and battery should be taken to an appropriate collection point for the disposal of electrical and electronic waste.
- Do not fill the tank with water at a temperature higher than 40°C.
- If the charging cable becomes damaged, it must be replaced by the manufacturer, an authorized service center, or a qualified professional to avoid hazards.
- Keep the power cord away from hot surfaces.
- Do not use the device if the nozzle is missing.
- Do not share your nozzle with other users.
- Do not immerse the device in water.
- Use the device only for oral hygiene. The device is designed for cleaning interdental spaces, gums, tongue. Do not use it in any other way.
- Keep the nozzle out of reach of small children. The nozzle is small enough to be swallowed and may cause choking.
- Some types of mouthwash may damage the device, causing cracks in the housing or water tank. Do not add mouthwash to the water tank.
- If you have recently undergone oral surgery or have been diagnosed with serious periodontal disease, consult your doctor before use.
- Do not press the nozzle too firmly against teeth or gums, as this may cause injury.

DANGER WARNING

To reduce the risk of electric shock:

- Always unplug the device from the power source after use.
- Do not use while bathing.
- For more detailed instructions regarding storage and usage of the device and included accessories, consult the manufacturer.

To reduce the risk of fire, short circuit, or injury:

- Unplug the charger from the power outlet once charging is complete.
- This product is not a toy. Exercise extra caution when the irrigator is used by or around children.
- Use the product only as intended. We recommend using accessories primarily approved or recommended by the manufacturer.
- Do not use the charger if the power cord or plug is damaged. If the irrigator has been dropped, exposed to the elements, or is not operating properly, return it to the place of purchase.
- Do not pull, pinch in doors, or wrap the charger's power cord around sharp edges. Do not operate the irrigator while charging. Keep all device components away from heated surfaces.
- Do not unplug the charger by pulling on the cord. Always hold the plug itself.
- Do not handle the plug with wet hands.
- Never use the product if you are feeling drowsy or having trouble concentrating.
- Do not use the device outdoors or in environments where it may be exposed to aerosols or pure oxygen.
- Plug the device into a power outlet. Check its specifications.

TWO-YEAR WARRANTY

The manufacturer provides a warranty for a period of two years from the date of purchase. The warranty is valid in Poland and within the territories of official distributors. The warranty covers damages related to manufacturing defects or improper workmanship of the product. Keep the proof of purchase for the duration of the warranty period. During the warranty period, you can file a claim if the product was used correctly and according to its intended purpose.

THE WARRANTY DOES NOT COVER:

1. Surface damages resulting from normal use.
2. Damages resulting from:
 - a. Accidental drops or impacts, improper use, abuse or commercial use, neglect, modification of the device,
 - b. Improper maintenance, lack of care for the product,
 - c. Connection to an incorrect power source,
 - d. Repairs or attempts to repair by a party not authorized by the manufacturer. Any attempts to open or disassemble the device, its parts, or accessories will void the warranty.

Warranty claims must be supported by evidence. The claim must be made within the warranty period. The material covering the handle is delicate, so care should be taken during use. Scratches or abrasions on the surface of the handle will be treated as surface damages resulting from normal use. Therefore, they will not be covered by the warranty. The warranty applies to components that affect the functioning of the irrigator. Proper care of the device significantly increases its lifespan. Failure to follow the recommendations may result in the loss of the warranty.

Manufacturer:

SHOP DENT Sp. z o.o.
Ul. Dziupli 12, 02-454
Warsaw, Poland
www.feelo.pl
Made in P.R.C

CONTACT

Customer service: bok@shop-dent.pl
Complaint service:
reklamacje@shop-dent.pl
Product safety: kontakt@feelo.pl

Model: FE.WF.010.WHT.1

Input Voltage: 5 V -- 1 A

Rated Input Power: 4.81 ± 1.85 W

Water Resistance (irrigator handle): IPX7. Protected against the effects of immersion in water up to 1 meter deep for 30 minutes.

The irrigators are equipped with a built-in lithium-ion battery.

The product meets the requirements of the so-called „New Approach” Directives of the European Union (EU).

The packaging and instruction manual are recyclable.

This device is marked in accordance with the Polish Act on Waste Electrical and Electronic Equipment and the WEEE Directive 2012/19/EU with the symbol of a crossed-out waste bin. This marking indicates that this equipment, after its service life, cannot be disposed of with other household waste. The user is required to hand it over to those conducting the collection of waste electrical and electronic equipment. Those conducting the collection, including local collection points, stores, and municipal units, create an appropriate system that enables the disposal of this equipment. Proper handling of waste electrical and electronic equipment helps to avoid harmful health and environmental consequences arising from the presence of hazardous components and improper storage and processing of such equipment.

For domestic use only.

Symbol of Class II protection designation for electrical equipment.

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