Red Light Therapy Panel


For Face and Body

Quick Start Guide

FR Guide de démarrage rapide DE Schnellstartanleitung IT Guida rapida all'uso ES Guía de inicio rápido PT Guia de início rápido BG Ръководство за бърз старт CS Průvodce rychlým startem DA Hurtigstartvejledning EL Οδηγός γρήγορης εκκίνησης ET Kiirjuhend FI Pika-aloitusopas HR Brzi vodič za početak HU Gyors üzembe helyezési útmutató LT Greito paleidimo vadovas LV Ātrās sākšanas ceļvedis NL Snelstartgids NO Hurtigstartveiledning PL Przewodnik szybkiego startu RO Ghid de pornire rapidă SK Stručný sprievodca spustením SV Snabbstartsguide SC 快速入門指南 TC 快速入門指南 KR 빠른 시작 가이드 JP クイックスタートガイド

What's in the Box?

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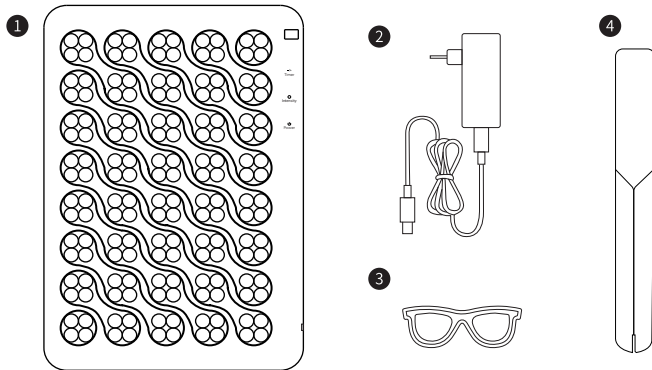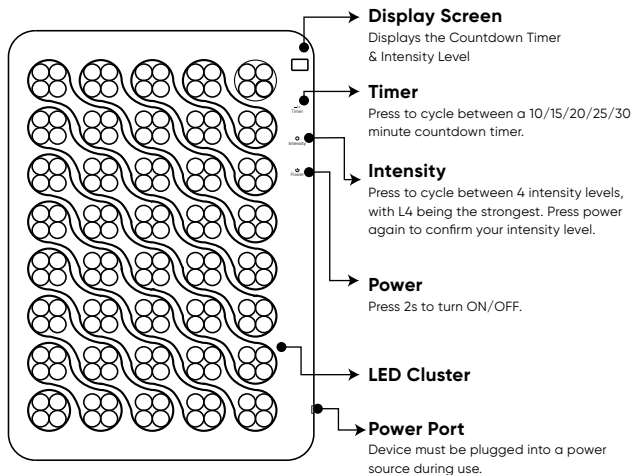
Displays the Countdown Timer & Intensity Level

Timer

Press to cycle between a 10/15/20/25/30 minute countdown timer.

Intensity

Press to cycle between 4 intensity levels, with L4 being the strongest. Press power again to confirm your intensity level.

Power

Press 2s to turn ON/OFF.

LED Cluster

Power Port

Device must be plugged into a power source during use.

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
Intensity Levels

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
L1

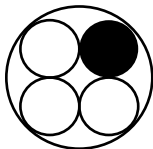
1 LED in each LED cluster is turned on.
1 LED allumée dans chaque groupe de LED.

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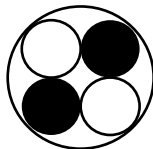
L3

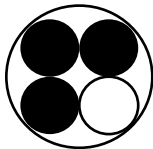3 LEDs in each LED cluster is turned on.
3 LEDs allumées dans chaque groupe de LED.

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
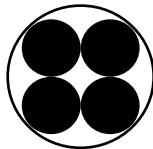
Red Light Therapy Panel must be connected to a power source during use.

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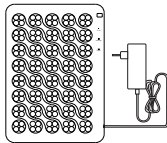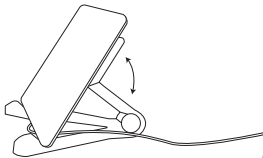Put on eye protection glasses, then press power button and select the intensity level and duration of your session.

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Important Things to Note


À savoir

This is not a medical device and cannot be used for professional treatment or as a substitute for professional medical advice. The actual effect of use varies from person to person.

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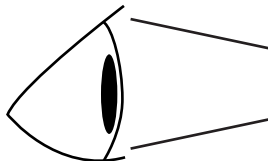
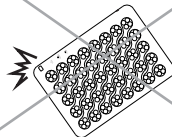
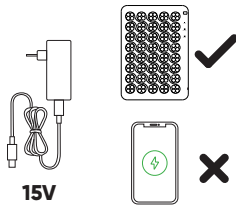
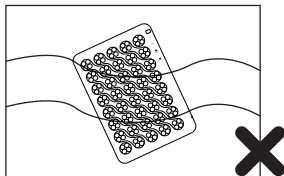
To clean, simply wipe with a damp cloth or with 70% alcohol.

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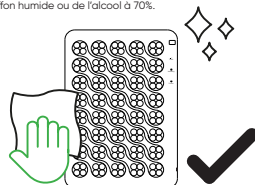EN WARNINGS This is not a medical device and cannot be used for professional treatment or as a substitute for professional medical advice. The actual effect of use varies from person to person. For first-time users, please start with a low intensity. If any abnormality occurs on your skin, please stop using it immediately and contact a professional doctor. This product contains an infrared spectrum. Do not use this device for more than 30 minutes at a time. Pregnant women, children, individuals with epilepsy, thyroid patients or those on thyroid medication should not use this product. Do not use with photoallergies, or while taking drugs that can cause photoallergies or photosensitivity. Avoid wounds, and areas with skin trauma. Individuals that utilize special medical implants, heart rate monitors or have special medical conditions should consult with a doctor before use. When this product is connected to a power source for charging, it can not be turned on. Please read this safety information carefully and keep this document for future reference. Help articles and customer support is always available at <https://nanoleaf.me/support>. Do not open or disassemble the device, as this will void the warranty and make the product unsafe for use. The flexible cable and power supply cannot be replaced by another third party product. If the cable is damaged, the entire controller and power supply should only be replaced with original equipment manufactured by Nanoleaf®. Using any non-approved replacement component will void the warranty and may present a serious safety hazard. This Nanoleaf® product is not designed, certified, nor intended for use in: emergency exits, emergency lighting, or other use cases involving critical systems for human safety; outdoor areas, regardless of weather conditions; regions or countries other than that in which the product was intended for and certified for sale in. When using this Nanoleaf® product: Only use the control provided with or specified by these instructions to control the lamp. The light sources or separate control gear are not intended for use in other applications combined. The light source of this device is not replaceable. When the light source reaches its end of life, the whole luminaire shall be replaced. Storage conditions; enclosed locations that limit the flow of air and/or trap heat which contributes to an environment temperature over 40 degrees Celsius.

CE INFORMATION NanoGrid Limited (Nanoleaf®) hereby declares that this device is in compliance with the essential requirements and other relevant provisions of Directive 2014/30/EU. The declaration of conformity is available at <https://nanoleaf.me/eudocs>. **CLEANING AND MAINTENANCE** Do not clean the device with liquid or sprays when it is connected to power. Clean with a soft dry cloth, and do not use harsh cleaning agents.

WARRANTY INFORMATION This product has a one year warranty from the date of purchase. Damages caused by time and use is not within the scope of the service guarantee. The warranty scope of this product does not include failures caused by any factors other than the defects of the product itself.

Disposal of your old product: Your product is designed and manufactured with high-quality materials and components, which can be recycled and reused. When this crossed-out wheeled bin symbol is attached to a product, it means that the product is covered by European Directive 2012/19/EU. Please find out about your local collection system for electrical and electronic products. Please act according to your local rules and do not dispose of your old products with your normal household waste. The correct disposal of your old product will help to prevent potential negative consequences for the environment and human health.

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Nanoleaf Europe SAS