




medivon.

**User manual**

PL | EN | DE | FR | ES | PT | IT | CZ | SK  
SL | EE | LT | LV | HU | RO | RU | NB | FI  
EL | HR | SV | DA | NL | TR | BG | AR

## Because You're Worth Every Moment

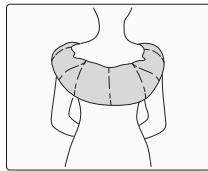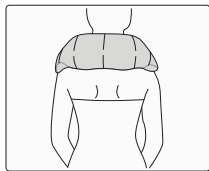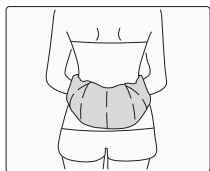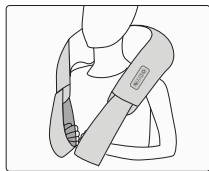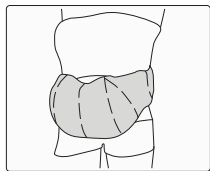One massager. Endless ways to use it.

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Before use, make sure the massager is fully charged. The device has a built-in battery that lasts for about 20 massage sessions. For charging, use only the included USB-C cable, connected to a proven power source. The device should not be used while charging.

The massager is operated using a single button:

Press and hold for 2 seconds - turn on

Short presses - change mode (2 intensity levels)

Press and hold for 3 seconds - turn off

The device is intended for use on the neck, but can also be used on other parts of the body according to the infographic at the beginning of this instruction manual. Use a maximum of several times a day, each time for no longer than 15 minutes (the device has a built-in timer that will turn the massager off after 15 minutes).

Contraindications:

Do not use the massager in the case of open wounds, skin inflammation or fresh injuries. Do not use if you have a pacemaker, circulation problems or other conditions that require medical consultation.

The device is not intended for children or people with limited ability to use it independently without supervision.

Do not use while sleeping or driving.

Store the massager in a dry place, away from sources of heat and moisture. Wipe with a soft, dry cloth when necessary.

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Used electrical and electronic equipment and batteries and accumulators must not be mixed with other waste. This is indicated by the crossed-out bin symbol on the product packaging. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution. Proper handling of used electrical and electronic equipment and batteries and accumulators in the household contributes positively to the reuse and recovery, including recycling, of used equipment. Used electrical and electronic equipment and batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and has performed the same functions as the equipment sold. A distributor, when delivering household equipment to a purchaser, is obliged to take back used household equipment free of charge at the place of delivery of such equipment, provided that the used equipment is of the same type and has performed the same functions as the equipment delivered. A distributor running a retail outlet with a sales area within the meaning of art. 2 item 19 of the Act of 27 March 2003 on spatial planning and development of at least 400 m<sup>2</sup> devoted to the sale of household equipment, is obliged to accept, free of charge, in that outlet or in its immediate vicinity, used household equipment, none of the external dimensions of which exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market after 13 August 2005

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The manufacturer hereby declares that the Medivon Motion X device complies with, Directive 2014/30/EU of the European Parliament and of the Council of February 26, 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65/EU of the European Parliament and of the Council of June 8, 2011 on restrictions on the use of certain hazardous substances in electrical and electronic equipment. The full content of the declaration of conformity is available at the following internet address: [www.medivon.pl/deklaracje](http://www.medivon.pl/deklaracje)

Importer: Armare S.A.  
ul Grzybowska 87, 00-844 Warszawa, Poland.  
Made in China.  
SN:06255904119286407  
Marka: Medivon  
Model: Medivon Motion X  
email: [serwis@medivon.pl](mailto:serwis@medivon.pl)